


AIR COMPRESSION RECOVERY SYSTEM

Model: FE-7221

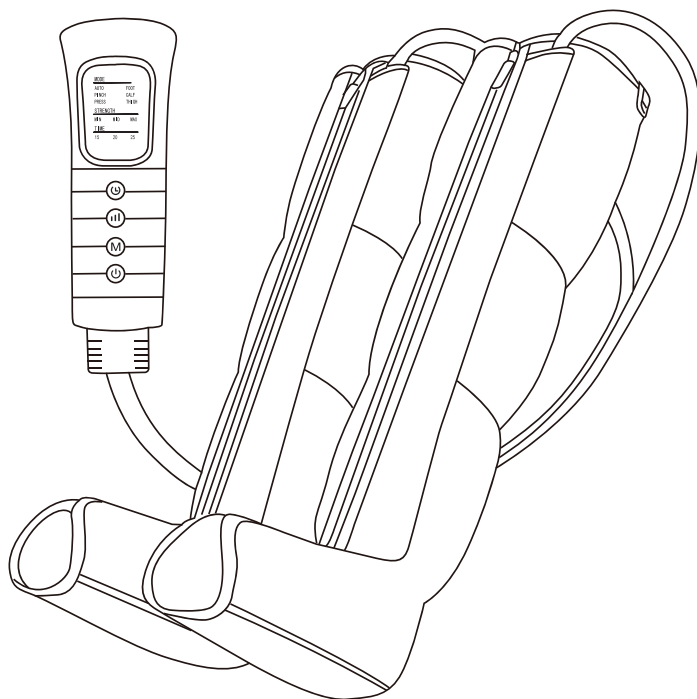

USER MANUAL

Please read the user manual carefully prior to use

Content


Safety Precautions	01-02
Package List	03
Control Unit	03
Button Introduction	04-05
Massage Sleeve	05
Operation Guide	06
Specifications	07
Product Cleaning and Maintenance	07-09
Troubleshooting and Treatment	09

Safety Precautions

- 01.** This product is for personal use only and is not a medical instrument or for commercial activities.
- 02.** Don't scratch, damage, process, excessively bend, pull or twist the powercord of power adapter. Otherwise, it may cause fire or electric shock. When unplugging, hold onto the plug and then unplug it.
- 03.** Don't use other power adapter but the original one.
- 04.** Don't plug or unplug the power adapter with wet hands.
- 05.** Don't touch the product with sharp objects to prevent damage to the product parts.
- 06.** Avoid operating the product under blankets or pillows and allow for proper heat dissipation to prevent overheating, which could lead to fire, electric shock, or injury.
- 07.** Don't use the product in areas with corrosive products or oxygen administration.
- 08.** Don't use the product in high humidity or high temperature environments that may cause electric shock or other safety hazards.
- 09.** Don't use the product while driving and obey driving laws for safety.
- 10.** Don't walk around when you are using this product or wearing the sleeves.
- 11.** Immediately stop using the product if you experience discomfort or severe pain during use.
- 12.** Don't use the product if it has a damaged cord or plug, is not working properly, has been dropped / damaged, or has come into contact with water.
- 13.** Stop using the product if it is not working properly.
- 14.** Supervision is required when children or individuals lacking knowledge use the product.

15. Keep it out of the reach of infants, children and people without the ability to use it independently.
16. Use the product only for its intended purpose as described in the manual.
17. It is forbidden to remodel, disassemble or repair the product without permission.
18. After use, please place the massager in a dry place, away from humidity and direct sunlight.
19. Unplug the power adapter when you are not using it.
20. Please keep the product clean and use a soft and dry towel during maintenance.
21. The product is designed to be used by only one person at a time.
22. Consult a doctor before using this product if you have any of the following conditions:
 - During pregnancy, menstruation, or shortly after delivery
 - Implanted medical electronic devices, such as a cardiac pacemaker
 - Blood circulation disorders or thrombosis
 - Osteoporosis
 - Skin infections
 - Heart disease patients
 - Wounds or surface diseases on the body
 - Suffering from malignant tumors
 - Physical disabilities, sensory or neurological deficits, or lack of experience and common sense (including children with mental illness or physical paralysis)
 - Receiving medical treatment, or have a contraindication requiring the use of any medical control unit.
 - Have bone fractures or dislocations at or near the site of application
 - Having serious peripheral neuropathy dysfunction or sensory disturbance caused by diabetes

Package List

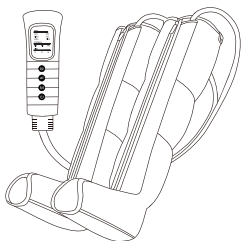

Air Compression
Leg Massager x 1




Adapter x 1



User Manual x 1

Control Unit

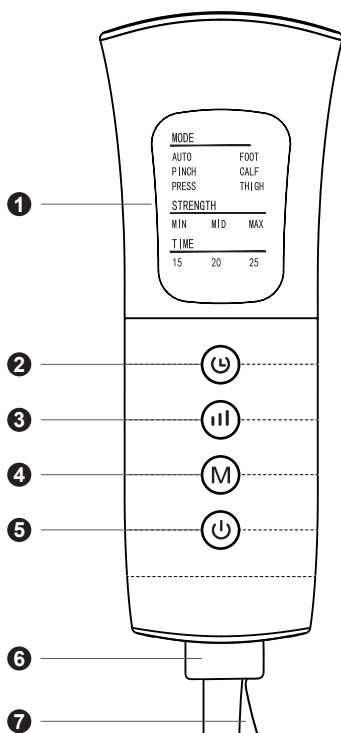

- ❶ Display Screen
- ❷ Timer Button
- ❸ Intensity Button
- ❹ Mode Button
- ❺ Power Button
- ❻ Air Hose Connect Port
- ❼ Power Port

Button Introduction


1. Power Button

Briefly press to turns on/off the device.


2. Mode Button

In the power-on state, short press to switch modes. There are 6 massag modes. please check the difference as below:

Mode 1 — Auto Mode

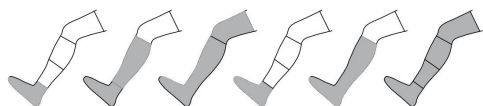

Mode 2 — Pinch Mode

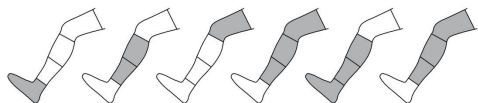


Mode 3 — Press Mode




Mode 4 — Foot Mode



Mode 5 — Calf Mode

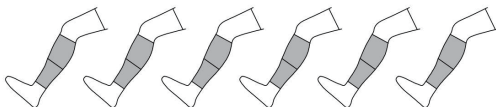

Mode 6 – Thigh Mode

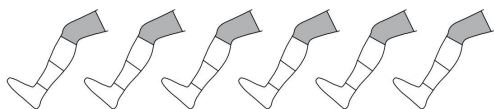


3. Intensity Button

There are 3 intensity levels, low (100 mmHg), medium (130 mmHg) and high (160 mmHg), short press to switch intensity mode.

4. Timer Button


Choose from three timers: 15 mins/20 mins/25 mins. Short press to switch timer button. Press the timer button to adjust the usage time, and the device will automatically power off when the timer elapses.

Massage Sleeve

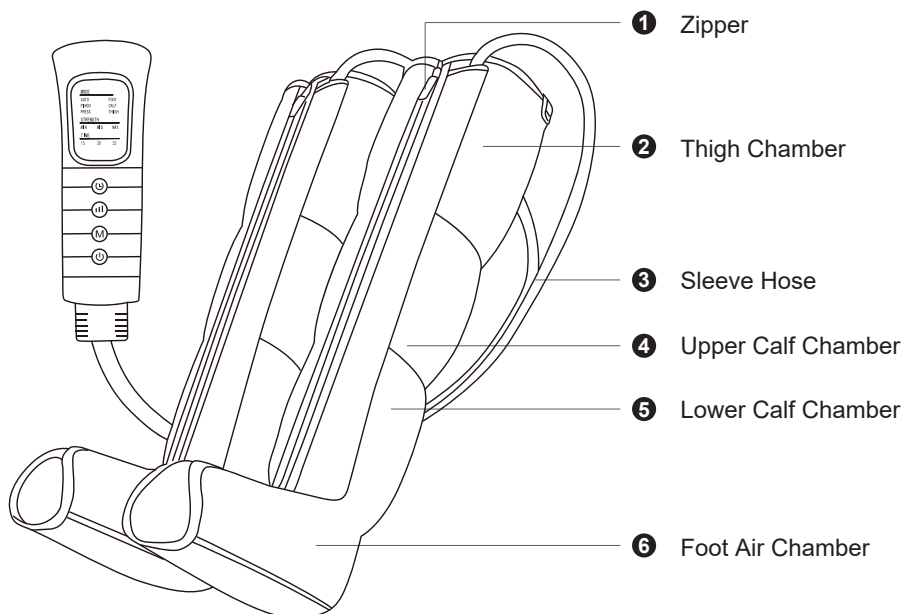

Operation Guide

Please read safety precautions carefully before using the product.
In order to get better massage, please use it according to the instruction.

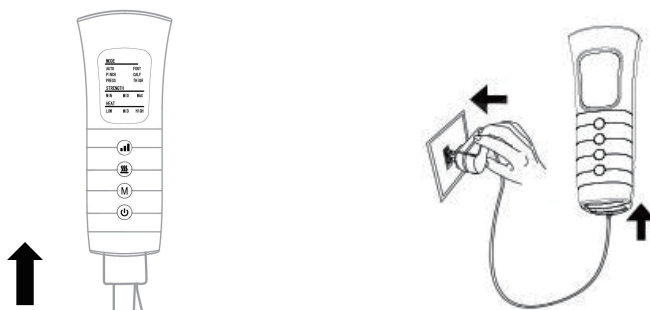

01. Connect the hose to the control unit.
* Two sleeves need to be connected to make it work.
02. Connect the adapter well with the outlet and the control unit, the screen will flash display for 1 second.
03. Put the sleeves on your leg properly.
04. Find a comfortable position sitting, reclined, or lying down. Press the power button to start.
05. You can casually switch the mode, intensity levels and timers to find the one that suits yourself, and then enjoy the massage (The product features a memory function that retains the last used pressure, mode and timing settings upon startup, preserving them as defaults for the next use.).
06. The device will automatically shut off after the treatment time is up, you can press the button if you wanna enjoy more or end the massage earlier at any time. Recommend to use 2 times a day (i.e.40 mins in one day)
07. After use, unplugging the air hose can lead to a rapid deflation.

Specifications

Product Name	Air Compression Recovery System
Mode	FE-7221
Color	Black
Pressure Levels	100-130-160 mmHg
Timer	15 / 20 / 25 Mins
Input	12V/1A
Rated Power	12W

Product Cleaning and Maintenance

Cleaning

Controller

Wipe the display with a soft cloth

Fabric

The product does not support washing, please use alcohol sterilization or clean wet towel wipe



Storage requirement

Wipe the display with a soft cloth



Discard

Dispose of discarded batteries in accordance with local regulations or environmental regulations

Use Environment

	Use Environment	Storage and Transportation
Temperature Range	41 – 104 °F	14 – 140.4 °F
Humidity Range	5% ~ 93% RH	5% ~ 90% RH
Barometric Range	86kPa ~ 105kPa	86kPa ~ 105kPa

Product Failure Analysis and Treatment

Phenimenon	Possible Causes	Solution
No Connection Adapter	No Connection Adapter	Connection Adapter
Product can not start	Connection Adapter	Connection Adapter
	If the timing is set, the machine will shut down antomatically when the time is up	Restart

Troubleshooting and Treatment

Phenomenon	Possible Causes	Solution
Doesn't work and screen is off	1 Power is not turned on. 2 Adapter is not securely connected. 3 Faulty wall outlet.	1 Press the power button to turn the control unit on. 2 Make sure the power adapter is connected well to the control unit and the wall outlet. 3 Check if the wall outlet works.
The screen is on but the sleeves doesn't work.	1 Only connect one of the sleeves with the air hose. 2 Air hoses are not securely connected.	1 It works only when 2 sleeves are connected to the air hose. 2 Please check if the air hoses are inserted correctly.
Break off during operation or stop pumping	1 Air hoses are not securely connected 2 The power adapter falls off. 3 The massage session time is up	1 Please check if the air hoses are inserted correctly. 2 Make sure the power adapter is connected well. 3 The device will shut off automatically when the pre-set timing is up press the power button to start a new session if needed.

