


SUNVIVI

ZRN017

WEIGHTED HEATING PAD USER MANUAL

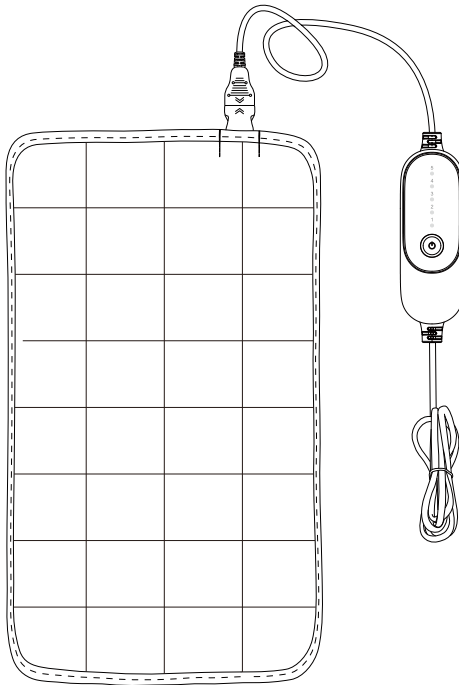

Model No.:NF10B

For your safety and correct use of the product,
please read the instructions carefully before use.



707-392-8405



www.sunvivi.com



service@sunvivi.com

INTRODUCTION

Thank you for purchasing the Heating Pad. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and order number.

INCLUDED













- 1 x Weighted Heating Pad
- 1 x Hand-held controller
- 1 x User manual

(Note: please check first to make sure the product is complete without any damage. If you have any doubt, please do not use and contact us.)

Product Specifications

Product Model	NF110B
SKU	ZRN017
Name	Weighted Heating Pad
Input voltage	120V 60Hz
Rated power consumption	110W
Heating level	5 Levels
Auto-off timing settings(hour)	2 Hours
Product dimensions	12in x 24in
Controller Model	RPHP-NFI-10B
Material	100% Polyester

Symbol Explanation

	Read all the instructions.		Gentle wash cycle in water at 30°C(86 °F).
	Do not puncture.		Do not bleach.
	Do not use if folded or creased.		Do not tumble dry.
	Not suitable for very young children (0-3 years).		Do not iron.
	Dispose of the packing in an eco-friendly manner.		Do not dry clean.
	WARNING: Warning of risks of injury or health hazards.		
	CAUTION: Safety information about possible damage to appliance / accessories.		

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER—TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.

01. READ ALL INSTRUCTIONS.
02. DO NOT USE WHILE SLEEPING.
03. DO NOT USE ON AN INFANT.
04. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.
05. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
06. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.
07. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
08. DO NOT SIT ON OR CRUSH PAD – AVOID SHARP FOLDS.
09. NEVER PULL THIS PAD BY THE SUPPLY CORD.
10. DO NOT USE THE CORD AS A HANDLE.
11. CAREFULLY EXAMINE COVER BEFORE EACH USE. DISCARD THE PAD IF COVERING SHOWS ANY SIGN OF DETERIORATION.
12. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
13. ENSURE THIS PAD IS USED ONLY ON A STANDARD 120V OUTLET.
14. NEVER LEAVE THIS PAD CONNECTED TO AN ELECTRICAL SUPPLY CIRCUIT WHEN NOT IN ACTUAL USE.
15. WHEN NOT IN USE THIS PAD SHOULD BE KEPT IN A COOL, DARK, DRY PLACE – PREFERABLY IN THE ORIGINAL CARTON.
16. THIS HEATING PAD HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER) AS A SAFETY FEATURE. THIS PLUG WILL FIT INTO A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY INTO THE OUTLET, REVERSE THE PLUG.

- IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN.
DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.
17. READ ALL INSTRUCTIONS BEFORE USING THE PRODUCT.
 18. THE PRODUCT IS DESIGNED ONLY TO WARM UP THE HUMAN BODY.
 19. DO NOT USE ON PETS.
 20. DO NOT USE ON PEOPLE INSENSITIVE TO HEAT OR VULNERABLE PEOPLE UNABLE TO REACT TO OVER-HEATING.
 21. THIS PRODUCT SHOULD NOT BE USED ON VERY YOUNG CHILDREN (0-3 YEARS OLD).
 22. THIS PRODUCT CAN BE USED ON CHILDREN (3-8 YEARS OLD) ONLY UNDER SUPERVISION FROM ADULTS AND IF THEY HAVE BEEN GIVEN ADEQUATE INSTRUCTIONS ON HOW TO USE IT SAFELY.
 23. THIS PRODUCT CAN BE USED BY CHILDREN OVER THE AGE OF 8, PEOPLE WITH REDUCED PHYSICAL, SENSORY OR MENTAL SKILLS AND PEOPLE LACKING EXPERIENCE OR KNOWLEDGE IF THEY ARE SUPERVISED AND HAVE BEEN INSTRUCTED ON HOW TO USE IT SAFELY, FULLY AWARE OF THE CONSEQUENT RISKS OF USE.
 24. CHILDREN SHOULD NOT PLAY WITH THE PRODUCT.
 25. CLEANING AND MAINTENANCE MUST NOT BE PERFORMED BY CHILDREN UNLESS SUPERVISED.
 26. THE PRODUCT IS NOT DESIGNED FOR MEDICAL USE IN HOSPITALS.
 27. DO NOT PUNCTURE; DO NOT FASTEN IN PLACE USING PINS OR OTHER METALLIC OBJECTS.
 28. DO NOT USE WHEN FOLDED OR BUNCHED UP.
 29. DO NOT USE IF WET.
 30. THE PRODUCT SHOULD ONLY BE USED IN CONJUNCTION WITH THE SPECIFIED CONTROLLER.
 31. THE ELECTRICAL AND MAGNETIC FIELDS EMITTED BY THE PRODUCT MAY INTERFERE WITH THE FUNCTIONS OF MEDICAL DEVICES LIKE A PACEMAKER. PLEASE CONSULT YOUR DOCTOR BEFORE USE IF YOU ARE EQUIPPED WITH A MEDICAL DEVICE.

32. DO NOT PULL, TWIST OR MAKE SHARP BENDS ON THE CORDS.
33. PLEASE CHECK THE PRODUCT FREQUENTLY FOR SIGNS OF DAMAGE OR WEAR AND TEAR. IF SUCH SIGNS ARE EVIDENT, IF THE PRODUCT HAS BEEN USED INCORRECTLY OR IF IT NO LONGER HEATS UP, IT MUST BE CHECKED BY THE MANUFACTURER BEFORE BEING SWITCHED ON AGAIN. THERE ARE NO USER SERVICEABLE PARTS.
34. IF THE MAINS CONNECTION CORD IS DAMAGED, IT MUST BE DISPOSED OF; IF THE CORD CANNOT BE REMOVED FROM THE PRODUCT, THE PRODUCT SHOULD BE DISPOSED OF ALONG WITH THE CORD.
35. THE ELECTRONIC COMPONENTS IN THE CONTROLLER WOULD WARM UP WHEN THE PRODUCT IS IN USE. FOR THIS REASON, THE CONTROLLER MUST NEVER BE COVERED WHILE IN USE.
36. DO NOT USE IN AN ENVIRONMENT WITH OVERLY ABUNDANT OXYGEN.
37. UNPLUG WHEN NOT IN USE AND NEVER LEAVE UNATTENDED WHILE USING, ESPECIALLY WHILE CHILDREN ARE PRESENT.
38. NEVER PULL BY THE CORD OR USE THE CORD AS A HANDLE; DO NOT PINCH THE CORD.
39. USE THE PRODUCT ONLY AT RATED VOLTAGE AS MARKED ON THE PRODUCT. IT'S RECOMMENDED NOT TO USE WITH GENERATORS AND POWER CONVERTERS.
40. DO NOT USE WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS.
41. IF THE PRODUCT IS PROVIDED WITH MAGNETIC CLOSURE, DO NOT USE IF YOU HAVE A PACEMAKER OR OTHER MEDICAL DEVICES SENSITIVE TO MAGNETS.
42. HOUSEHOLD USE ONLY.
43. IF THE SUPPLY CORD IS DAMAGED, IT MUST BE REPLACED BY THE MANUFACTURER, ITS SERVICE AGENT OR SIMILARLY QUALIFIED PERSONS IN ORDER TO AVOID A HAZARD.


44. THE INSTRUCTIONS FOR WASHABLE APPLIANCES FITTED WITH A NON-DETACHABLE CORD SHALL STATE THAT THE SWITCH OR CONTROL UNIT MUST NOT BE ALLOWED TO GET WET DURING WASHING AND THAT DURING DRYING THE CORD MUST BE POSITIONED TO ENSURE THAT WATER DOES NOT FLOW INTO THE SWITCH OR CONTROL.
45. WHEN THE APPLIANCE IS LIKELY TO BE USED FOR AN EXTENDED PERIOD, CONTROLLER SHOULD BE ADJUSTED TO THE LOWEST SETTING FOR CONTINUOUS USE.
46. DO NOT DRY-CLEAN THIS HEATING PAD. CLEANING SOLVENTS MAY HAVE A DETERIORATING EFFECT ON THE INSULATION OF THE HEATING ELEMENT.
47. WHEN NOT IN USE, STORE AS FOLLOWS: STORE THE PRODUCT ONLY AFTER IT COOLS DOWN COMPLETELY AND IT IS DRY COMPLETELY. LOOP THE CORD LOOSELY. STORE IN AN AIRY, DRY ENVIRONMENT.
48. WHEN STORING THE APPLIANCE, ALLOW IT TO COOL DOWN BEFORE FOLDING;
49. DO NOT CREASE THE APPLIANCE BY PLACING ITEMS ON TOP OF IT DURING STORAGE;
50. EXAMINE THE APPLIANCE FREQUENTLY FOR SIGNS OF WEAR OR DAMAGE. IF THERE ARE SUCH SIGNS, IF THE APPLIANCE HAS BEEN MISUSED OR DOES NOT WORK, RETURN IT TO THE SUPPLIER BEFORE SWITCHING IT ON AGAIN;
51. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS;
52. DO NOT USE IF WET;
53. THIS APPLIANCE MUST NOT BE USED BY PERSONS INSENSITIVE TO HEAT AND OTHER VERY VULNERABLE PERSONS WHO ARE UNABLE TO REACT TO OVERHEATING;
54. CHILDREN UNDER THE AGE OF THREE ARE NOT TO USE THIS APPLIANCE DUE TO THEIR INABILITY TO REACT TO OVERHEATING.
55. THE INSTRUCTIONS SHALL STATE THAT THE APPLIANCE IS NOT TO BE USED BY YOUNG CHILDREN OVER THE AGE OF THREE UNLESS THE CONTROLS HAVE BEEN PRE-SET BY A PARENT OR GUARDIAN, AND UNLESS THE CHILD HAS BEEN ADEQUATELY INSTRUCTED ON HOW TO OPERATE THE CONTROLS SAFELY.

56. THIS APPLIANCE CAN BE USED BY CHILDREN AGED FROM 8 YEARS AND ABOVE AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.
57. IF THE SUPPLY CORD IS DAMAGED, IT MUST BE REPLACED BY THE MANUFACTURER, ITS SERVICE AGENT OR SIMILARLY QUALIFIED PERSONS IN ORDER TO AVOID A HAZARD.
58. THE INSTRUCTIONS FOR WASHABLE APPLIANCES FITTED WITH A NON-DETACHABLE CORD SHALL STATE THAT THE SWITCH OR CONTROL UNIT MUST NOT BE ALLOWED TO GET WET DURING WASHING AND THAT DURING DRYING THE CORD MUST BE POSITIONED TO ENSURE THAT WATER DOES NOT FLOW INTO THE SWITCH OR CONTROL.
59. THE INSTRUCTIONS FOR HEATING PAD SHALL STATE THE SUBSTANCE OF THE FOLLOWING:
 - THE APPLIANCE SHOULD NOT BE USED ON AN ADJUSTABLE BED; OR.
 - IF THE APPLIANCE IS USED ON AN ADJUSTABLE BED, CHECK THAT THE BLANKET AND CORD DO NOT BECOME TRAPPED OR RUCKED, FOR EXAMPLE IN HINGES.
60. WHEN THE APPLIANCE IS LIKELY TO BE USED FOR AN EXTENDED PERIOD, CONTROLLER SHOULD BE ADJUSTED TO THE LOWEST SETTING NO. 1 FOR CONTINUOUS USE.
61. THE PRODUCT CAN ONLY USE WITH RPHP-NFI-10B CONTROLLER. SAVE THESE INSTRUCTIONS.

Note: Do Not Plug Into the Wall Outlet Until You Have Connected the Controller Cord to the Blanket to Make Sure It Works Properly.


PRODUCT OPERATION GUIDE

STEP 1

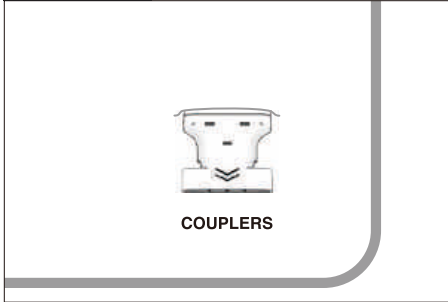


Locate the controller cord coupler on the Heating Pad.


STEP 2

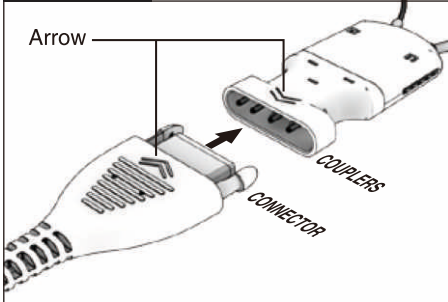

Connect Power Cord to the heated Heating Pad by matching the arrow on the connectors to the arrow on the couplers. Once the arrows are aligned, push power cord connectors into the couplers until they are securely connected.

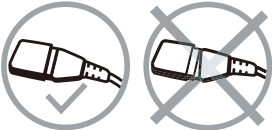


STEP 3

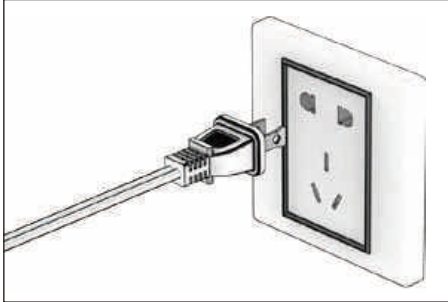

Plug Controller into the wall outlet(120V,60Hz)Heating Pad has one hand held controller.

Controller Diagram

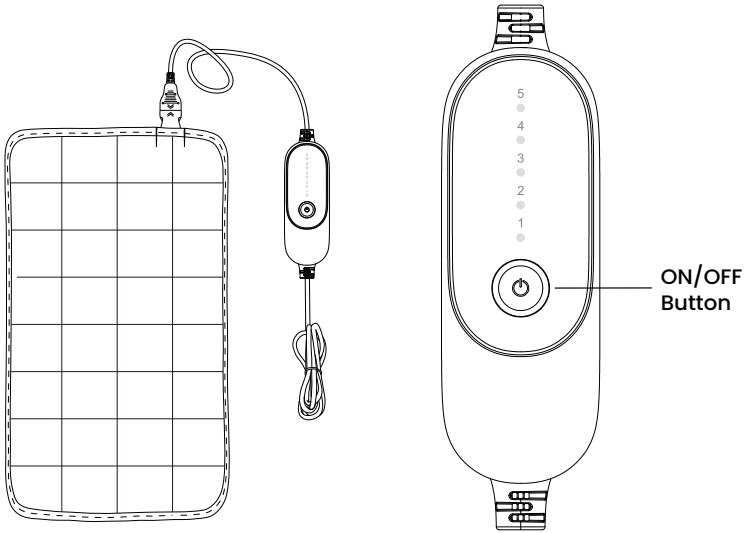

Product Operation Guide

1. Connect the product to the power supply. Press the power button to turn on the product and heating level 1 will be activated by default, which would be displayed on the screen.
2. Press the heating button to choose heating level from level 1 to level 5 by turn.
3. The heating pad will switch off automatically after 2 hours of consistent use.
4. Press the power button to turn off the product, and disconnect it from the power supply.

CLEANING AND MAINTENANCE WARNING

1. Before cleaning, always remove the mains plug from the outlet first.
2. Disconnect the controller from the product.

CLEANING

1. Soak the product for 15 minutes in mild soap and cool water.
2. Machine wash on gentle cycle in water at 30°C/85°F for 5 minutes or hand wash.
3. Rinse with fresh, cool water.
4. Machine dry on gentle, cool cycle or hang to dry.

(Immediately after washing, reshape the product to its original dimensions and spread it out flat on a clothes horse or clothesline.)

WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER! Do not iron your heating pad.


1. Disconnect the cord from pad by firmly pulling the control connector away from the pad connector. Do not pull by the cord.
2. Machine wash pad in cold water on gentle cycle or hand wash; tumble dry on low heat or line dry. Remove promptly from the dryer.
3. Do not plug in pad and pad connector until they are completely dry.

Questions & Answers

Problem	Cause	Solution
<p>The product does not produce enough heat or does not work.</p>	<p>Heating products provide safe heat and protect consumers' skin within temperature range based on international safety standards, which helps avoid low temperature burns. Therefore, the temperature setting is for keeping warm instead of getting hot.</p>	<p>If you don't feel heat, you can fold the product and test it for 30 minutes to see whether it works or not (But please make sure that the product is placed flat without being folded for normal use). If you still can't feel the heat after the test, please contact the customer service.</p>
<p>The product does not work upon arrival.</p>	<p>There exist product malfunctions.</p>	<p>When receiving the product, please take it out from the packing and place it flat (do not fold). Check the coupler and mains plug have been connected properly, then switch on to check whether the product works normally or not. For products unable to work normally, please contact our customer service.</p>
<p>The product does not heat uniformly.</p>	<p>The product is folded or the heating wire is coiled during the use;</p>	<p>Please place the heating pad flat for using. The temperature of the body areas that contact the heating zones is higher than the areas do not, which is a normal phenomenon.</p>
<p>The product no longer works after being used for a period of time.</p>	<p>The product is not used or washed properly.</p>	<p>Please contact our customer service.</p>
<p>If the above suggestions could not solve the problems, this might result from the permanent switch-off protection via the safety system. In this case, please contact our customer service.</p>		

POLARIZED PLUG:

Plug power cord into a 110–120 volt wall outlet. Do not attempt to force the plug. This appliance has a polarized plug (one blade is wider than the other). To reduce the risk of electrical shock, this plug is intended to fit only one way into an outlet. If the plug does not fit in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way. If the plug fits loosely into the AC outlet or if the ac outlet feels warm do not use that outlet.

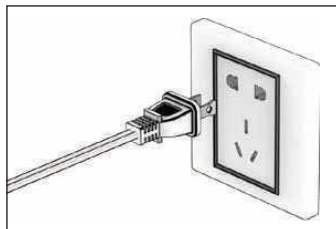

FCC CERTIFIED:

This product complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This product may not cause harmful interference
2. This product must accept any interference received, including interference that may cause undesired operation

Note: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or televisions reception, which can be determined by turning the product off and on, the user is encouraged to correct the interference by one or more of the following measures:

- ***Reorient or relocate the receiving antenna.**
- ***Increase the separation between the product and receiver.**
- *Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- *Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

WARRANTY AND CUSTOMER SUPPORT

If you need help regarding product issues, after-sales service, or other relevant problems, please provide the original order number through your purchase channel for further communication. We provide a 1-year limited warranty.

PRODUCED BY:

Ningbo Royal Peace Household Products Co.,Ltd
No. 269 Yindong South Road, Yinzhou Economic Development ZONE,
Yinzhou District,
Ningbo, Zhejiang, China 315145



Made in China



SUNVIVI[®]



Any Problem? Let Us Know.

 707-392-8405  www.sunvivi.com  support@sunvivi.com