

Portable Toilet User Manual



Read Carefully Before Use
Keep For Future Reference

Safety Information


Warning!

- **DO NOT** use this toilet on inclined or uneven ground.
- **DO NOT** use this toilet if it is damaged or malfunctions. Repair or replace worn or broken parts before further use.
- **ALWAYS** use the seat when sitting on this toilet. Sitting off-center or directly on the edge of the flush tank may cause damage or personal injury.
- **DO NOT** stand on this toilet at any time.
- **DO NOT** lean on the backrest of this toilet.
- **DO NOT** place more than 150 kg (330 lb.) on this toilet at any one time.

Specifications

Overall Dimensions	35×43×29.5 cm	13.8×16.9×11.6 in.
Seat Inner Diameter	25.5×22 cm	10×8.7 in.
Clean Tank Volume	12 L	3.17 gal.
Waste Tank Volume	12 L	3.17 gal.
Weight Capacity	150 kg	330 lb.
Net Weight	4 kg	8.8 lb.

Product Diagram

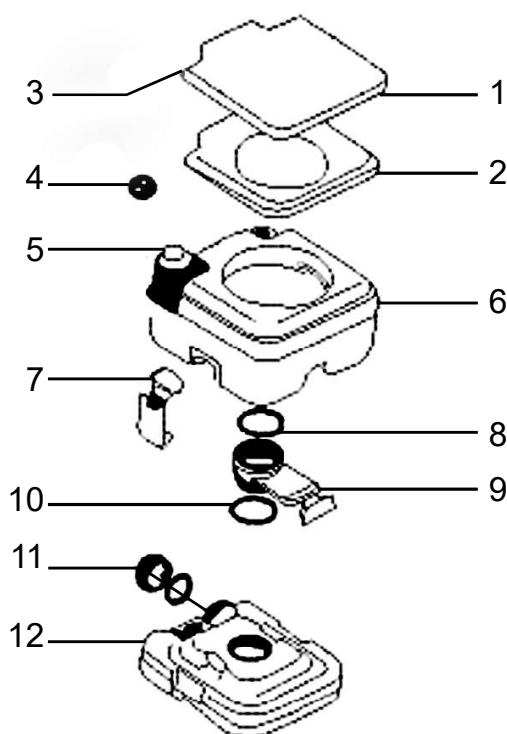

No.	Name
1	Seat Cover
2	Seat
3	Hinge
4	Water Storage Cap
5	Flush Button or Handle
6	Flush Tank (Bowl)
7	Detachable Latch
8	O-Ring
9	Slide Valve
10	O-Ring
11	Outlet Cap
12	Waste Tank

Assembly

1. Place the waste tank on a firm and level surface.
2. Mount the main body onto the waste tank's slide valve.
3. Secure the main body by fitting its latches into the corresponding slots of the waste tank.

Operation

1. Remove the water storage cap and fill it with clean water. **DO NOT** exceed the "Max. Water Level" line. Replace the cap, screwing it into place.
2. Open the slide valve by pulling the rod at the back of the waste tank.
3. Use the toilet in standard fashion.
4. Rinse the bowl. For models equipped with a flush button, press it as many times as needed. For models equipped with a flush handle, pull it up to flush the bowl and push it down completely to stop flushing.
5. Close the slide valve by pushing its rod back in.

Maintenance

Emptying the Waste Tank

The waste tank should be emptied when waste accumulates to about 2 cm (0.8 in.) below the slide valve.

1. Close the slide valve completely by pushing its rod in.
2. Remove the latches and remove the main body.
3. Carry the waste tank to a permanent toilet or other location designated as a safe disposal site for human waste. Unscrew its outlet cap and pour out accumulated waste into the facility.
4. Rinse the waste tank with clean water. Pour out this water into the same disposal facility as in Step 3.
5. Replace the outlet cap of the waste tank, screwing it into place.
6. Reassemble the toilet as before.

Cleaning

The exterior of this toilet should be cleaned periodically. Do not use abrasive cleaners or caustic chemicals.

- Clean the external surfaces of this toilet with water and any mild detergent. Rinse clean.
- To clean the toilet bowl, first make sure it is well attached so nothing will leak. Open the slide valve and clean all surfaces with any toilet cleaner and a toilet brush. Spray it clean or rinse it using the flush button/handle. After rinsing is done, empty the waste tank as described above.

Check and Storage

- Check the parts of this toilet for any signs of looseness, wear, or damage between uses. Tighten, repair, or replace any problematic parts before further use.
- To store this toilet, empty both the waste tank and the water storage completely. Clean all tanks and surfaces as described above. Store in a well-ventilated place away from direct sunlight and rain.