

USER MANUAL

3D Massage Chair (With Air Pressure)

Many thanks for purchasing our massage items. Before using this machine, please read this manual carefully to make sure fully understand the exact method of using, especially the points for attention. Please keep this manual carefully for your reference!

Functions

◆ NECK PART:

1. Four kneading heads working by reverse&forward.
2. It's adjustable height for met different height people.

◆ BACK PART:

- 1.Kneading & Tapping together work
- 2.Up & down massage full back
- 3.Can be fixed massage for Higer back, Lower back for your choice

◆ HIP PART:

Three vibration intensities supply comfortable vibration on the butt& thigh

◆ LEG PART:


Five angles to adjust different position for relaxing you leg.

- ◆ It's a collapsibel chair for package,no need assemble,Once you get it can used directlly.
- ◆ 15 mins time limited.

◆ With 5 pcs air bags supplies pressure for the back part & hip part

◆ Specially pressure in Mid-part of hip

Controller

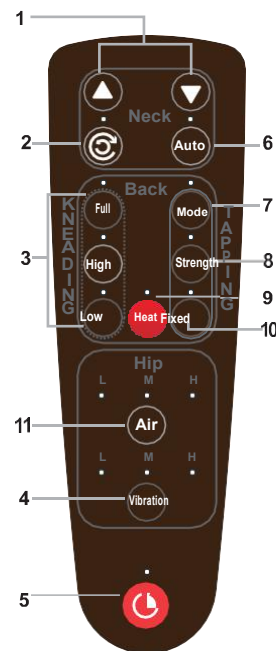

1. Neck Massage Adjustment

Adjust your desired neck massage level.

2.3D Neck Massage Programme

Change the direction of kneading heads.

3. Back Zone Selection

Select your desied back massage area (FULL,UPPER,LOWER).

4.Vibration Massage

To relax your buttocks and thighs.

5.Power ON/OFF

Switch ON or OFF the machine.

6. Auto Programme

Press this button by starting working automatically.

7. Tapping, Kneading & Rolling Massage

Provides 3 different massage techniques to relief sore muscles.

8. Massage Intensity

Select your desired massage intensity levels.

9.Heat Therapy

Improve metabolism rate.

10. Spot Massage

Aunique feature which allows you to manually focus on desired area of your body

11. Air Pressure Massage System

Total 9 built-in air bags from back to waist, buttocks and thighs possition to provide relaxing massage.

How to operate

1. Open the packing ,place the massage chair in a flat place& open
2. Connect the DC with the adaptor
3. Connect the adaptor with the power
4. Press the button on controller to start working.

Notes

Before Use	Please read carefully the instruction before use.
Practice	- Working 20 minutes every day at the beginning stage and more time step by step; do not exceed 20 minutes for safety. - Could place one folded towel or blanket on back in case of too strong massage pressure.
Longer Life	Please let the machine be rest for 20 minutes after 40 minutes operation so to ensure of longer life of machine.
After Use	Please turn-off power and then pull out the plug after usage; please do pull the plug from the socket and do not pull the wire.
Place to put	- Please do not operate the massager where with much dust, oil crust, fog. - Please do not put the massager near to electric appliance. - Please do not put the massager where under direct sunshine - Please do not use the massager where with much humidity (like hot spring, swimming pool, bathroom etc). - Please do not use the massager outdoor. - Please do not operate the massager where with fog and lack of oxygen.

Notes

Suitable for the crowd	- Those who feel uncomfortable, been operation or be injured on back, neck, please consult with doctor before using. - Do not use within one hour after lunch/dinner. - The following people please use the massager after consulting with doctor or specialist : a) with cancer b) heart disease - Pay much attention when the massager is running, especially for patients or those without ability of operation. - You are suggested to wear baggy clothes when using the massager; Do not use it nakedly which might cause hurt to your skin; Too thick or thin clothes or those with collar, belt or other hard objects should also be avoided because they might affect massage effects. - Prohibit using for kids and those pregnant women
------------------------	--

FAQ

Phenomenon	How to check	How to handle
Sounds like squeak, ping, click,	Sound is caused by motor; normal; no need to check	N/A
Noise during massager	Has it been working for a long time or continue?	Please turn off power and be rest over one hour
Feeling hot on the surface	Has the working time exceed time limit?	Please operate within time limit
No working or stop suddenly during running	1. Power off 2. Temperature protection device	1. insert the plug and turn on the power 2. turn off power and be rest over one hour
Please consult the distributor if no use after the above handling		

DANGER

To reduce the risk of electric shock

1. Always unplug the appliance from the electrical outlet immediately after use and before cleaning.
2. Never use pins or other metallic fasteners with this appliance.
3. Keep dry-Do not operate in a water or moist condition.

Specifications

Product Name	3D Massage Chair
Voltage	12V $\equiv$ 5A
Rated Power	60 W
Time Limited	15 min
Product Size	130 X 66.5 X 92cm
N.W.	22 KG
G.W.	23 KG

CAUTION

- ◆ TO REDUCE THE RISK OF ELECTRIC SHOCK,DO NOT REMOVE COVER,THERE ARE NO SERVICEABLE PARTS INSIDE.
- ◆ TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS ITEM TO RAIN OR MOISTURE.




The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated "dangerous voltage" within the unit's enclosure that may be of sufficient magnitude to constitute a risk of electric shock.



The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit


ASSEMBLY INSTRUCTIONS

►STEP1

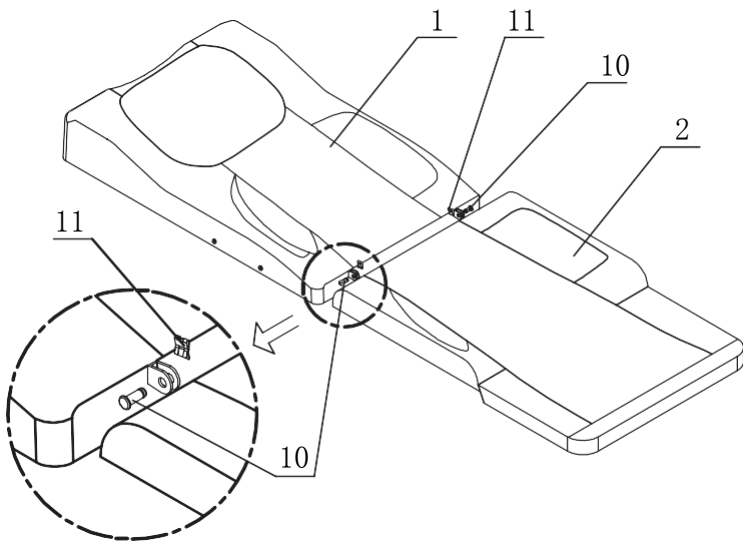

No.	Description	Qty
1	Chair Back	1 pc
2	Chair Cushion	1 pc
10	Coupling shaft	2 pcs
11	Spring piece	2 pcs

►STEP2

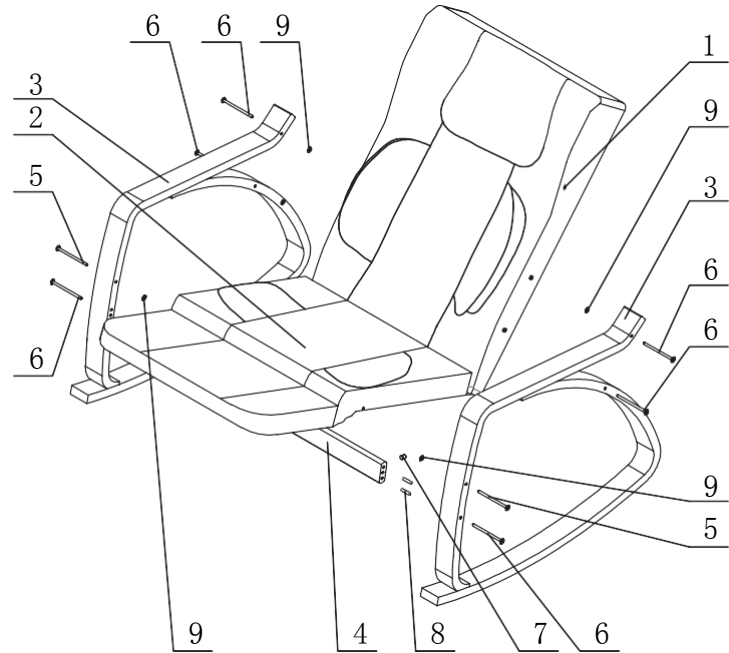

No.	Description	Qty
1	Chair Back	1 pc
2	Chair Cushion	1 pc
3	Armrest	2 pcs
4	Connecting boards	1 pc
5	Screw bolt M6 X 105mm	2 pcs
6	Screw bolt M6 X 95mm	6 pcs
7	Screw nut M6	2 pcs
8	Wooden pin	4 pcs
9	Gasket	6 pcs