


MASSAGE CHAIR

USER MANUAL

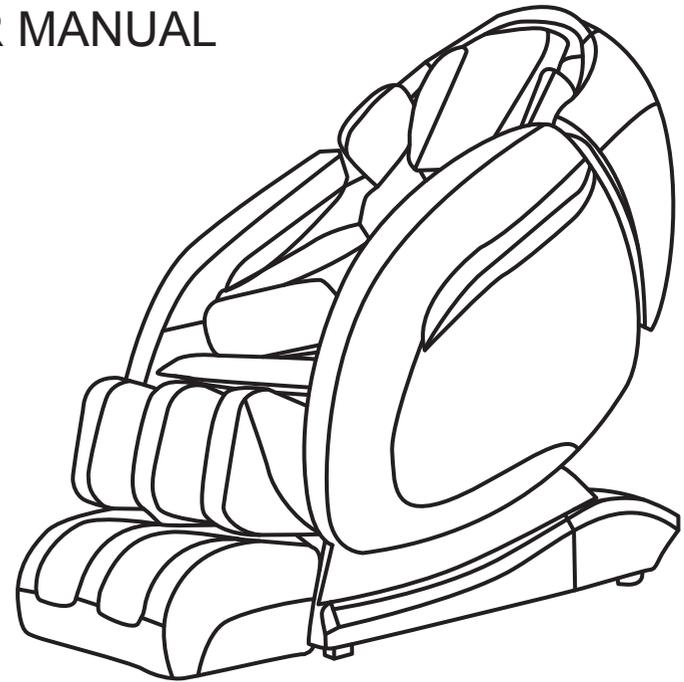

Operating manual

Before using this product, be sure to read this manual carefully to avoid malfunction and danger caused by improper use.

Thank you for purchasing this product. Please read this manual carefully before use so as to understand the correct using method of this product. Please notice the important safety warnings while reading and keep this manual well for your reference at any time.

The company reserves the right to make changes and explain to the design of the product, please forgive that it also subjects to change without notice, and the color of the product should be subject to the real object.

Contents

Safety and maintenance issues.....	1
Function.....	4
Part name and function.....	5
Massage chair installation.....	6
The using method of the massage chair.....	8
Common phenomena and malfunction elimination.....	12
Product's technical parameters.....	13

Product specification

Product name	Zero gravity multifunctional massage chair
Rated voltage	110V/60Hz
Rated power	170W
Safety classification:	Class I
Rated time	20min

Common phenomenon processing

- Plug in the power, and the massage chair cannot work. Check whether the power plug is well plugged in and connected to power:
- If the power connection is intact, but the massage chair can't work, and the switch light won't respond, please replace the fuse and refer to page 5 for replacement method;
- If it suddenly stopped while massaging, it may be that the timing time has expired, or the shutdown button is pressed accidentally on the control panel;
- If the motor sound emits during operation, it would be a normal running sound;
- If the machine stops running, it means the working time is up, please repress the button to operate;
- If the power button light on the hand controller flashes, but it cannot be turned on. Please check whether the manual controller button is not popped up, or the line core of the hand controller connector falls off;
- If the massager doesn't work normally, try to turn off the main power and then turn it on again.

Prompt

- If you encounter any problems when using this product, please refer to the operating instructions and solutions in this manual. If the product is still in a fault state or it cannot be used, please turn off the power and pull out the plug, and contact the manufacturer or dealer;
- Do not disassemble or repair this product yourself. Damage to the product due to disassembly or repair by itself will void the product warranty.

Warnings!

1. Do not use the massage chair at voltages that do not meet the standards;
2. Always use a socket that matches the plug;
3. Please use this product according to this operation manual;
4. Do not use this product outdoors;
5. This appliance is not intended to be used by persons (including children) who are physically or mentally handicapped or lacking experience and common sense, unless they are supervised and directed by those who are responsible for their safety regarding how to use this massage chair;
6. Children must be supervised to ensure that they won't use it for playing;
7. If the power cord is damaged, it must be replaced by professionals of the manufacturer, its maintenance department or similar departments in order to avoid danger.
8. This unit has a heating surface, those who are not sensitive to heat must be cautious while using.

Operating environment

1. Please do not use it in the shower room or other damp environments;
2. Please do not use it when the environment temperature changes sharply.

Those who are not suitable to use this product

1. Patients with osteoporosis;
2. Patients who have heart diseases and those who are wearing electronic medical appliance such as cardiac pacemaker;
3. Women who are pregnant or in menstrual period;
4. Those who have physical wounds or have physical appearance illness;
5. Those who are allergic to heat;
6. Those who should take rests as the doctor instructed or those who feel not well.

Safety precautions

1. Please check whether the voltage meets the specification requirements of this product;
2. Please do not unplug the power plug with wet hands;
3. Do not use wet cloth to clean the power switch, plug or the like which may carry electricity;
4. Please stop using it immediately if you detect any abnormality during using this product.

Cut the power and stop massaging

- During use, press the on-off  key to close all massage functions and enable the massage chair to reset automatically;
- When the timing time of massage is up, it will automatically turn off all massage functions and the backrest, shanks frame and seat will automatically reset and keep the original shape.

Key function of the hand controller

Manual mode

It has a manual massage technique button, press the button to exit the automatic massage technique, and the massage technique switches between seven manual massage techniques;

Automatic mode

Switch the button in automatic massage technique mode, press this button to switch among six automatic modes;

Automatic air bag

A button for switching to the automatic mode and off mode of the whole body air bag;

Heating

Waist heating function on/off button;

Roller

Foot roller function on/off button;

Shift up

The massaging part of the robot moves upwards. When release the button under the manual mode, it can massage at the current position;

The massaging part of the robot moves downwards. When release the button under the manual mode, it can massage at the current position:

Product maintenance and repair

1. This product can only be repaired by the service center designated by the manufacturer, and the user cannot disassemble and repair it by himself;
2. Please turn off the main power switch after the action;
3. If you do not use it for a long time, please wind up the power cord of this product and store it in a dry and dust-free environment;
4. Do not store it adjacent to high temperature or open flame, do not expose to the sun for a long time;
5. Please use a dry cloth to clean the massage chair other than thinner, benzene or alcohol;
6. Please use it intermittently, should not let this product run continuously for a long time.

Function

1. L-shaped guide rail long massage stroke, covering the head and neck - shoulder - spine - waist-hip;
2. A group of freely walking, four-wheel drive, silent design massage manipulator on the backrest;
3. Seven kinds of massage techniques including finger massage, kneading, flapping, tapping, kneading & tapping synchronization, kneading & flapping synchronization and squeezing;
4. Six sets of full-automatic massage modes including soothing, vitality, whole-body, sleep, neck and waist, experience
5. Pressure massage function for shoulders, arms, buttocks, legs, feet;
6. Back heating function;
7. The semi-automatic telescopic function of the shank;
8. Finger pressure of the foot roller cooperates with the foot air bag to wrap the foot and extrude;
9. Connect the bluetooth to play music, and the music synchronization function will allow the intensity of massage to move with the music frequency;

Key function of the hand controller



Turn on/off the massage chair. The massage chair will automatically tilt at an angle when turn it on; When turn off the unit, the TFT will display that it is in the off state, and the massage chair will be automatically reset. After the reset is completed, the screen will be off;



The massage chair will be automatically adjusted to zero gravity state when it is turned on; When it is turned off, the posture of the massage chair will be restored to the initial state of starting;



When the button is pressed, the shanks will rise and the posture adjustment will stop after it is released;



When the button is pressed, the shanks will fall and the posture adjustment will stop after it is released;



When the button is pressed, the backrest will rise and the posture adjustment will stop after it is released;




When the button is pressed, the backrest will fall and the posture adjustment will stop after it is released;

Button layout of the hand controller




Turn on the power of the whole machine

Refer to page 5 to find the power cord jack, plug in the power cord, turn the power switch to 1 position, and the switch indicator light will be on, indicating that the power is on.

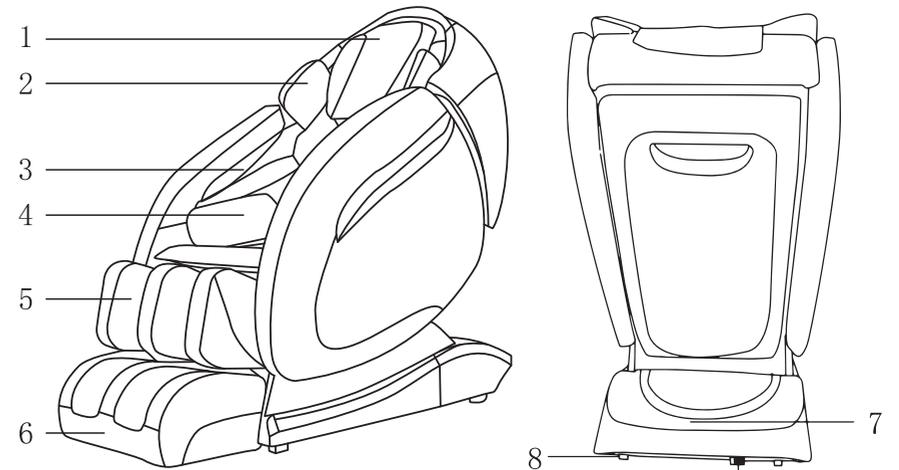


- 1. Loudspeaker
- 2. Shoulder air bag
- 3. Arm air bag
- 4. Hip-fixing air bag
- 5. Leg air bag
- 6. Foot roller
- 7. Electrical cabinet
- 8. Moving caster

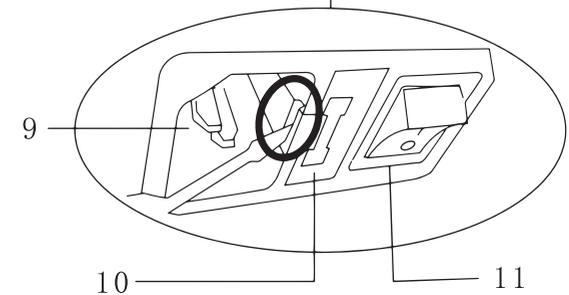

Notes: When replacing the fuse, use a screwdriver to open the circle indication position and remove the fuse box.

Unpacking the massage chair

Lift the outer box up, remove the fixing post, and take out the instruction manuals and the power cord and other accessories.

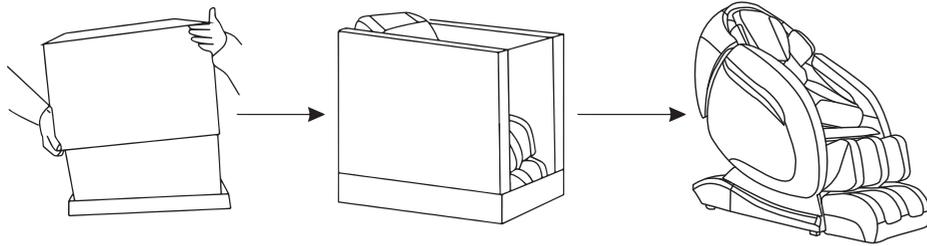


Move the massage chair

As shown in the image on the right, tilt the chair backwards at about 45°. Use the casters to carefully push forward or backward to the appropriate position and place it down.

Note:

- * When someone is lying on the massage chair, please don't move it.
- * When moving the massage chair, the floor may be damaged. It is recommended to put a cushion under the product.
- * In order to avoid falling, please hold the massage chair firmly when you move it.

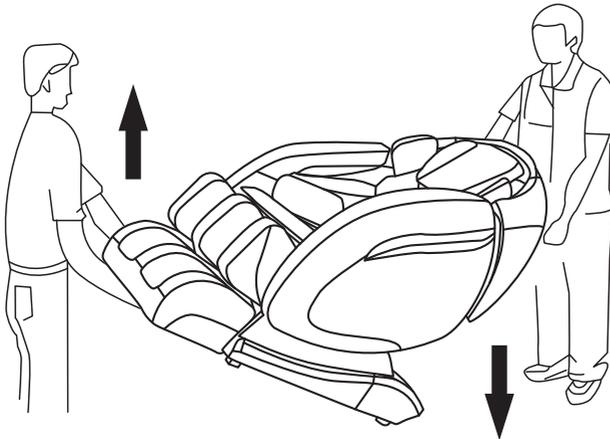

Place the massage chair

When placing the massage chair, be sure to leave enough space around the chair so that the legs and backrest of the chair can rise and fall freely. Arrange the space as suggested in the right picture.

