

尺寸：120*80mm 双面印

材质：

120mm

LED Skin Care Mask

Product Description


This Light Therapy Mask includes eight modes, each mode is effective for skin improvement.

The effects on the skin are as follows

- 1.Red light: Collagen regeneration, stimulate metabolism, enhance elasticity, resist wrinkles and fine lines;
- 2.NIR infrared: Promotes blood circulation and tissue repair, enhance skin radiance and helps soothe sensitive and damaged skin;
- 3.Blue light: Anti-acne and cleansing,shrinking pores,repairing sensitive skin, relieving various discomforts of skin,and insist on using blue light can improve the balance of skin water and oil, improve seasonal affective disorder;
- 4.Orange light: Soothe the skin, and repair pigmentation issues
- 5.Purple light: Control oil secretion,brighten and fade dark spots,stimulate collagen production for firmer and more radiant skin;
- 6.Indigo light: Soothe the skin,enhance cell energy gradually,and has a good metabolism facilitating effect;
- 7.Green light: Normalizes sebaceous glands to reduce blackhead formation and moderate oil production, suppresses excessive melanin production,reduces pigmentation,reduces intraocular pressure;
- 8.Yellow light: Improve rough skin and wrinkles,alleviate redness and sun-burn,soothe the skin and gently promote collagen production;

1

Product Specification

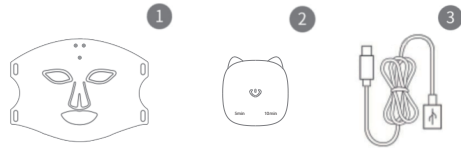

① Silicone Flexible Mask ② Controller ③ USB Charging

④ Goggles ⑤ Fixed Strap

Product Parameters

Product model: RB-008.J

Charging voltage: 5V

Rated voltage: 3.7V

Charging current: 600MA

Rated power: 4W

2

Mask Operation Instructions

1. Assemble eye goggles. The face of the mask is facing down, the eye goggle with "L" assemble on the left eye, with the "L" position facing outward. The eye goggle with "R" assemble on the right eye, with the "R" position facing outward.
2. Assemble fixed straps. Fixed straps do not distinguish between left and right. With the flannel facing inward, place two straps with small velcro strips through the two holes on the left and right of the mask. The velcro is attached to the tape and the position can be adjusted to the left or right according to your needs.
- After wearing the mask, two large velcro pieces can be attached together.
3. Connect the mask to the remote control.

Controller Operating Instructions



①Power key: Long press the power key to turn on/off. One click to adjust the light modes. Double-click to switch the timing.

②Power key/Time display: When it' s power on, the default time is 5 minutes (the left white LED indicator is on and the buzzer beeps). Double-click to switchon the right is on, and the buzzer beeps twice), if switch modes, recalculate the timing for 10 minutes (the white indicator on the leftis off,the white indicator the time during the timing process (cycle in turn)

3

③Power key/Mode display: Long press the power key, it's on, boot default red light, buzzer beeps, the mask enters into red light mode, automatically shut down after 5 minutes.

After first click, buzzer beeps,the mask enters into infrared + red light mode, automatically shut down after 5 minutes.

After second click, buzzer beeps, the mask enters into blue light mode, automatically shut down after 5 minutes.

After third click, buzzer beeps, the mask enters into orange light mode, automatically shut down after 5 minutes.

After forth click, buzzer beeps, the mask enters into purple light mode, automatically shut down after 5 minutes.

After fifth click, buzzer beeps, the mask enters into indigo light mode, automatically shut down after 5 minutes.

After sixth click, buzzer beeps, the mask enters into green light mode, automatically shut down after 5 minutes.


After seventh click, buzzer beeps, the mask enters into yellow light mode, automatically shut down after 5 minutes.(Cycle in turn)

④Power key/power level: When the remote control is working, the corresponding mask LED light is steady on, and when the battery voltage is lower than 3.4V, the red LED indicator flashes. When battery voltage is 3.2V, the Buzzer beeps three times, then the machine shuts down.

When charging, the red LED indicator is breathing, the green LED indicator is on. After fully charged, and the charging power indicator is off.

Instruction for use

4

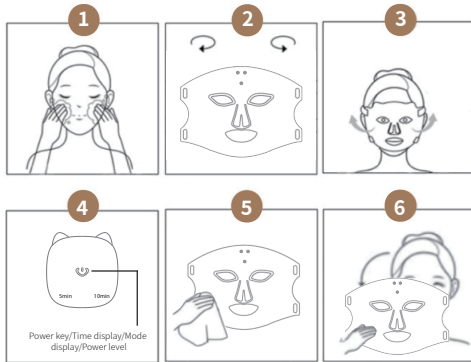

Notice:

1. After using the care, use a soft towel or wet towel to wipe the silicone part, please dry it in ventilation and store it in a dust-proof bag
2. Do not directly contact the water source.
3. Recommended 3-5 times a week, 4 consecutive weeks as a cycle.
4. For the first use, please fully charge the device, and do not use it while charging.

The following people should not use this product	1. People who are pregnant and breastfeeding; 2. People who are photosensitive or who are taking any drugs that may cause photosensitivity; 3. People with eye genetic diseases or skin diseases (such as lupus erythematosus, photosensitive eczema or albinism, etc.)
Charging instructions	1. This product contains a lithium-ion battery, and the charging of this product does not support power-on work 2. When the control brightness indicator flashes, the device is in the charging state. When the three control brightness indicators are all on and stop flashing, the charging is complete.

5

Light Sensitivity Test

Before using this product for the first time, please conduct a light sensitivity test:

First step	Put the device on your arm and turn on the switch to irradiate for 10 minutes. If the skin becomes red or itchy. Symptoms, indicating that you may have light-sensitive skin. It is recommended that you stop using it, or before using it again, Consult a dermatologist
Second step	Before use, please cleanse your face and keep it dry (without makeup) to prevent cosmetics that may contain Some light-sensitive ingredients affect skin care efficacy
Third step	After use, the skin may appear reddened, which is a normal reaction of the skin, and will occur within 1-2 days. It will return to normal within hours, if you feel burning or other discomfort during the process, please stop using it immediately

How to use and steps

First step	Clean face thoroughly
Second step	Connect the straps to the mask, place the mask on the face, adjust and fasten the straps to secure the mask On the face, long press the controller button to turn on the device
Third step	According to the different light energy, this device sets eight working modes respectively, which can be adjusted according to the needs. After starting up, the device defaults to red light mode - RL, can change light modes by pressing the light mode button. Each light mode can adjust time 5mins/10mins by press power button.
Fourth step	The device will automatically turn off within the set time after it is turned on. After the device is turned off, take off the mask and end the irradiation, Follow your daily skin care regimen.

6

Equipment maintenance and maintenance

Deposit	Please store this product in a normal temperature and dry environment, avoiding dust, direct sunlight or humid environment.
Clean	Before and after use, the mask can be wiped with a water-based rag or a damp cleansing cloth. When cleaning, make sure that the device Turn off and disconnect the power supply
Recycling	This product contains a lithium battery, and please send the device to a designated local. The collection point for recycling processing is used for recycling processing. It is recommended to separate them from other household wastes to contribute to environmental protection. Do your part.

Warranty Service

Under normal use, if there is any quality problem with the product itself, this product provides a one-year warranty for the main unit and a half-year warranty for accessories. Man-made damage and non-quality problems are not covered by the warranty.

If you have any questions or issues regarding the product, please feel free to contact us via email at supportsir@163.com, and we will assist you promptly.

7

80mm