

Massage Chair

An instruction manual

Pictures are for reference only

Living Art

Thank you for purchasing our massager products.

This massage chair adopts ergonomic design to massage neck, arm, back, waist, buttocks, legs, feet and other parts, with various massage methods, warm compress, air bag extrusion massage and other functions. It is a luxury high-end massage chair.

Before use, please read the manual carefully and pay special attention to safety precautions to ensure the correct operation and use of the product.

Take good care of the manual for future reference.

Catalog

I. Product Characteristics

II. Usage

III. Maintenance and Maintenance

IV. Matters needing attention

V. Safety Protection

VI. Warm Tips

VII. Technical parameter

VIII. Circuit schematic diagram

IX. Names and contents of harmful substances in products

X. Product warranty instructions

XI. After-sales Service Guarantee Card

XII. Analysis and Removal of Common Faults

XII. Analysis and Removal of Common Faults

Failure phenomenon	Possible cause	Processing Method
Sound (working sound) *inflatable sound of airbags *bleed sound of airbags *changing-over sound on solenoid valve *air pump working sound	The sound produced by the product at work Sound of mechanical operation	This sound is normal No need to deal with it
The whole machine is not working at all	*The power plug is not inserted into the socket *Socket switch not open *The massage button is not pressed	Check the power plug and socket, ensure that the plug is in good contact with the socket, press the massage button.
The massage stopped suddenly in the process	This product has just completed 15 minutes of rated working time and automatically switched to standby mode.	If you want to massage again, please press the key again and open the massage function.
	Overheating protection for motors exceeding rated service time	Turn off , Cut the power, use the unit again after the motor is cool down

*The company reserves the right to improve products and instructions.

If alteration is made without notice, the appearance of the product shall be in accordance with material object.

*This volume has been carefully checked, but it does not exclude the possibility of misprints or misunderstandings in the content.

We are sorry to inform you that we reserve the right of interpretation.

♻️ Improper disposal of waste causes environmental pollution.

Certificate of Conformity

Name: _____ Model No.: _____

Inspection No.: _____ Production Date: _____

XI. After-sales Service Guarantee Card

Warranty Card

With this card warranty for one year,
users will keep it properly.

User Union

Buyer ID: _____
Order No.: _____
Product Name: _____
Model No.: _____
Store: _____
Date: _____ Price: _____
User name: _____
Telephone: _____
Address: _____

Warranty Card

With this card warranty for one year,
users will keep it properly.

Company Stubs

Buyer ID: _____
Order No.: _____
Product Name: _____
Model No.: _____
Store: _____
Date: _____ Price: _____
User name: _____
Telephone: _____
Address: _____

Maintenance record:

Date	Fault Condition	Address	Result

I. Product Characteristics

- The whole machine has beautiful lines, bright colors, humanized design and green environmental protection materials.
- This product is designed with 2 groups of movements and 4 massage heads, which can knead and massage the human neck, back, waist and buttocks.
- The whole machine is equipped with about 20 airbags in feet, legs, arms and shoulders, which can freely display the massage function of air bags.
- It can manually adjust the massage movement of neck, back, waist and buttock. Various massage techniques can be changed freely, and three-level adjustment can be realized.
- Air bag massage one key switch, can manually adjust the shoulder, arm, leg air bag massage, massage strength to achieve three gear adjustment.
- The recline angle of the back can be adjusted according to the needs, and the free adjustment of 115 ° to 155 ° can be realized. When the angle of the back is adjusted to the minimum, the zero gravity function can be realized.
- One button heating of waist and legs, constant temperature 35 ° to 45 ° C.
- The massage time can be adjusted according to the need, with 5 minutes plus or minus echelon, you can set 5 minutes, 10 minutes, 15 minutes, 20 minutes, 25 minutes, 30 minutes to adjust a massage journey.
- Special intelligent body detection, according to the shoulder height set massage itinerary.
- Special SL type extended guide rail.
- Special sole roller.
- Special hifi Bluetooth audio, massage at the same time, enjoy the fun of music.

II. Usage

One-in-one installation-free

No tedious installation operation is needed to open the box.

1. Installation Method of Connection Line



Connect the type plug of the power cord to the plug of the electrical box (As shown in the figure) and the other end to the power socket.

- People stand on the top of the back of the massage chair and press down to lift their legs off the ground.
- Tilt the chair back 40 ° to facilitate free movement of the massage chair indoors.
- Do not press down too much to avoid damage to the fuselage due to the complete collapse of the fuselage.
- Please do not move the massage chair when someone is sitting on it.
- It is suggested that the massage cushion may be damaged when the product is laid on the floor.

IX. Names and contents of harmful substances in products

Part name	Hazardous substances					
	(Pb)	(Hg)	(Cd)	(Cr (VI))	(PBB)	(PBDE)
PCB board	X	○	○	○	○	○
Leather cover	○	○	○	○	○	○
Connecting line	○	○	○	○	○	○
Electric machinery	X	○	○	○	○	○
Sponge	○	○	○	○	○	○
Packing materials	○	○	○	○	○	○
Metal frame	X	○	○	○	○	○
Plastic parts	X	○	○	○	○	○

This form is prepared in accordance with SJ / T 11364
○Indicates that the content of the hazardous substance in all mean materials of the part is below the limit requirements specified in GB / T 26572.
×Indicates that the content of the hazardous substance in at least one average material of the component exceeds the limit requirements specified in GB / T 26572.

The "harmful substances" shown in this table will not cause any harm to human body and environment under normal use of the product. The "hazardous substances" shown in the table and their existing components provide consumers and recycling practitioners with information about the existence of relevant substances, which is helpful for proper disposal when the products are discarded.

Suggestive notes to regulations on electronic recycling and treatment of waste electrical appliances



The company advocates environmental protection and low-carbon consumption experience. It is recommended that users who own the product should abide by the national laws and regulations on recycling of waste electrical and electronic products when they no longer need the product, and hand it over to local manufacturers with national recognized recycling qualification for recycling treatment, so as to reduce resource waste and protect the earth's natural environment.


X. Product warranty instructions

1. In order to protect the interests of consumers, our products are strictly implemented in accordance with the three-package regulations of the state. Consumers can enjoy a one-year warranty service. If the product is damaged due to product quality problems during the warranty period, the company can repair and replace spare parts free of charge.
2. During the warranty period, the company will charge the cost of materials as appropriate if the service is required due to failure caused by the following circumstances:
 - Used improperly or repaired by oneself;
 - Users who do not comply with the prescribed power supply voltage;
 - Self-moving, transporting or refitting after installation;
 - Those whose appearance is naturally dirty due to use;
 - Damage caused by natural and man-made disasters and insect pests.
3. In case of the following circumstances, the company charges according to the price:
 - Beyond the warranty period;
 - Failure to show the warranty card;


VII. Technical parameter

Product Name:	Massage chair
Rated voltage:	110V~
Rated frequency:	50Hz
Rated power:	110W
Rated working time:	15mins.
Security level:	Class I
Implementation criteria:	GB4706.1-2005 GB4706.10-2008

VIII. Circuit schematic diagram

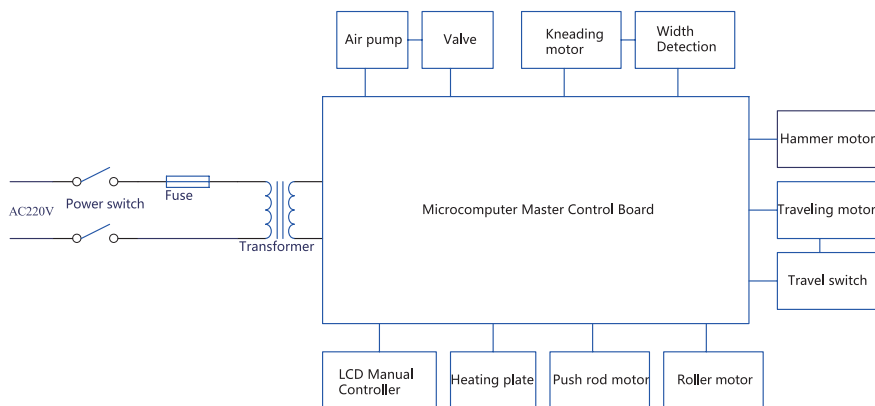

2. Manual operation instructions

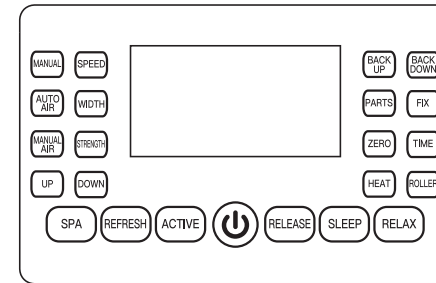

: on or off.

Manual: Adjust different massage ways

Speed: Adjust the massage speed (only available in some manual modes)

Auto air: open the whole body air bag with one button.

Width: adjust the massage head "wide medium narrow" (only available in some manual modes)

Manual air: manually open the shoulder arm and leg airbags.

Strength: air pressure and force adjustment.

Up: long press the movement to walk upward (when manual fixed-point massage).

Down: long press the movement to walk downward (when manual fixed-point massage).

Back up: press the back for a long time.

Back down: press the back down for a long time.

Parts: Press it to adjust massage in different parts of body (only available in manual mode)

Fix: press the fixed point key to massage (only available in manual mode).

Zero: It is a comfortable reclining degree between the body and the legs


Time: timing setting.

Heat: heating the waist and legs.

Roller: foot roller massage.

Automatic massage mode: SPA, REFRESH, ACTIVE, RELEASE, SLEEP, RELAX

3. Shortcut Key Operation Instructions

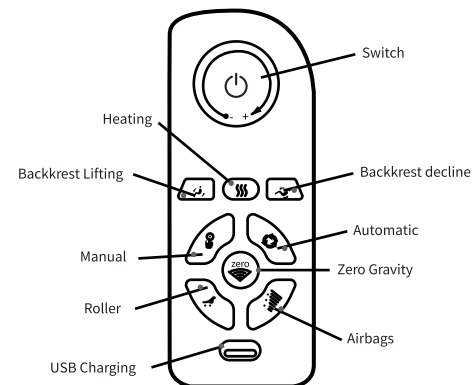

Switch: Turn on or off. Rotate the switch to slow down or accelerate the speed.

Backrest lift: Press and hold the backrest for a long time to lift it up.

Heating: Turn on or off the heating of the waist and calves.

Backrest decline: Long press the backrest to lower

Manual: H1, H2, H3 different massage.

Automatic: Different massage modes with A1, A2, A3.

Zero Gravity: The massage chair is in the zero gravity position.

Roller: Turn on or off the foot rollers.

Airbags: Turn on or off the full body airbag.

USB Charging: For phone charging.

III. Maintenance and Maintenance

1. Please use neutral detergent and soft cloth for daily cleaning, and do not use corrosive detergent.
2. If the power cord is damaged, it must be replaced by the full-time personnel of the maintenance shop designated by the manufacturer to avoid danger.
3. When not using this product or cleaning, please unplug the power supply to avoid danger.
4. Do not use any sharp objects to cut the product surface.
5. Do not put the product in a humid environment, but in a dry and ventilated place.
6. Do not press or press the product hard.

IV. Matters needing attention

In order to ensure the safety and correct use of the massage chair, please refer to the operation manual in detail before using this product.

To reduce the risk of fire, burns, electric shock, or injury to others, follow these instructions:

- 1) Power supply other than AC 220V shall not be used.
- 2) When no one is using the machine, please be sure to remove the power plug of the massage chair from the socket. When not using the machine or before removing the device / parts, please also remove the power plug from the socket in time.
- 3) In order to avoid fire, electric shock or injury caused by overheating, please do not put the product under pillows, blankets and other covers.
- 4) If the product is used for children or the disabled, adult guardians are required to supervise and guide it.
- 5) Please be sure to use the product according to the usage described in the user manual. Do not use parts that are not recommended by the manufacturer.
- 6) Do not use this product in case of any of the following conditions: the power cord or plug is damaged, the operation is not normal, if it falls to the ground or other damages, or falls into water, please send the product back to the service center for inspection and maintenance.
- 7) Do not drag the power cord to move the product or use the power cord as a handle.
- 8) Keep a safe distance between the power cord and the heating surface.
- 9) Do not operate the product when the air outlet is blocked. Please keep the air outlet unobstructed to avoid being blocked by cotton chips, hair and other things.
- 10) Do not insert or drop objects into any openings of the product.
- 11) Do not use this product outdoors.
- 12) please do not operate the product in the place of spray product or oxygen filled place.
- 13) Cut off the power supply, please turn off all operation controls, and then unplug the power plug from the socket.
- 14) Do not cover the massage chair so as not to cause malfunction or fire due to excessive temperature.

V. Safety Protection

- Do not use the massage chair in excessively humid or dusty environment to avoid electric shock or mechanical failure.
- Do not use the massage chair when the room temperature exceeds 40 °C.
- Do not place the massage chair near the electric heater, stove or direct sunlight.
- The massage chair is designed for personal use and is not suitable for any commercial use.

Careful

- Massage chair is only for indoor use.
- Use massage chair at low room temperature, do not suddenly adjust the temperature. It is recommended to gradually adjust the temperature to the normal range before use.
- If the massage chair is moved from a low temperature storage place to a warm environment for use, it may affect the normal operation of the machine due to the accumulation of moisture in the parts, or even cause parts failure. It is recommended that the product be placed at normal room temperature for 1 hour before use.

VI. Warm Tips

- The power switch must be turned on before other functions can be selected.
- Children should use this massage chair under the supervision of adults to avoid equipment damage and children's injury due to improper operation.
- Please do not let children climb this massage chair to avoid children falling.
- Use the massage chair to adjust the strength of the chair massage, because the massage force is too small to play its due role in stimulating, too large prone to fatigue, and damage the skin, massage force and times from less to more, from light to heavy.
- Use the massage chair to relax the body and mind. In addition to thinking at home and abroad, you should be calm and relaxed.
- Massage chairs are not suitable for patients with heart disease, hypertension, osteoporosis, local skin damage, ulcer, bleeding, tuberculosis, tumor, people who are too hungry, full, drunk or tired.
- In order to avoid the reduction of massage effect or injury, please do not use the massage chair together with other effective equipment such as electric blanket.
- Do not use the massage chair within 1 hour after eating to avoid adverse effects.
- Do not use the massage chair when your hands or other parts of your body are wet. Do not plug or unplug the power plug.
- Do not insert any articles or inflammables into the massage chair. This may cause electric shock, malfunction or fire.
- If you feel any discomfort during massage, please stop immediately.
- If the product is not used for self diagnosis or treatment and needs diagnosis or treatment, it shall be handled by professional medical personnel.
- Do not use this product if you are in medical condition and are pregnant or in physiological period.
- Use the massage chair to master the massage time, 15 minutes for each massage, and no more than half an hour for each person (time) (when it reaches half an hour, let the massage chair rest for ten minutes before use). If it is used continuously for a long time, the overheat protection device may be started and the power will be cut off automatically.
- The surface of the appliance is hot. People who are not sensitive to the heat must pay attention to it when using it.
- Do not directly massage elbow joint, knee, chest, abdomen, etc.

Not suitable for massage chair:

In case of any of the following conditions, please turn off the switch immediately, cut off the power supply and contact our company.

- accidentally spilling liquid on the massage chair may cause electric shock.
- the cloth cover of massage chair is damaged and the internal parts are exposed.
- in case of abnormal or severe pain during massage, please stop massage immediately and seek medical consultation.
- please stop using immediately if any mechanical failure or abnormal condition is detected during use.
- in case of current interruption in use, attention shall be paid to avoid user injury due to sudden power recovery.
- in case of lightning.