READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

⚠ Warning

As with all electrical appliances, there is an inherent danger of overheating, fire, or personal injury if the product is not used correctly. Please read all instructions carefully before use.

⚠ Do Not Use This Heating Pad if You Are

*This appliance is not to be used by children.

Close supervision is necessary when any appliance is used near children. Children should not play with the appliance.

*An incapacitated or immobile person.

*Diabetic or insensitive to heat.

*A paraplegic, quadriplegic or amputee person.

*Having poor blood circulation.

*Over cuts, skin lesions, sores or areas of damaged skin.

*Anyone who cannot clearly understand the instructions and/or operate the controls.

⚠ Caution

1. Do not use while sleeping.
2. Do not use on an infant or on an animal.
3. Do not use pins or other metallic fasteners.
4. Do not sit on, or crush pad - avoid sharp folds.
5. Do not use on insensitive skin areas.
6. Always unplug when not in use. Never leave unattended, especially with children present.
7. Not for use by invalids, unconscious persons, or those with poor circulation without supervision.
8. Burns may occur at any heat setting - check skin frequently during use.
9. Never pull this pad by the supply cord. Do not use the cord as a handle, do not pinch cord.
10. Loop cord loosely when storing. Tight wrapping may damage cord and internal parts. Allow to cool completely before storing. Do not place items on top that may cause creasing. Store in original packaging in a cool, dark, dry place.

11. Do not dry clean this heating pad. Cleaning solvents can damage the heating element's insulation.

⚠ Caution


12. This Heated pad is for household use only. Not for medical/hospital use.
13. Always turn off and unplug when not in use.
14. Do not use with other heating products.
15. Never open up the controller case. Doing so could prove potentially hazardous and will void YOUR WARRANTY!
16. Keep all animals away from this heated product. They could damage the wiring, leading to a potentially harmful situation.
17. Do not use this heated product with analgesic creams, liniments, salve or ointment preparations. Skin burns could result.
18. If the heating pad malfunctions, stop use immediately and contact the manufacturer for repair.
19. In the event of abnormal operation, immediately disconnect the power and remove the device.
20. Do not plug or unplug the USB charging cable with wet hands, as this may cause electric shock or malfunction.
21. Do not use in bathrooms or high-humidity areas. Risk of severe electric shock and product damage.
22. When unplugging the DC charging cable from the device, hold the plug instead of pulling the wire. Pulling the wire may cause electric shock or damage the product.
23. The product may become warm during extended use—this is normal.
24. Save these instructions.

Operating Instructions for the Switch

After connecting the mobile power supply, press and hold the switch for two seconds to turn it on/off. Adjust three heating settings (high 149°F, medium 131°F, low 113°F) with just a simple press of the button.

Product Introduction

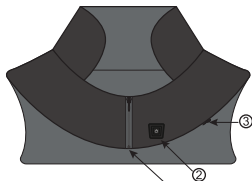- ① KCC Zipper
- ② Silicone Controller
- ③ KCC Pocket Zipper

Weighted Design
Filled with micro-glass beads
to relax your body

Capacity:

15000mAh, 3.6V, 54Wh

USB-C input:

5V=3A, 9V=2A, 12V=1.5A

USB-C output:

5V=3A, 9V=2.22A, 12V=1.67A

USB-A output:


5V=3A, 9V=2A, 12V=1.67A

DC output:

9V~12.6V=3A

4 power level indicator lights

Please charge the power bank for 4 hours before first use for best performance.

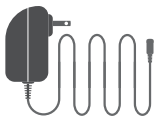12V3A 36W

Can be charged via adapter

TYPE-C / DC5521


Usage Instructions

Thank you for purchasing the Weighted Heating Pad. We designed this heating pad to provide soothing heat therapy to the neck and shoulder area.

The benefits of heat have been clinically proven to alleviate muscle pain, reduce stiffness, and promote healing. This product is designed with two usage options:


You can power it directly with a **power bank** or **adapter** for heating, making it more portable during use.

Power Bank Heating

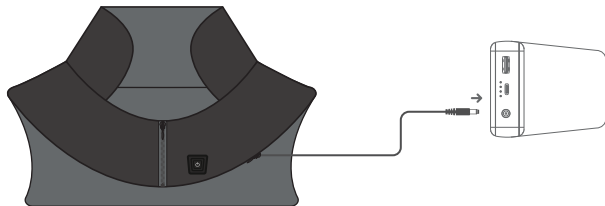
- Connect the power bank and press and hold the "⏻" button for 3 seconds. The indicator light will turn red.
- After turning on, short press the "⏻" button for 1 second to cycle through the temperature settings: low, medium, and high.
- By default, the device will automatically turn off after 2 hours of operation. If you wish to continue using it, you will need to turn it on again.

Adapter Heating

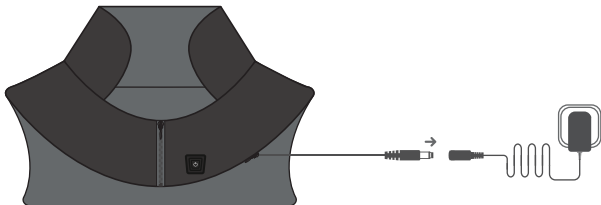- Connect the adapter and press and hold the "⏻" button for 3 seconds. The indicator light will turn red.
- After turning on, short press the "⏻" button for 1 second to cycle through the temperature settings: low, medium, and high.
- By default, the device will automatically turn off after 2 hours of operation. If you wish to continue using it, you will need to turn it on again.

Available for Power Bank Charging

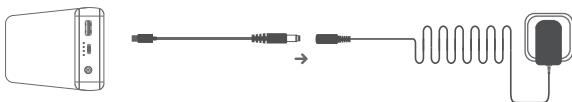Power bank + adapter cable + adapter connection = Charging

Portable Power Bank Charging Guide

I. Connection Steps:

Prepare the original accessories:

Use the included original power adapter

Use the provided DC-to-Type-C cable

Connection sequence:

- ① First, insert the Type-C end into the power bank's input port.
- ② Then, connect the other end of the cable to the adapter.
- ③ Finally, plug the adapter into a power outlet.

II. Safety Precautions:

During charging:

Do not unplug any cables while charging.

Avoid moving or shaking the charging device.

Keep away from high temperatures and humid environments.

After charging is complete:

First, unplug the adapter from the power outlet. Then, disconnect the cable from the power bank.

Power Bank Safety Instructions

I. Charging Requirements

Only use the included original adapter and cable for charging.

Do not use universal or incompatible chargers, as they may cause battery failure, capacity degradation, or reduced lifespan.

The product is not responsible for damage or injury caused by using non-original accessories.

II. Safety Warnings

Never use the power bank beyond the rated voltage—overloading may cause fire hazards.

Do not submerge the product in water or expose it to rain or snow.

If the battery becomes excessively hot, emits unusual odors, appears dented, or shows abnormal deformation, stop using it immediately.

III. Charging Management

Fully charge the power bank before first use.

A completely drained battery requires 4 hours to fully charge. Disconnect the charger once the indicator light stays on to prevent overcharging damage.

Charge in a well-ventilated area—avoid high temperatures (e.g., direct sunlight), humidity, or enclosed spaces (e.g., inside a car).

IV. Long-Term Storage Recommendations

If storing for extended periods, maintain ~50% battery level—avoid full or empty charge.

Perform a full charge-discharge cycle every 3 months to preserve battery health.

V. Prohibited Actions

Never charge the power bank while simultaneously using it (discharging).

Avoid drops, impacts, or compression to prevent internal short circuits and fire risks.

Do not disassemble, modify, or attempt to repair any part of the product.

Dispose of used batteries properly—do not discard in household trash or incinerate.

Keep away from metal objects during storage to prevent short-circuiting.

Product Features

The Weighted Heating Pad for Neck features a secure, ergonomic design that contours perfectly to the body. Its weighted construction allows for deeper heat penetration, while adjustable settings provide soothing warmth to relieve muscle tension and discomfort.

Usage Areas

This weighted heating pad is designed to fit the neck and shoulder area, offering soothing warmth to alleviate muscle tension, stiffness, and discomfort. The ergonomic design ensures a perfect fit, allowing the heat to penetrate deeply for fast relief. It can also be used on the abdomen and legs, providing versatile heat therapy for full-body relaxation.

Shoulder and Neck

Abdomen and Legs

Frequent Questions And Answers

Q:What should I do when my battery won't start?

A: Try charging the power bank with an adapter and a conversion cable for at least 5 seconds to check if the battery starts functioning normally again.

Q: My heating product isn't heating properly.

A: Try the following tests to see if it works normally:

*Check the power bank's charge level and ensure the power cable is properly inserted, press and hold the power button to turn it on and see if it can be used normally.

*Use the included adapter to connect to the heating pad, then press the power button again.If it still doesn't light up, unplug the power and contact customer support.

Q: Why does the heating pad still not feel hot enough after heating for a long time?

A: Our products are tested according to safety standards, under normal room temperature (25°C),the heating element temperature meets the product's design specifications. In different usage environments,people's perception of temperature can vary depending on the environment. Placing the heating pad directly against your skin will still provide warmth.

WASHING INSTRUCTIONS :

NEVER submerge cords, power supply, or control unit.

NEVER use bleach.

Do not use wringer!

NEVER dry clean this heating pad for neck.

Cleaning solvents may damage the wiring insulation.

For hand washing, wash the heating pad for neck by hand in cold water using a small amount of mild laundry soap or detergent.

NEVER rub or scrub.

Rinse thoroughly in cool water until all traces of soap or detergent are gone.

Squeeze excess water out by hand.

NEVER iron your heating pad for neck.

For line drying, dry your heating pad for neck by simply draping it over two parallel clotheslines or a shower rod.

IMPORTANT: DO NOT USE HEATING PADS FOR NECK WHEN WET.

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