

# 21x14cm



## FOOT MASSAGER

360° ALL DIRECTIONAL FOOT MASSAGING

### Contact us

#### THANK YOU FOR PLACING YOUR TRUST



**SATISFIED**

Thanks for choosing us. We're so glad that you're satisfied. Here we sincerely invite you to write an Amazon comment about your happy experience.

Your satisfaction always stimulates us to do better.



**NOT SATISFIED**

Please put your trust in us. Just send us an email with any problem that you have been confronted with.

Our professional service team will never fail your expectations.

12

### Table of contents

Safety Instructions	2
Warranty	4
Product Components	5
Features and Use Diagram	6
Settings	7
Operating	8
Specifications	9
Troubleshooting	10
Maintenance	11

1

### Safety Instructions

#### ► Important Safety Instructions

- Be sure to read the user manual in a detailed way prior to use to ensure safety and proper use of the product
  - Make sure to use the product by following the instructions in this manual
  - After using or before cleaning the product, please disconnect to the power.
  - Always cautiously use the heating function to avoid serious burns.
  - Do not apply the product to too insensitive and over sensitive skin or to body parts of poor blood circulation.
  - Never allow children or those, of physical and mental disorders to use the heating function or improper operation that may be dangerous.
  - To prevent fire, electric shock or injuries caused by overheating, never use it under a pillow or a blanket.
  - Do not cover it during use so as to avoid malfunction or damage resulting from motor overheating.
  - Children, those of physical and mental disorders and the disabled should be under guidance when using the product
  - Please stop using it immediately if you feel uncomfortable.
  - When the power cords or plugs are damaged, inappropriate operation is performed and the product has been dropped or damaged or dropped into water, do not use this product please contact the service center for inspection and repair.
  - Do not use spare parts that are not recommended by the manufactures.
  - Do not insert or drop objects into any opening of the product.
  - Please close all operational controls before disconnecting the product to the power if you want to cut off the power.
  - This product cannot be used for self - diagnosis or treatment and please visit professional medical care if in need of diagnosis or treatment.
- #### ► Use Environment
- Do not use the product in a wet or dusty environment in order to avoid product failure.
  - Do not use this product at room temperature over 40 °C
  - Do not use the product outdoors.
  - Do not use the product next to an electric heater or scores or in direct sunlight.

2

### Maintenance

#### ► Storage

- The product should be put into the color box and stored under safe, dry well-ventilated and cool condition when it is not in use.
- Do not twist switching the power button and the cord.
- Do not hang the product with the cord
- Avoid scratching by sharp objects on the surface of the product.
- Do not put the product indirect sunlight or under condition of high temperature

#### ► Cleaning

- Make sure the power button is off and unplug the power cord before cleaning
- Use a dry cloth to clean the surface of the product
- Do not hang the product with the cord

#### Note

- Do not use corrosive detergents, such as benzene or thinner, to clean the product.

#### ► Maintenance

- Keep the product clean and away from high humidity conditions.
- Please use a dust cover (not include in the box) to cover the product to prevent dust from coming into it if it is not used for a long period of time

11

### Troubleshooting

Please refer to the following several tips and possible solutions if problems occur during usage of this product. The product may be returned to the service center or undergone inspection and repair by authorized technicians.

Problems	Possible Causes	Remedies
Not working	1. Unplugged 2. Power button not pushed	1. Connect the product to power or check whether or not the product is well-connected 2. Push the power button
Suddenly stops running	1. Power disconnection 2. Operation duration over 15 minutes 3. Product overrun and overheat 4. Too much operating pressure that the massage balls stop rotation	1. Check whether or not the product is 2. Restart the product 3. Restart the product after cooling. 4. Restart the product

10

### Safety Instructions

#### ► People Not recommended Using the Product

- Please consult your doctor before using if you have the following situations
- Receiving medical treatment or experiencing physical discomfort
- Suffering from hypertension, heart disease, cerebral arteriosclerosis, and cerebral thrombosis
- Suffering from cancer
- Pregnant or in a period
- Suffering from osteoporosis
- Using a pacemaker or implantable medical devices
- Suffering from acute diseases
- Minors

#### ► Buttons

- Do not operate the button with wet hands.
- Keep the buttons from water or other liquid
- Do not fall asleep when using the product
- Do not use this product when you are drunk or feeling uncomfortable
- Do not directly contact the massager pad with your hands, feet joints, knees chest, and abdomen
- Do not use this product over thirty minutes in a row
- Do not use if the product's cloth cover has been damaged.
- Do not use the product one hour before and after eating so as to avoid adverse effects
- Do not enable the heating function when the product is being used in a low room temperature. It is recommended that you should gradually adjust the temperature to a normal range before enabling heating.
- If the product is used immediately after being transferred from a cold environment to a warm environment, the accumulation of water vapor inside the product may affect the normal operation of the product or even cause mechanical failure. You are suggested to put the product in normal room temperature for 1 hour before starting to use it after the transfer

3

### Warranty

- Warranty is not transferable.
- Users can enjoy free warranty from the date of purchase during the warranty period (The surface is not covered under warranty.)
- The product warranty is valid for 12 months from the invoice date.
- If we write that supplier shall provide a warranty card.
- The purchasing date shall be based on the effective receipt or, invoice issued by the retailers
- This warranty shall not apply to leasing, trade, commerce, institutions and other non-buyer users.
- The following situations are not covered by the warranty
  - Destruction caused by human factors, including the non-normal using environment, usage and storage that are not compliant to the instructions, and damage caused by improper use, etc.
  - Disassemble, repair and modification of the product without our authorization
  - Damage caused by users poor transportation after Purchase;
  - Damage resulting from any other force majeure (such as flood, fire, lightning and earthquakes)

4

### Specifications

Name	Airbag Heating Foot Massager
Power supply	110V~220V /50Hz
Power consumption	50w
Dimension(L*W*H)	42*36*24cm
Weight approx	4.7kg

9

### Features and Use Diagram

#### ► Product Features


- Egg-shaped design, fashionable and beautiful, lightweight and portable;
- A combination of a variety of function including scratching, roller kneading and push kneading acupressure, air pressure, and heating;
- Massager Pad that provides all-around massages and in -deep acupressure to clear the body's meridians and foot reaction zones;
- Removable and washable cover, allowing you to keep clean at any moment.;
- Cord wrap, bringing convenience to life

#### Use Diagram



6

### Settings

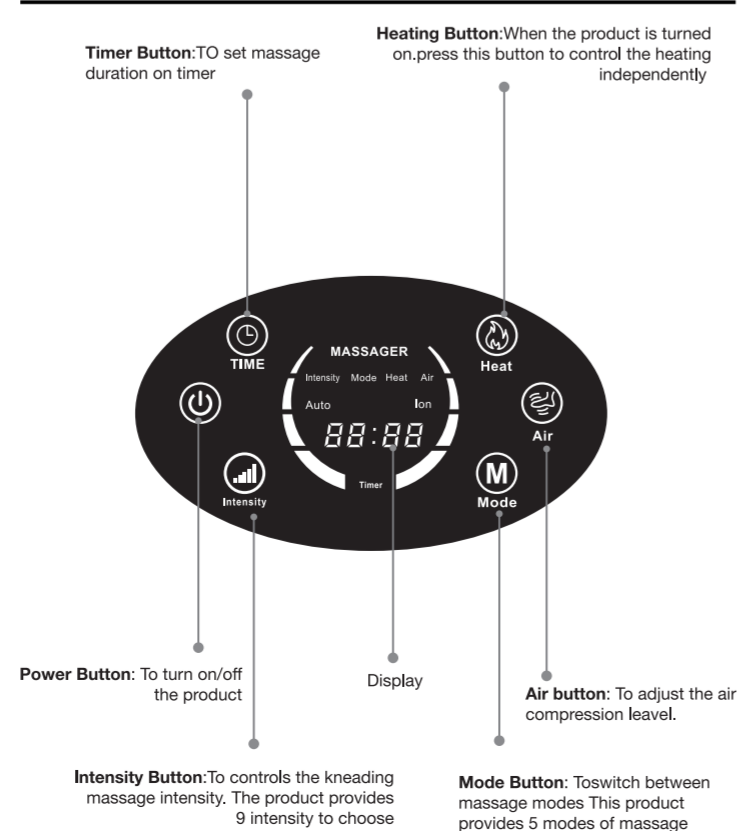


7

### Operating

- Place the foot massager on the floor in front of a comfortable chair and plug it into the mains outlet. Sit down in a comfortable position.


- Power on the unit by turning on the power switch
- Set time by timer key, can be adjusted to 10-15-20-25-30 minutes. The default time is 15 minutes.
- Set intensity by intensity key. There are 9 levels to meet different demands. The default intensity is level 5.
- Set air compression level by air key. There are 2 levels.
- Set heat on/off by heat key.
- Set mode by mode key. There are 5 modes.

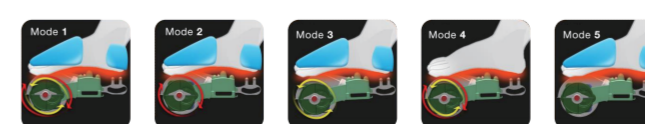

- Mode1: Air on +heating on + motor two direction turns
  - Mode2: Air on +heating on + motor clockwise turns intermittently
  - Mode3: Air on + heating on + motor Counterclockwise turns intermittently
  - Mode4: Air off + heating on + motor two direction turns
  - Mode5: Air on + heating on
- Turn off the power switch to switch off the appliance completely and unplug the appliance from the wall outlet after use.

8

### Product Components

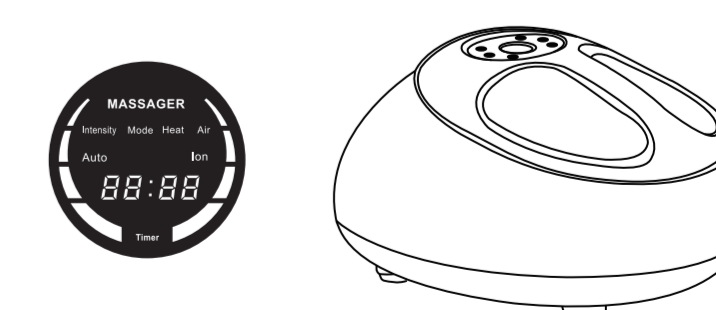

- Central Surface: Six independent buttons - controlling the product's settings and functions
- Central Massager Pad: Scratching-Acupressure on the arches with heating
- Front Massager Pad: Roller Kneading-Acupressure on the front soles
- Rear Massager Pad: Push Kneading-Acupressure on the heels
- Interior airbag-Air pressure on the insteps
- Removable and washable fabric cover-cleaning
- Bottom: Cord Wrap-keeps things tidy

5