Upgrade · The Home Creator Way

Tub Transfer Benches

Model: FBL630264

MODEL:FBL630264

<Picture Only For Reference >

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

IMPORTANT SAFEGUARDS

WARNING: Read and understand this entire manual before operating or servicing this product. Failure to follow these warnings and instructions can cause personal injury or damage to valuable property.

ASSEMBLY PRECAUTIONS

1. Assemble only according to these instructions. Improper assembly can create hazards. The manufacturer declines all responsibility in case of damage due to a bad installation.
2. Caution with the risks of pinching or cuts during assembly.
3. Avoid children using the product. And this product is not a toy. Do not allow children to play.
4. Wear ANSI-approved safety goggles and heavy-duty work gloves during assembly and use.
5. Keep the assembly area clean and well-lit.
6. Keep bystanders out of the area during assembly.
7. Do not assemble when tired or when under the influence of alcohol, drugs, or medication.
8. Product capabilities apply to a properly and completely assembled product only.
9. The product complies with relevant standards and regulations. Installation of this device must comply with all applicable regulations, and it should only be operated in well-ventilated areas.
10. This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
11. Please properly keep product. Check whether the product is abnormal before use. If it is abnormal, please stop using it.
12. Please use the product according to the intended use, and warn of potential injury from misuse.
13. For additional information regarding the parts listed in the following pages, please refer to the Assembly Diagram of this manual. Unwrap and separate all parts in a clean work area. Please keep small spare parts out of children's reach.

SAVE THESE INSTRUCTIONS

WARNING:


- Read the manual completely to avoid harm, product damage, or warranty avoidance. Seek professional help if you encounter discomfort, pain or worsening symptoms. If necessary, use supervision in order to assemble and use this product correctly and safely. Do not exceed the weight limit.
- Regularly check for secure assembly and contact support if in doubt. Ensure the product's legs are even and locked before use, and discontinue usage if there are signs of wear, damage, or missing parts. The product is to be used on flat, stable surfaces. Do not disassemble if not necessary, as it might compromise the longevity of the product. Do not modify the product unless specified in the manual.
- Exercise caution during operation and assembly. Avoid extreme temperatures or open flames. If an allergic reaction occurs, cease use immediately and consult a healthcare professional.
- Clean regularly using a mild detergent and warm water, rinse off any soap residue. After cleaning, dry thoroughly
- Not suitable for environments with high humidity, explosive atmospheres, salt spray, acidic or alkaline conditions.

IMPORTANT


- 1) Don't exceed the maximum weight capacity.
- 2) Only for use on bathtubs that are stable and have flat surfaces on both sides.
- 3) Regularly inspect for signs of wear, damage, or loose parts.
- 4) Before use, confirm with your healthcare provider that the product meets your specific needs. The following guidelines are not a substitute for medical advice.

PART LIST

When unpacking this product, check to make sure the parts listed below are included and carefully inspect for any damage that may have occurred during transit. If any part is missing or damaged, do not attempt to assemble or use the product.

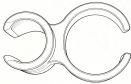NO	PART NUMBER	PICTURE	QTY	NO	PART NUMBER	PICTURE	QTY
#1	Backrest		1	#6	Chair foot B		2
#2	Backrest connect tube		2	#7	Armrest		1
#3	Trestle A		1	#8	Shower buckle		1
#4	Trestle B		1	#9	Manual		1
#5	Chair foot A		2				

DIMENSIONAL PARAMETERS

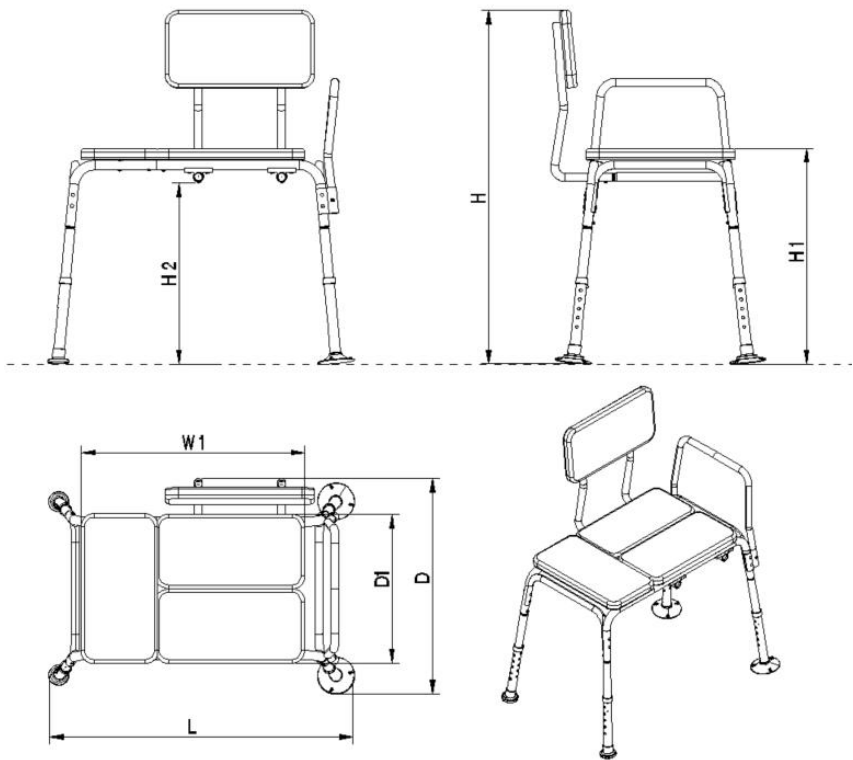Unit:(mm/in)			
L	D	H	D1
800-810 (31.5-31.9)	530-540 (20.9-21.3)	820-930 (32.3-36.6)	400(15.7)
W1	H1	H2	
600(23.6)	485-600 (19.1-23.6)	385-500 (15.2-19.7)	


※There may be a ±20mm error in the measurement of some dimensions, which is normal and not a product quality issue.

PAREMETERS



Model	FBL630264
Main Materials	Aluminum
Seat Rated Load	227kg/500lbs

INSTALL INSTRUCTIONS

Step1.Place trestle A and trestle B upside down on a plane, and insert trestle A into trestle B.

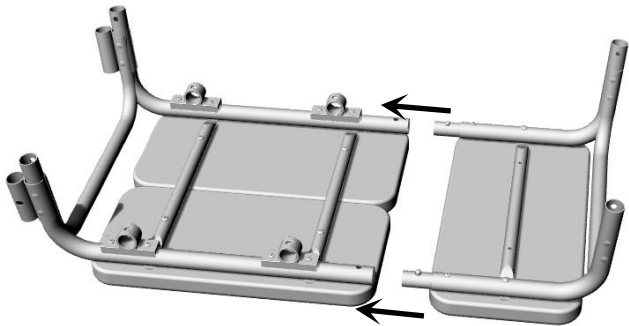
Step2.Insert the Chair foot(A&B) assembly into the trestle A and the trestle B.

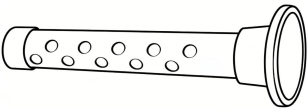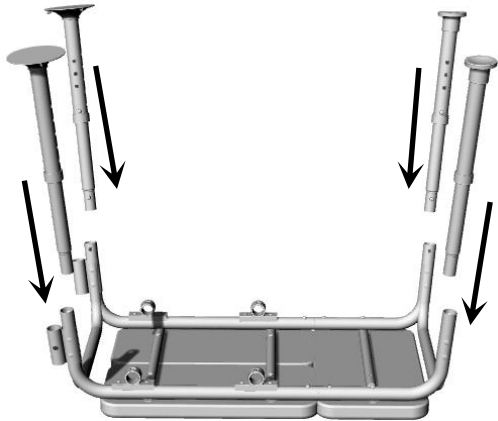
The chair legs are adjustable in 10 gears, and users can adjust the position of the chair leg ball hole to meet the height requirements of the seat.

Step3. Insert the backrest assembly backrest connection pipe into the backrest assembly.

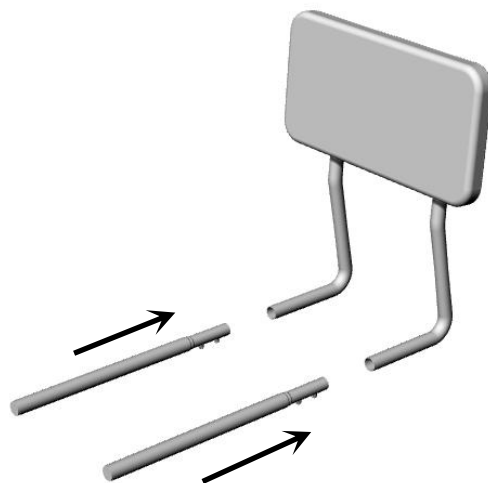
Step4. Insert the backrest assembly into the support hole of connection pipe under the seat.

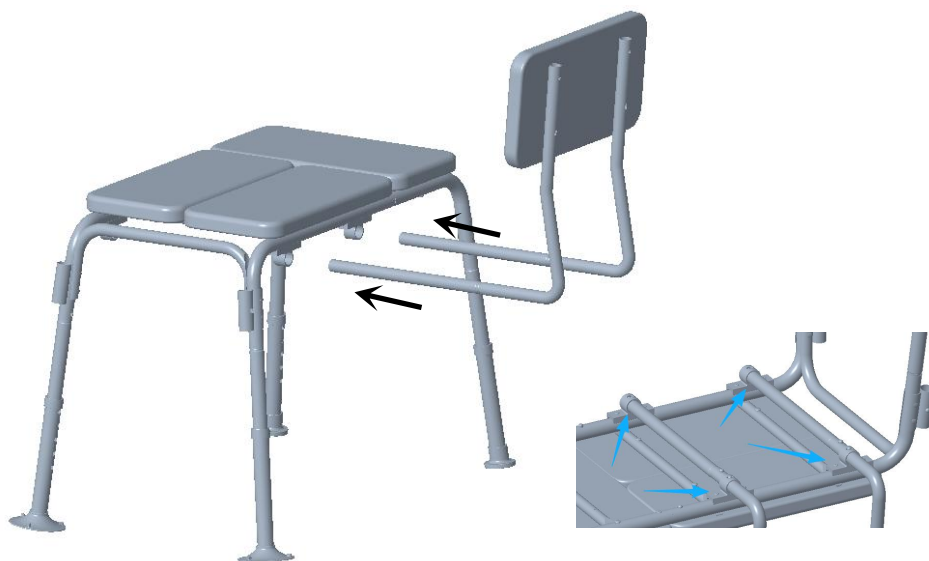
Step5. Insert the armrest into the pipe of trestle B.

Step6. Fix the shower buckle to the armrest.

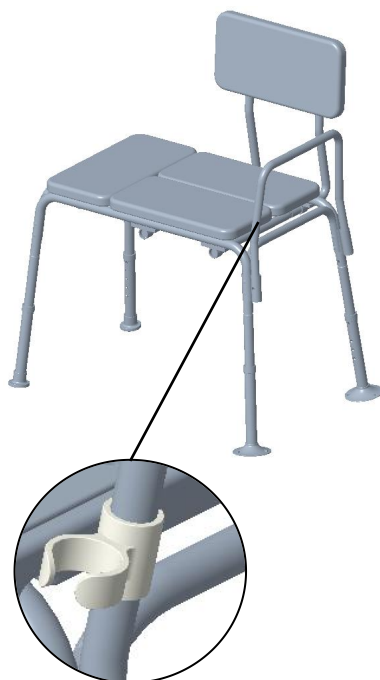Leveling

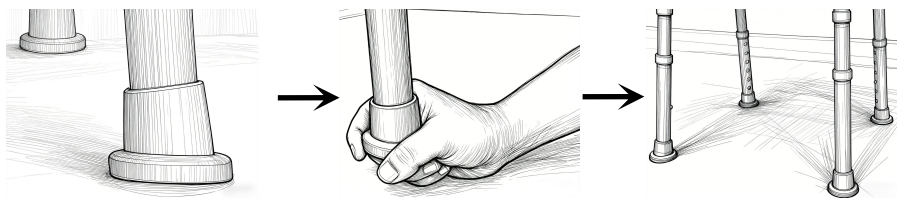If the foot mat is not placed flat on the ground, rotate it until it is fully in contact with the ground.

USE INSTRUCTIONS

USING WARNING:

- (1) Put the shower bench stably on a flat floor.
- (2) Before use, please make sure that the screw has been tightened and the ball pin of the legs has been ejected from the slot.
- (3) Before use, adjust the shower bench to an appropriate height so that when the two legs are supported in the bathtub, the seat is smooth.
- (4) Check the stability of the shower bench before using it every time. Stop using it if the legs are shaking or the shower bench is unstable.
- (5) Avoid the overload.
- (6) The shower bench should not be used to push the user while seated.
- (7) Shower bench must be used ONLY for taking a shower inside bathroom, not for any other purposes.
- (8) Adjust the shower bench legs to a suitable height so that the user's feet from toes to heels comfortably touch the floor when sitting on it.
- (9) Stay away from open flames.
- (10) Periodically check whether screws have been tightened.
- (11) Ensure that rubber tips at the bottom of the shower bench's legs are not worn out. If not, replace them immediately.
- (12) Please read carefully before use and do not tear off the label.

MAINTENANCE

- Wash with clean or mild soapy water and air dry after use.
- When not in use, please store in a cool and dry place
- Check the shower bench regularly.

STORAGE CONDITIONS

* The fabric should be stored completely dry, in a clean and dry place to avoid the risks of mold.

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