Upgrade · The Home Creator Way


VEVOR Support Center

FLAT THRESHOLD RAMP

MODEL: JJLXMKPD001

FLAT THRESHOLD RAMP

MODEL: JJLXMKPD001

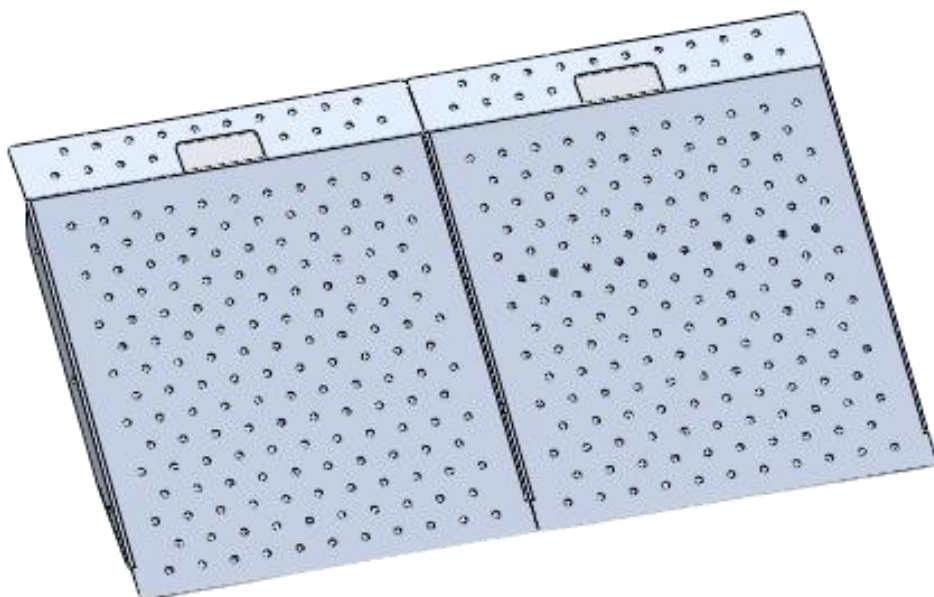(The picture is for reference only, please refer to the actual object)


This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

Warning-To reduce the risk of injury, user must read instructions manual carefully.

PRODUCT DESCRIPTION


Thank you very much for choosing this Flat threshold ramp. Flat threshold ramp is primarily used for accessing thresholds and stairs during the transportation of wheelchairs, carts, and large items. They can be fixed to the ground, and some models are foldable for storage and easy portability. Please read all of the instructions before using it. The information will help you achieve the best possible results.

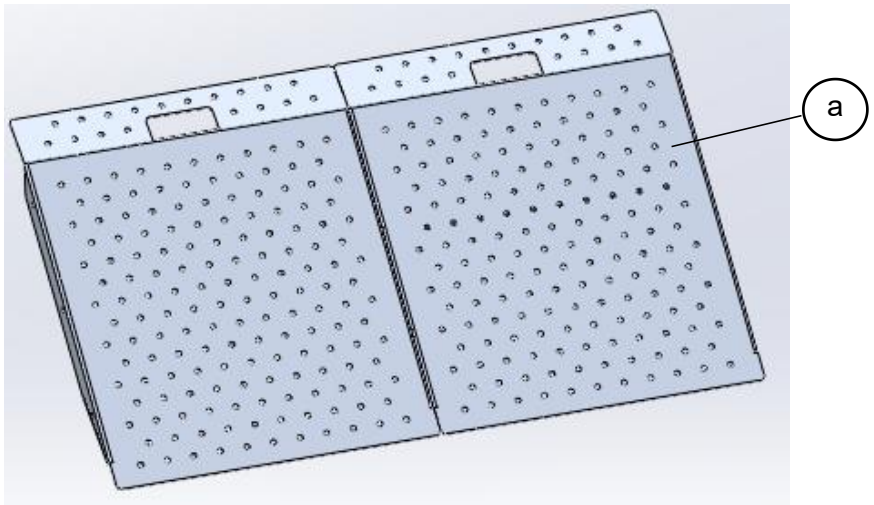
A: Flat threshold ramp

SAVE THIS MANUAL

Flat threshold ramp	
Model	JJLXMKPD001
Max Weight	750lbs
size	750*617mm

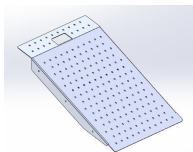×2

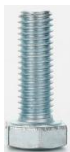×3


×3


OPERATION

Installation:

Step 1: Connect the two together with screws.

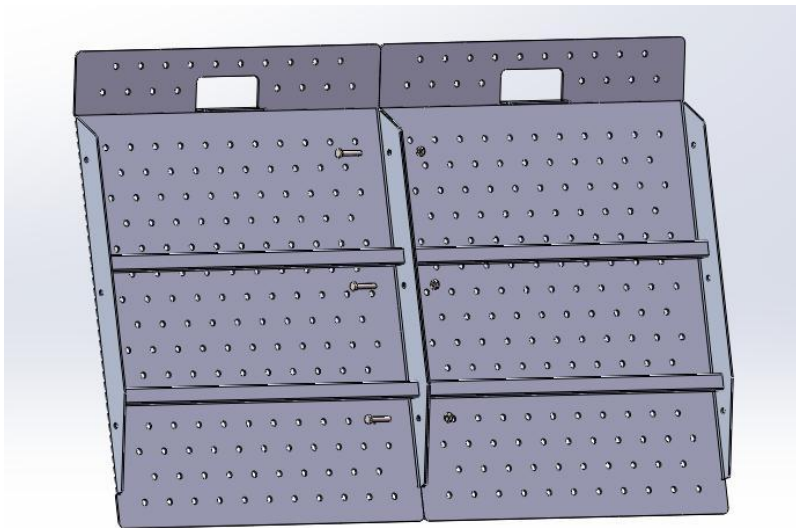Step 2: You can start using it.

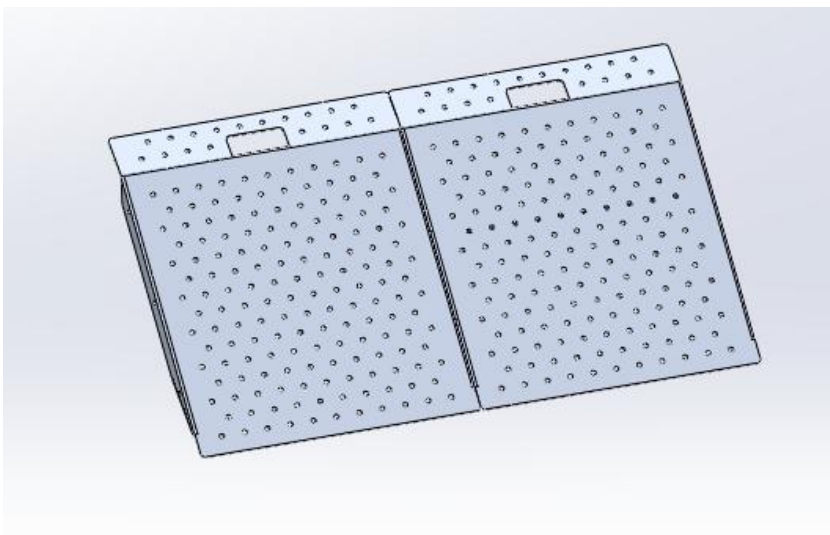SAFETY NOTES AND WARNINGS

Important Safety Instructions:

WARNING:

(1) The product complies with relevant standards and regulations.

Installation of this device must comply with all applicable regulations, and it should only be operated in well-ventilated areas. Please read the user manual before installing and using this equipment.

- (2) This product is mainly used by users to push wheelchairs, strollers, and carry large objects across thresholds.

1.AVOID SERIOUS INJURY OR DEATH

2.DO NOT Exceed the weight capacity or load height of the ramp.

3. Use Ramps on a stable and level surface only.

4.It is strictly prohibited for users to drive the wheelchair alone on the ramp.

5.Only use the ramp with a qualified helper and always use a lap belt.

6.When carrying heavy objects through a ramp, be sure to have enough auxiliary personnel to assist in ensuring that people or goods do not tip over and cause damage to people or goods when using the ramp.

- (3) Before use, refer to your mobility equipment's owner guide for acceptable ramp usage, including proper slope/incline, chair direction; never exceed its recommendations.

- (4) Do not use the ramp if it is damaged or unstable.

- (5) Do not tamper with, attempt to repair, or modify any portion of the ramp.

- (6) Use caution when folding and unfolding the ramp to avoid all potential pinch points.

- (7) Use it only if the secure, unobstructed, level landings can support, at minimum, the same weight capacity as the ramp (800lbs.)

- (8) The ramp surface may be slippery in wet conditions; use caution if the surface is wet as traction may be reduced

- (9) Do not use the ramp if its surface is covered with ice, snow, dirt, leaves, or other debris. Remove all accumulations and sweep the surface clean before use.

- (10) Improper use of the ramp could result in serious injury.

- (11) Always exercise caution when handling, installing, and or using the ramp.

- (12) Casual placement may result in accidental injury due to children's climbing or playing; or the product may fall off and cause major accidents. During walking, the product may shift or dislocate due to the

reaction of your walking. Please be careful.

- (13) Please do not disassemble and reinstall this product, Otherwise, it may cause serious accidents.
- (14) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- (15) Children shall not play with the appliance.
- (16) This appliance shall not be used by children. Keep the appliance and its cord out of reach of children.
- (17) Please properly keep product. Check whether the product is abnormal before use. If it is abnormal, please stop using it.
- (18) Please use the product according to the intended use, warn of potential injury from misuse.

MAINTENANCE AND TROUBLE SHOOTING

- (1) Place the product on a sturdy, secure, accessible, and level starting point. The minimum weight capacity of the starting point must be equal to that of the ramp (750 lbs). Do not exceed the maximum weight capacity of the product.
- (2) Ensure that the product is securely flush against the starting point. If a safe, accessible, and level starting point cannot be found, or if the product does not fit securely against a higher starting point, do not use the ramp.
- (3) Do not tamper with, attempt to repair, or modify any part of the product. Exercise caution to avoid pinch points during assembly and handling.
- (4) The ramp surface may become slippery in wet conditions. Walk with caution and pay attention to your footing to prevent slipping.
- (5) If using the product with other mobility equipment, refer to the user guide of the equipment to understand the acceptable usage

guidelines, including appropriate slope, incline, and chair direction. Do not exceed the recommendations of the mobility equipment.

- (6) When using the ramp, ensure it is securely flush against the starting point to ensure safety.
- (7) Before each use, clean the ramp surface of ice, snow, dust, leaves, and other debris. Remove accumulated ice and snow, and sweep the surface clean. For dust and debris, use a soft bristle brush and non-alkaline detergent to clean the ramp surface, then rinse thoroughly.
- (8) If the ramp will not be used for an extended period, store it in a dry and well-ventilated area away from direct sunlight.
- (9) Do not use thinner, gasoline or other volatile organic solvents to clean the product.
- (10) Do not scrub this product with a hard brush or hard sponge.

Upgrade · The Home Creator Way

Centre d'assistance VEVOR

FLAT THRESHOLD RAMP
MODEL: JJLXMKPD001

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