


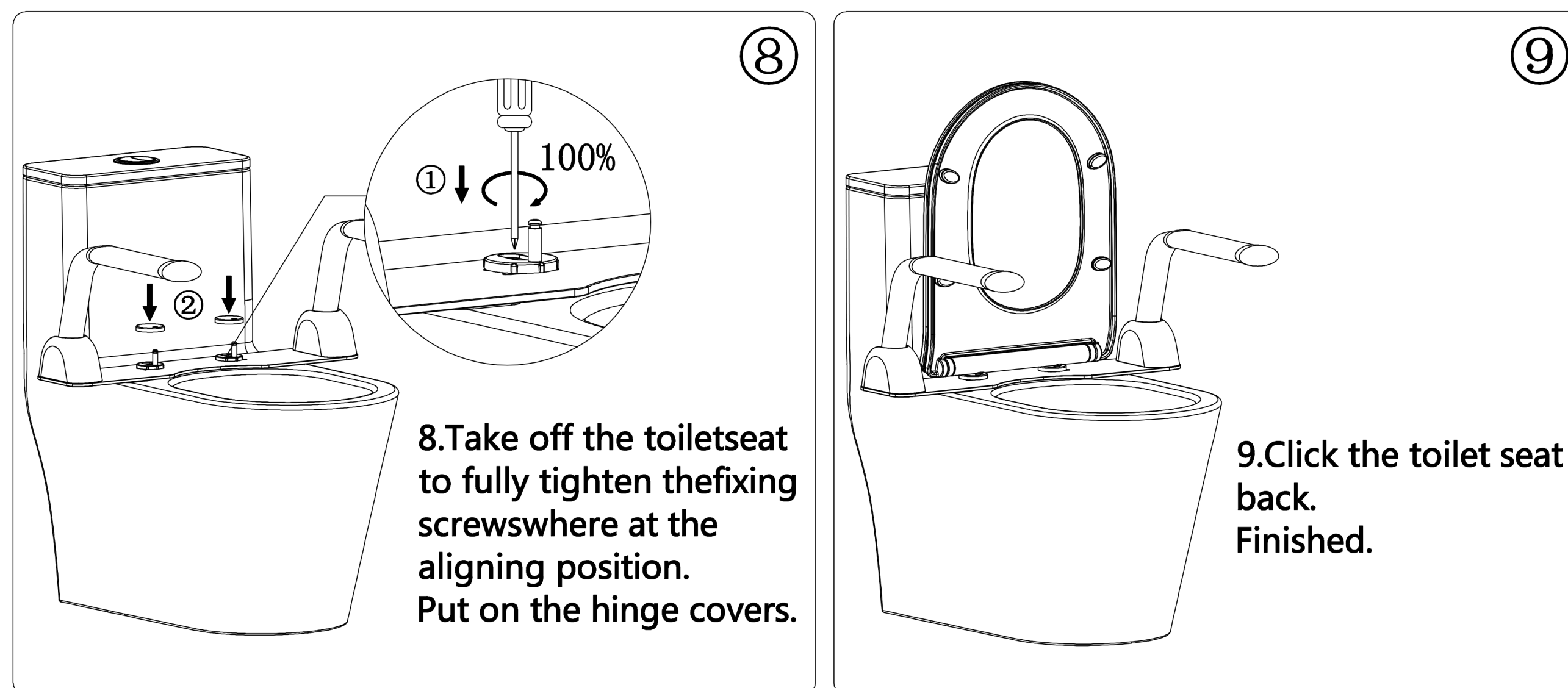


INSTALLATION INSTRUCTION

NOTE

- Please read the instructions carefully to understand our toilet handrails and prepare the necessary tools for installation in advance. Correctly install the toilet handrail according to the instructions. This will help avoid hurting you, your family and your property.
- Please keep the instructions available after reading it carefully.
- Please subject to the real products when the instructions are different from the real products. If there is any new update to products and its accessories, we will give no further notice.


NOTES

- To prevent the lock nut from being pulled out or loosened during use of the arm rest, try to secure the hinge base fixation kit as tightly as possible without causing the nut to be stripped.
- To avoid the Toilet Seat from falling to one side, or being tilted up significantly, do not forcefully shake or put excessive pressure ($F1 \leq 80\text{kg}$) onto one side of the arm rest after the arm rest is correctly installed in accordance herewith.
- To prevent damage to the lifting handrails, the dual handrails simultaneously bear ($F2 \leq 120\text{kg}$).
- The figures above are for reference only.

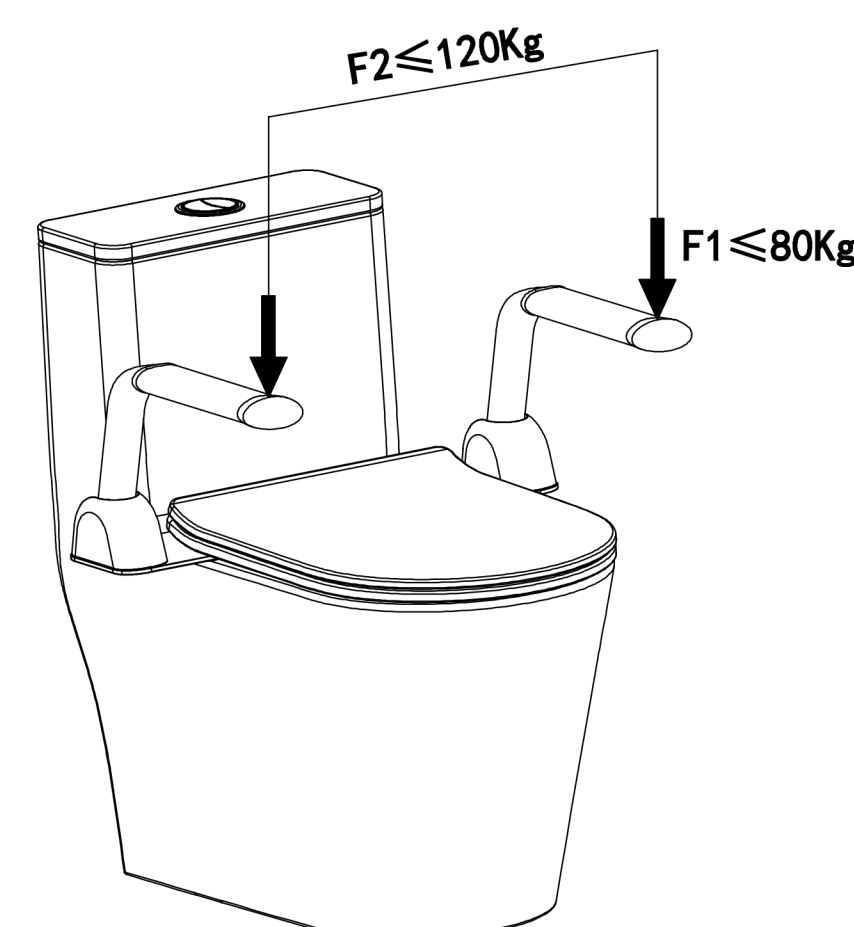


CARE AND CLEANING

- Use a mild detergent such as liquid dishwashing soap and warm water for cleaning. Do not use abrasive or high density cleaners that may scratch or dull the surface.
- Carefully read the cleaner product label to ensure the cleaner is safe for use on the material.
- Always test your cleaning solution on an inconspicuous area before applying to the entire surface.
- Wipe surfaces clean and rinse completely with water immediately after cleaner application. Rinse and dry any overspray that lands on nearby surfaces.
- Use a soft, dampened sponge or cloth. Never use an abrasive material such as a brush or scouring pad to clean surfaces.

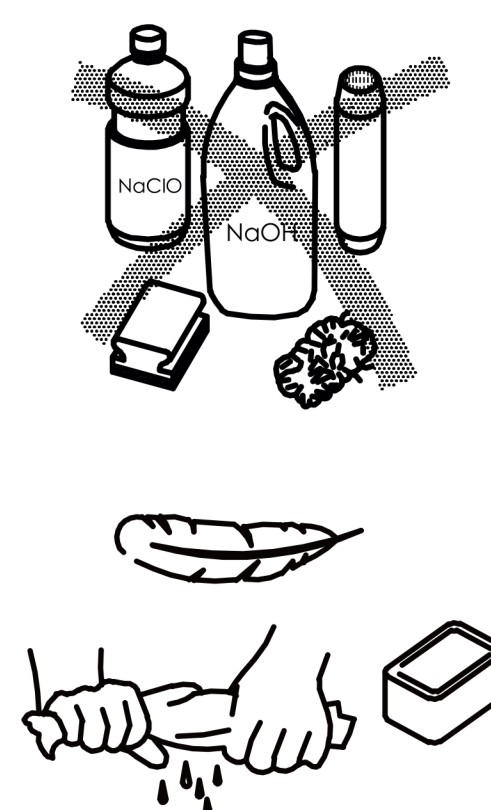

Toilet handrail and Accessory set

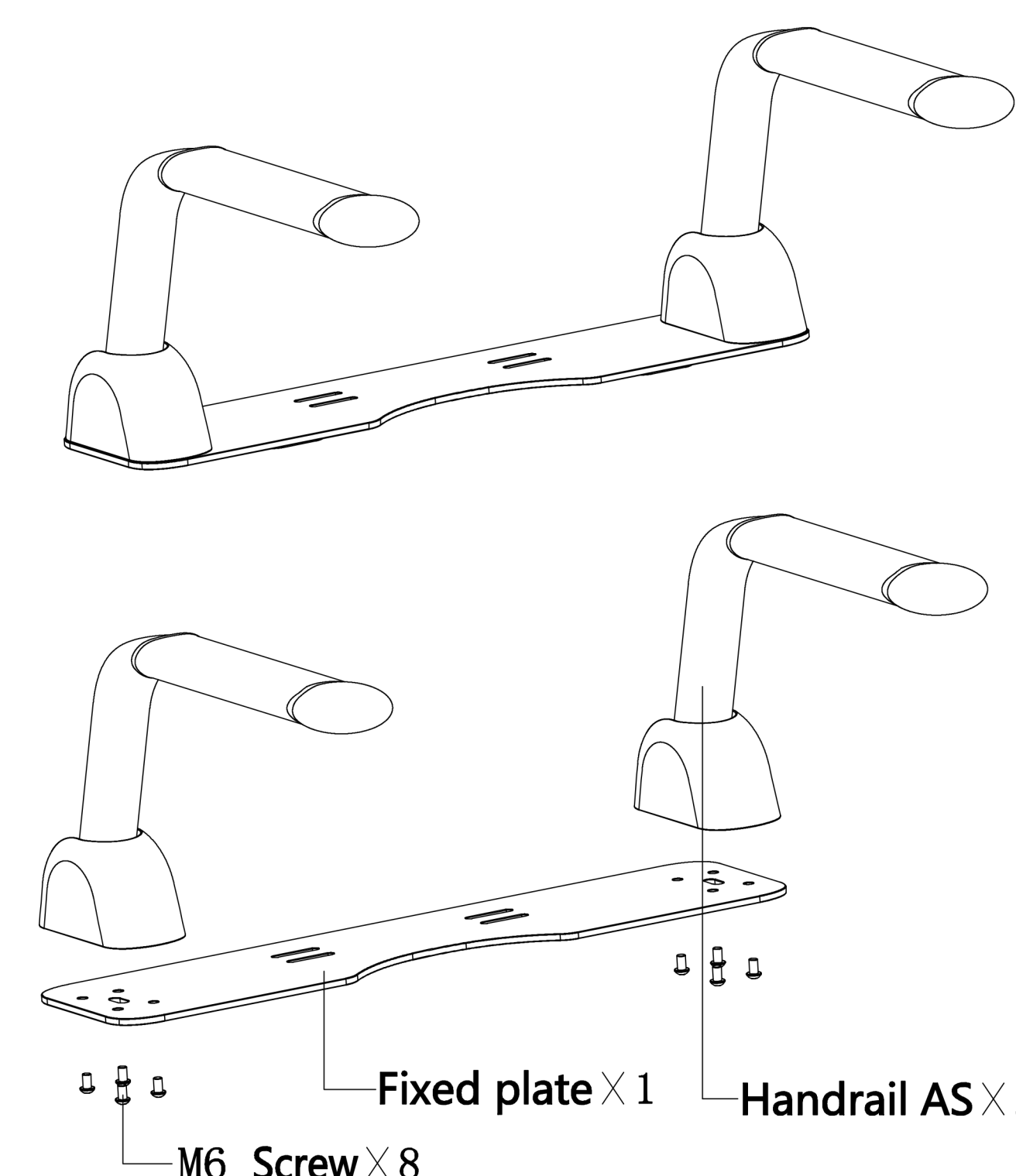

Parts	
	M6 Screw × 8
	Allen Key × 1
	Instructions × 1
	Shim × 2
	Integrated Screw Kit × 2


Intallation Tools Needed

Screwdriver	Soft wet cloth

Cautions

As shown in the following figure, a toilet seat suitable for handrails must meet the following requirements the requirements below:

1. To prevent interference between the arm rest end and the toilet tank, the area of $A \geq 35\text{mm}$ shall be flat.
2. Distance between the installation holes in the toilet seat:
 $130\text{mm} \leq B \leq 200\text{mm}$.
3. Space distance required for installation of toilet armrest: $C > 570\text{mm}$.

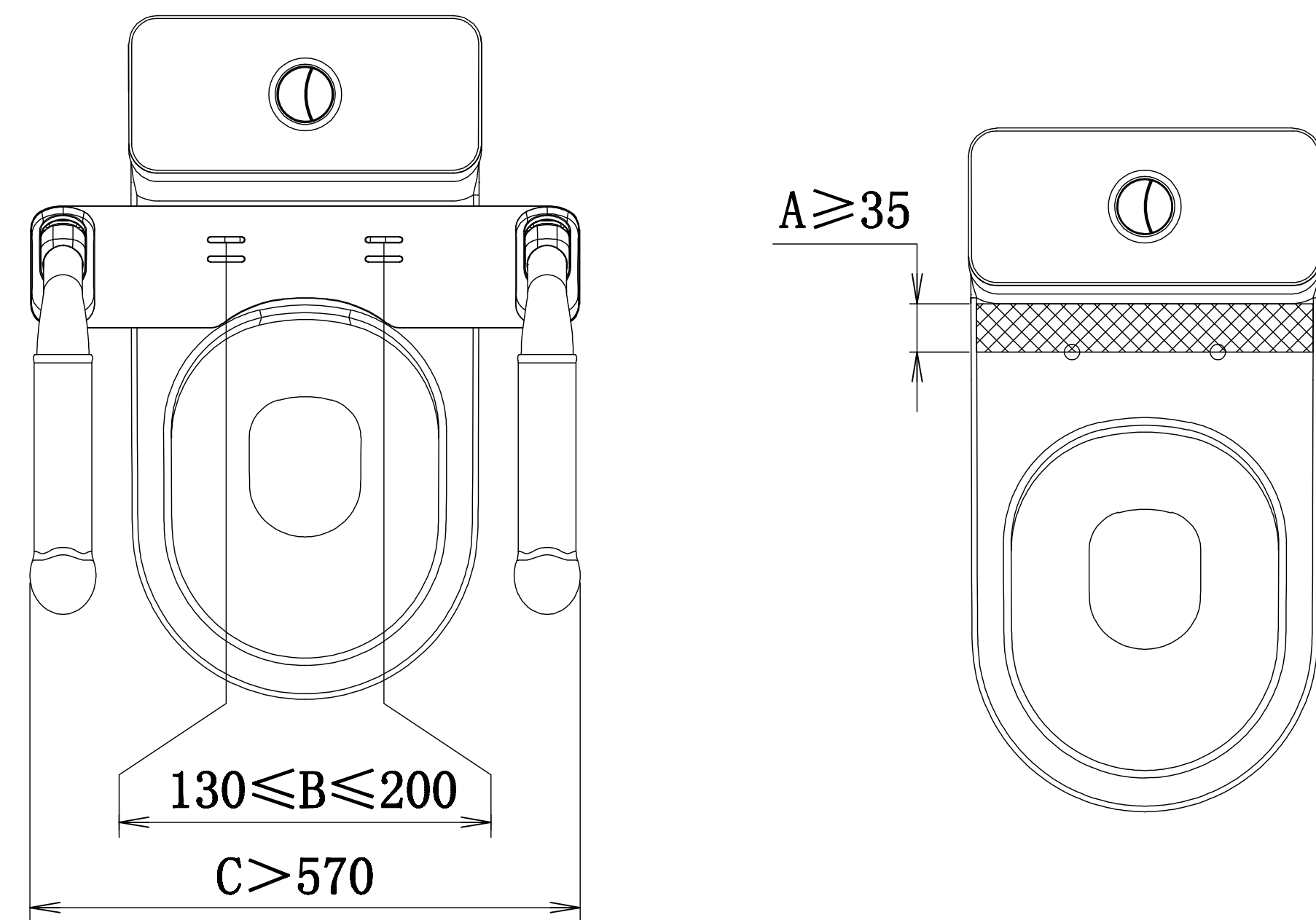

INSTALLATION

