

USER'S MANUAL

Folding Seat Crutch

THIS INSTRUCTION BOOKLET CONTAINS **IMPORTANT** SAFETY INFORMATION. PLEASE READ AND KEEP FOR FUTURE REFERENCE.

Please give us a chance to make it right and do better!

Contact our friendly customer service department for help first.

Replacements for missing or damaged parts will be shipped ASAP!

Before You Start

- ⚠ Please read all instructions carefully.
- ⚠ Retain instructions for future reference.
- ⚠ Seperate and count all parts and hardware.
- ⚠ Read through each step carefully and follow the proper order.
- ⚠ We recommend that, where possible, all items are assembled near to the area in which they will be placed in use, to avoid moving the product unnecessarily once assembled.
- ⚠ Always place the product on a flat, steady and stable surface.
- ⚠ Keep all small parts and packaging materials for this product away from babies and children as they potentially pose a serious choking hazard.


Introduction

Dear customers:


- Thank you for using our seat crutch.
- Before using it, please read the instruction manual carefully and use it correctly.
- Keep this manual in a safe place for future reference.

Product Introduction

- Seat crutch is a multipurpose and convenient crutch with a seat. When open, it is a triangular chair. When folded, it becomes a sturdy crutch. The three legs with rubber feet are more stable than other crutches. The main structure is aluminum alloy, with durable ergonomic ABS seat board, ergonomic handles, and the seat crutch is adjustable in height so that users can adjust it according to their height.
- Seat crutches bring the crutches a new definition, bring the leisure chair a turning innovation. The perfect combination of the characteristics of crutches and leisure stool, the most significant increase in the safety of crutches and leisure stool, which is the most trusted leisure chair and walker crutches to family.
- Seat crutches suits for leisure-loving old people, elder who wait in line, can help to release their burden, and also fit for people who like to go fishing, and people who want auxiliary functions.

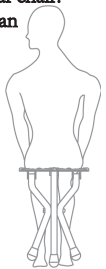

Product Structure

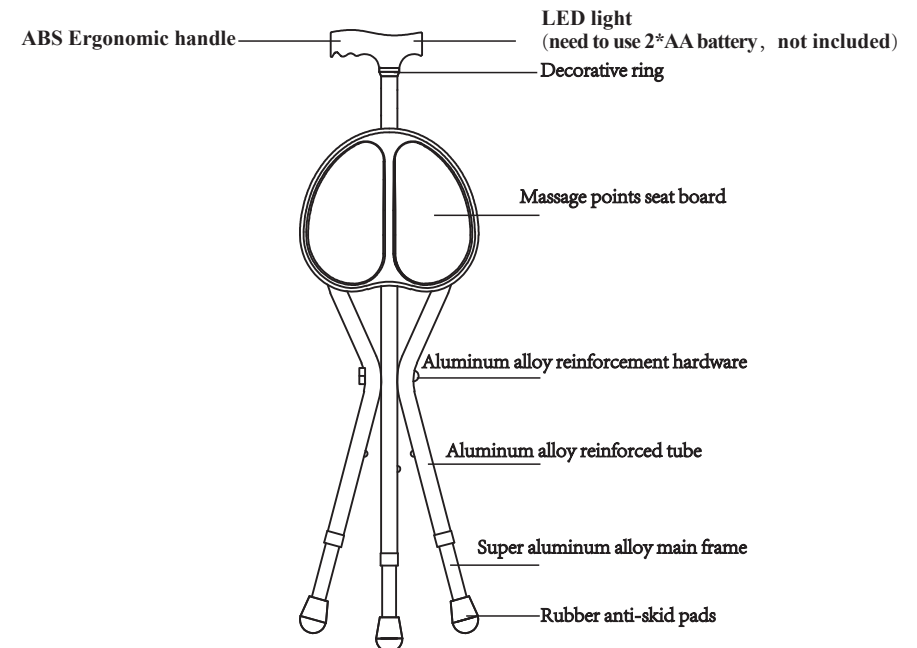

Steps for Usage

How to Use When Folded



First Method :
Hold the handle forward
This kind of usage can massage a lot of acupoints of the hand to achieve a variety of beneficial effect to human body.



Second Method:
Hand grip backhand
This method can massage Hoku, this one is more comfortable than the method one.


Usage When Open the Seat Crutch

Open the Seat Board
When open, pull downward the seat board to the fixed position.
When folding, pull back and remember to fix the seat board to the buckle.




X

Error sitting method

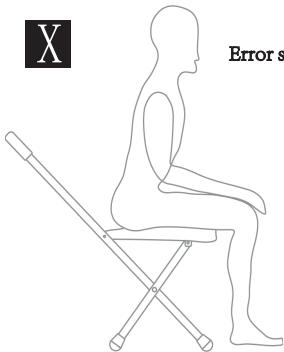

✓

The right usage:
After opening, it is a beautiful heart-shaped seat board. It comes with massage points of which the height is in accordance with the ergonomics. There are two big dots designed to massage Huantiao (GB30) point, to prevent movement series diseases: sciatica, lower limb paralysis, cerebrovascular disease sequelae, lumbago and leg pain, hip and surrounding soft tissue diseases, beriberi, etc.

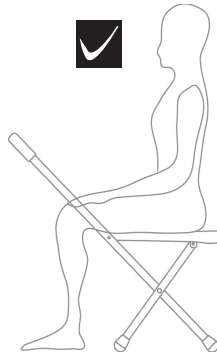

Precautions

Please note the following:

1. The wire on the floor
2. The liquid on the floor
3. Easy sliding carpet
4. The stairs
5. The courtyard of the door
6. The gap of the floor

- * Reward Points
- * Exclusive Customer Service
- * Personalized Recommendations
- * Permanent Shopping Cart
- * Order History

5 Star Rating



We would like to extend our heartfelt thanks to all of our customers for taking time to assemble this product and giving us valuable feedbacks.

THANK
YOU
😊