



■ Help stand function

■ Massage and fever

■ Reclining function

INSTALLATION BOOK

Simple installation manual for heating sofa



We're excited for you to regain your independence with Functional sofa!

To ensure the best experience, please read this guide carefully before using your new Functional sofa.

Be sure to keep this manual close by, preferably in your sofa pocket so that you can always refer to it as needed.

| DIRECTORY |

Safety Rules.....	2&3
Heat/Massage Remote	4
How to use the massage features.....	4
Massage Areas & Assembly Parts.....	5
Assembly Instructions.....	6&7
Help stand function.....	8
Recline Function.....	9
Date Sheet&Troubleshooting.....	10

SAFETY RULES

- ★ Never set the lift chair in wet/damp areas or where dust can accumulate as these could cause an electrical fault or a mechanical dysfunction.
- ★ Never use the lift chair inside your home if the room temperature is over 90°F.
- ★ Never place the lift chair near a heat source(oven,chimney,fernace etc...) or in direct sunlight.

Do Not use the armchair if:

- ① The massage chair is wet.
- ② Your hands are wet.
- ③ Only clean the chair with a dry cloth.
- ④ It's advisable not to use the heat or massage feature for more than 30 minutes at a time and that you drink 2 glasses of water after using it.

Situations where you should not use the massage chair:

Switch the massage chair's massage system off if it gets wet or seems not to be working properly,if any part of the char is damaged or malfunctioning,if any part of the chair appears to be damaged or wiring is exposed,if you feel unwell or feel any pain during the massage,or if ther has been a general electrical fault.



Usage:

- Indoors, avoid high temperatures and wet or damp areas.

Precaution:

- Never sit on the remote or place heavy objects on top of it.
- Do not drop the remote on the floor.
- Do not tug on the electric cable.
- Do not use the remote with wet fingers.
- Do not press or touch the buttons with sharp objects.
- Never unplug the remote without first turning the power off.
- The chair should rest for 30 minutes after it has been used for 30 minutes.

When the massage system of the chair is not in use:

- Make sure it is switched off.
- Unplug the chair's massage system from the wall outlet when you are not going to use it for long periods of time.
- Keep the chair in a clean and dry place.

HEAT/MASSAGE REMOTE



Power On/Power Off

Heating function

Intensity

Massage mode

Select massage area

Back

Lumbar region

Thighs

Legs

Timer

How to Use the massage features

- **POWER:** Press the Power button to start, press it again to Stop.
- **MODE:** With the Mode button, you can choose from 9 different massage systems.


We recommend the default setting when the massage feature is powered On. To select the massage area, just press the button you want to select: back, lumbar region, thighs and legs.

- **INT:** The INTensity of the massage can be controlled, select from high or low.
- **Timer:** You can choose 15Min, 30Min or 60 Min timer for massage and heating.



Massage Areas

The functional sofa has 8 vibrating message points for 2 points in 4 areas: Back, Lumbar, Thighs, Legs. Each Region can be controlled by the massage remote as well as the type of massage pattern and intensity.

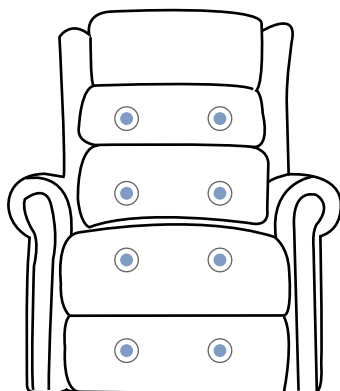

Please check you have all hardware before assembling.


Please ensure you have all the parts and hardware as the instruction shown.

Please contact us by mail, if you have any questions, and we will reply you ASAP.


Enjoy your new purchase! .

Assembly Parts:


A.chair base

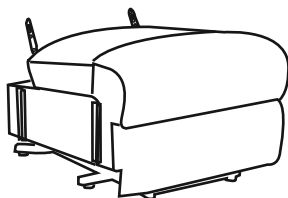


B.backrest

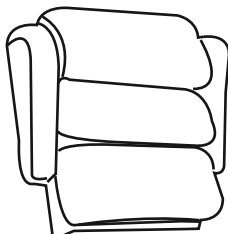

C.armrest

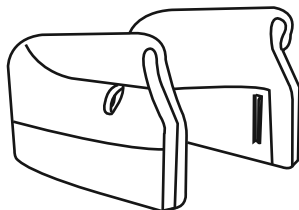


Supplies:

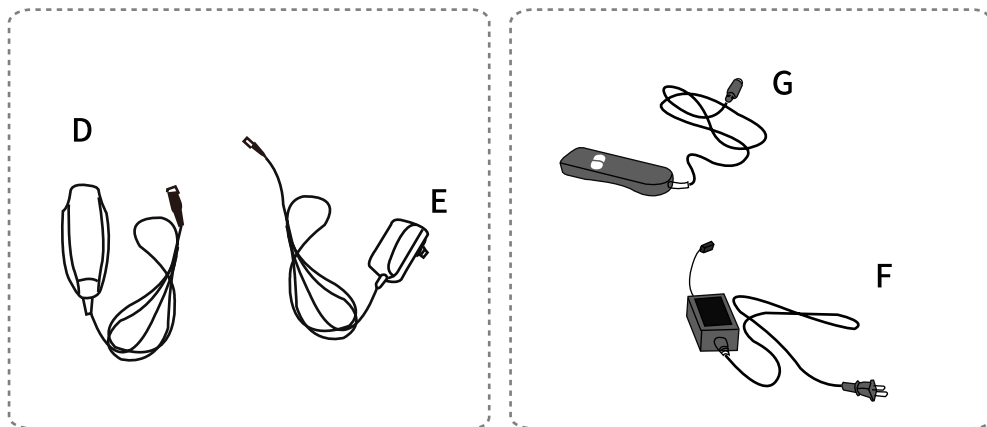

1. Position the armrests


Put the armrests fit into the seat, you will hear a "click" when they are in place.

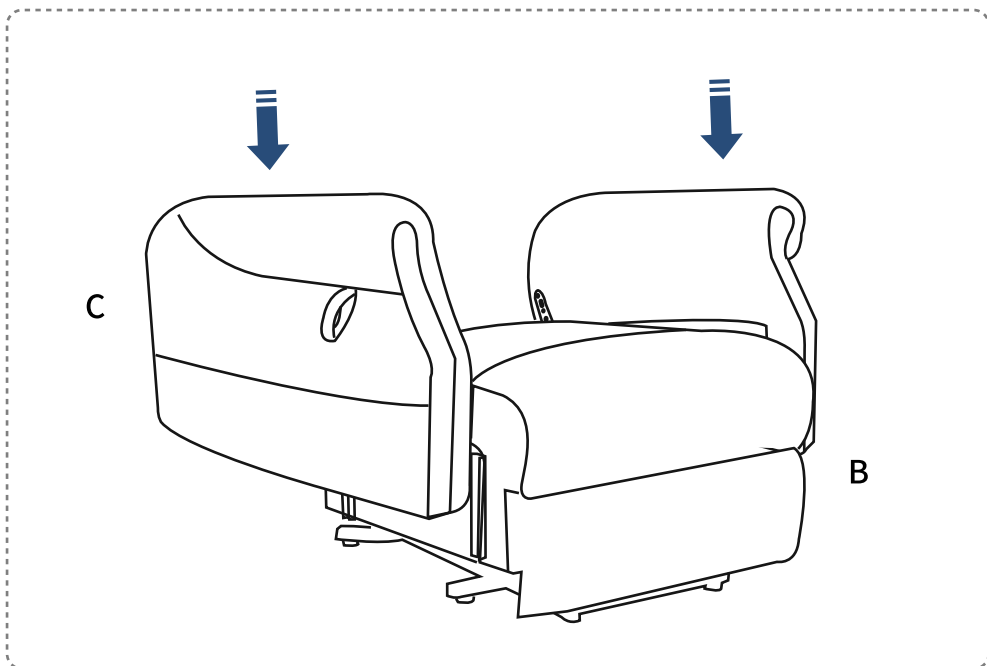






2. Position the backrest

Put the back fits into the seat, you will hear a "click" when it is in place.

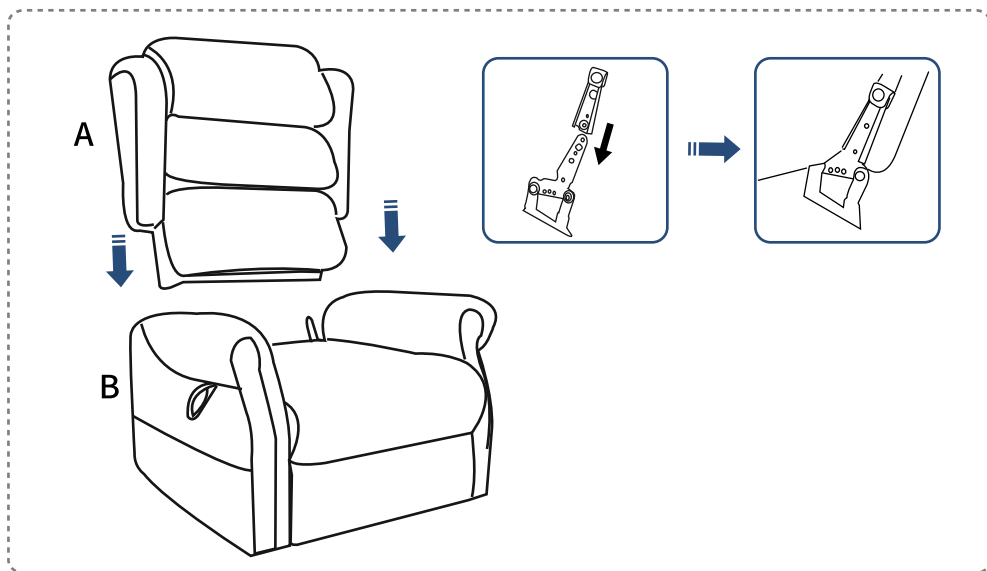



3. Connect the supply cable pins

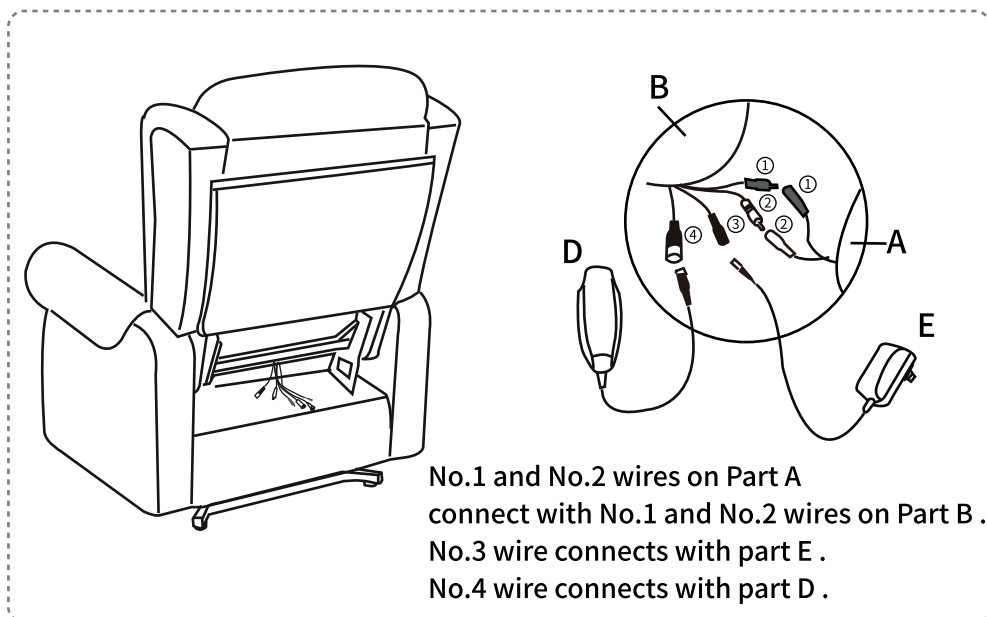

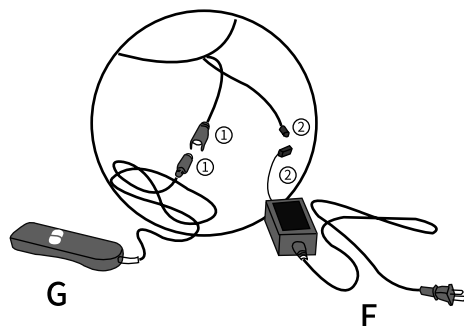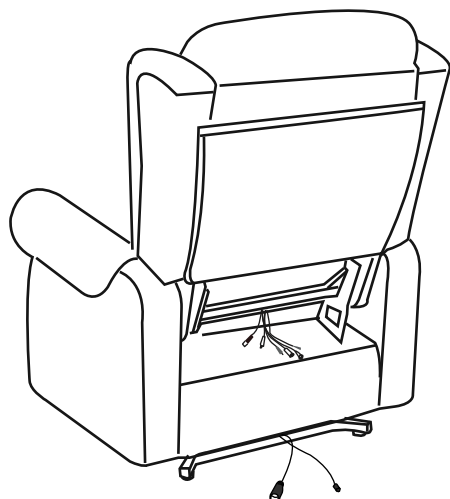


No.1 wires on Part G
connect with No.1 wires on Part A .
No.2wires on Part F
connect with No.2 wires on Part A .

How to help stand?

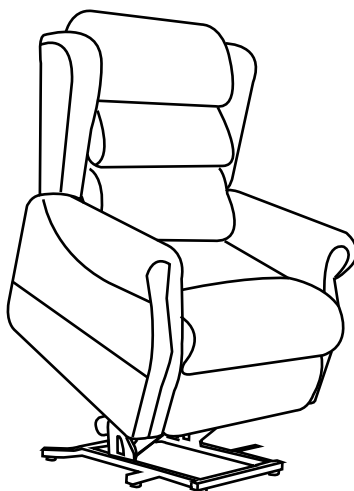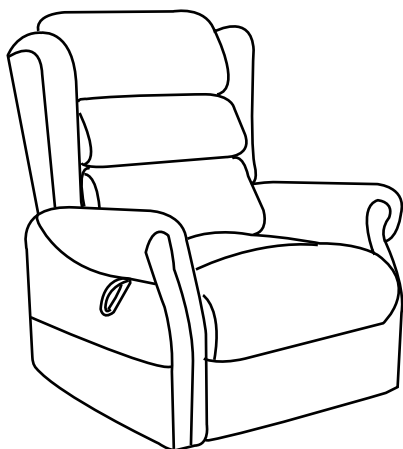




How to recline?

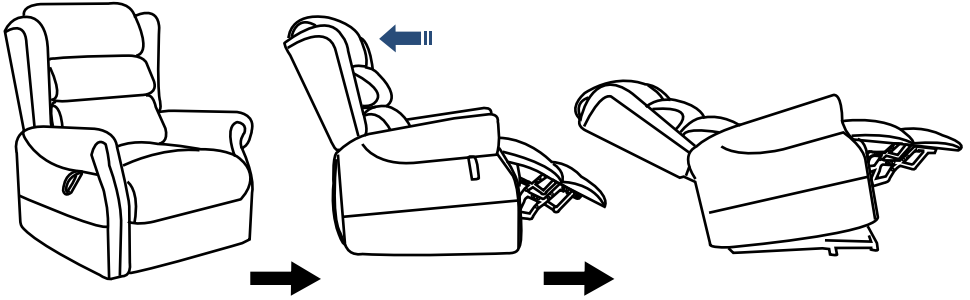

1. Lift the handle up to release footrest.
2. Sit deeper on the chair, put hands on armrests, lie down and use both head and back strength to recline.
3. Have a relax on chair

How to get footrest latch down?

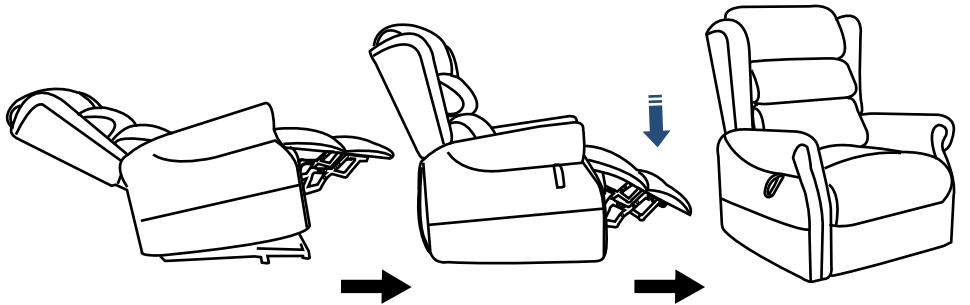

1. Sit on the chair (not recline), backrest will be back.
NOTE : footrest will not close in reclined position.
2. Press the bottom of footrest with some strength (with feet) ,try some times to get footrest latch down.
3. The chair will be in seated position

Data sheet

Packaging size	76cm x 60cm x 58cm 76cm x 58cm x 40cm
Armchair measurements	90cm x 95cm x 106cm
Weight	55 Kg
Maximum weight supported	180 kg
Max Heating Temperature	50-60°C(90°F Max)
Material	FABRIC

Troubleshooting:

If the lift chair is not working, check the following chart. Make sure you turn the power off and unplug from the wall outlet. If you have any questions, please contact us. Do not disassemble or repair any part of the lift chair yourself.

PROBLEM

There is a noise during
the massage
.....

The chair does not function
when the button on the
control unit is activated
.....

POSSIBLE CAUSE

The noise is caused by
the massage motors
and is normal.
.....

Check that the plug is
pushed into the wall
outlet properly.
Check whether the power
supply button is on.
.....

Thank you for using!