

# 2 IN 1 ROLLATOR WALKER INSTRUCTION MANUAL



Scan the QR code to  
view the video tutorial >




## MAINTENANCE

1. The exterior of this rollator walker should be cleaned periodically with a soft damp cloth. Do not use abrasive cleaners or caustic chemicals.
2. Check the parts for any looseness, wear, or damage periodically. If any problems are found, adjust them first before use.
3. Please keep the walker in a cool, dry and safe place when not in use for long periods.

## ⚠ WARNING

1. Before use, ensure the rollator is fully unfolded. Confirm that all screws are tightened and that both the brakes and especially the seat are securely locked in place.
2. Please lock the brakes before sitting down to prevent the walker from moving and causing falls.
3. Never lean forward onto the handles. This can cause the walker to tip forward and result in a fall.
4. Do not recline fully on the backrest, or it may cause the walker to tip over, resulting in a fall and serious injury.

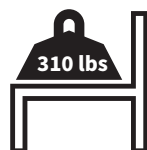

Weight limited  
(140kg/310lbs)

Any problems with the product, please  
contact us through:  
[support@ablecareglobal.com](mailto:support@ablecareglobal.com)

## PARTS LIST

|                                           |                                             |                                                       |                                                     |
|-------------------------------------------|---------------------------------------------|-------------------------------------------------------|-----------------------------------------------------|
| <b>1</b><br><b>Main frame</b><br><br>1 PC | <b>2</b><br><b>Front wheel</b><br><br>2 PCS | <b>3</b><br><b>Rear wheels handrails</b><br><br>2 PCS | <b>4</b><br><b>Cane/umbrella holder</b><br><br>1 PC |
| <b>5</b><br><b>Foot rest</b><br><br>2 PCS | <b>6</b><br><b>Backrest</b><br><br>1 PC     | <b>7</b><br><b>Velcro I</b><br><br>4 PCS              | <b>8</b><br><b>Velcro II</b><br><br>1 PC            |

## ASSEMBLY

### STEP 1

1 ×1



Open the main frame and unfold the armrests on both sides so that the marbles fit into the slots.

### STEP 2

6 ×1




Align both clips of the backrest with the screws and insert, then pull backward to secure firmly.

### STEP 3

3 x2

7 x4


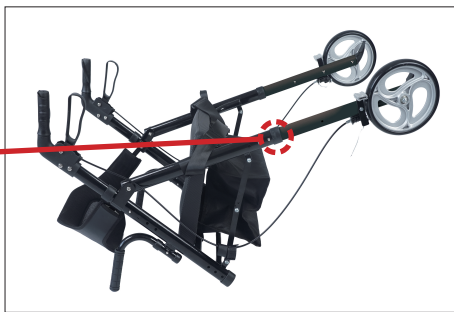


1. Install the handrail tubes and rear wheels into the main frame, ensuring the spring-loaded balls snap into the corresponding holes according to your desired height.

2. Secure the brake cable with velcro.

### STEP 4

2 x2

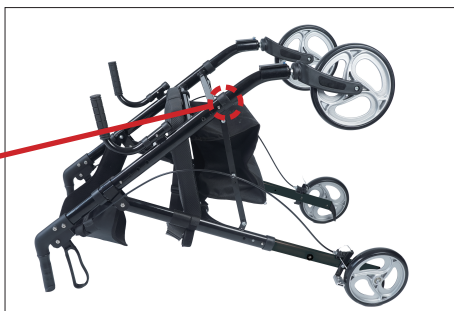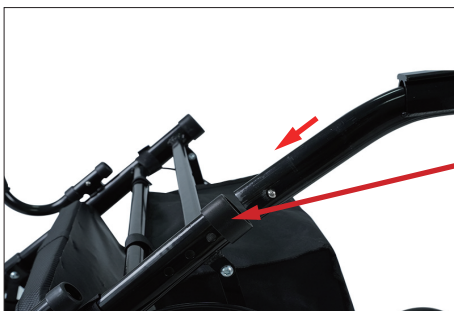

Install the front wheel tubes into the main frame until the spring-loaded balls snap into the corresponding holes.

### STEP 5

5 x2



Install the left and right footrests on the corresponding front wheel tubes, sliding them in from top to bottom and ensuring they are securely in place.

### STEP 6


4 x1

8 x1



1. Install the cane/umbrella holder to the rear wheel tube, secure it from top to bottom.  
2. Secure the cane/umbrella with velcro.

### Brake Adjustment

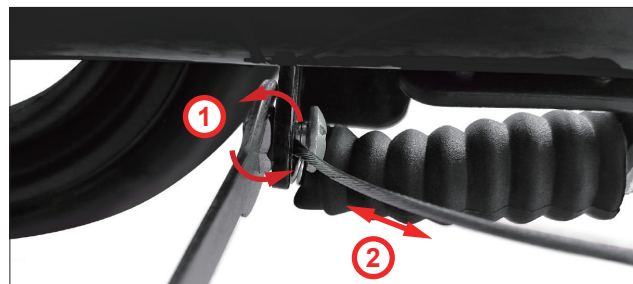

Use a wrench to turn the nut clockwise to loosen it. If the brake is too loose, pull the brake line outward. If it is too tight, pull it inward a little. After adjustment, tighten the nut counterclockwise.



Pull up to slow down

Push down to park