
ADJUSTABLE BATHTUB GRAB BAR

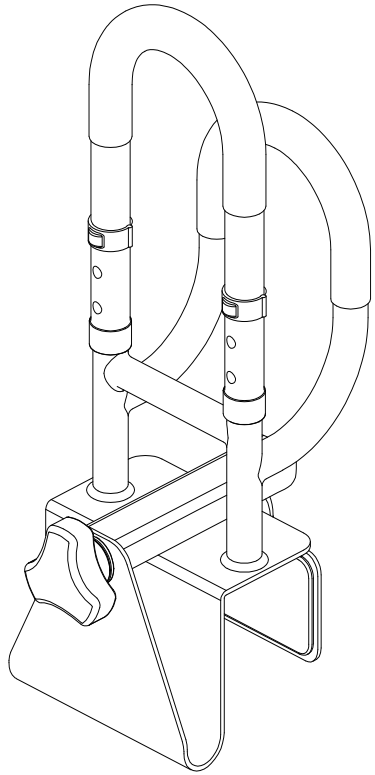


Product Name	Adjustable Bathtub Grab Bar
Model	202603-Bathtub Safe Rail
Main Material	Iron
Overall Dimensions (L×W×H)	20 1/16" x 8 1/4" x 5 7/8" (51 x 21 x 15cm)

ATTENTION

- We will provide our best customer support.
- If you have any questions or issues during installation, please send an email with your telephone number and address.
- We will provide installation support and offer replacement parts for any defective parts.
- Customer satisfaction is crucial for us.
- Customer support email address: **support@ablecareglobal.com**

NOTE

- Before use, verify that the Bathtub Bar is secured by pulling it from side to side and front to back.
- The Bathtub Bar should be checked tightening every time before usage. Check periodically and tighten again when needed.

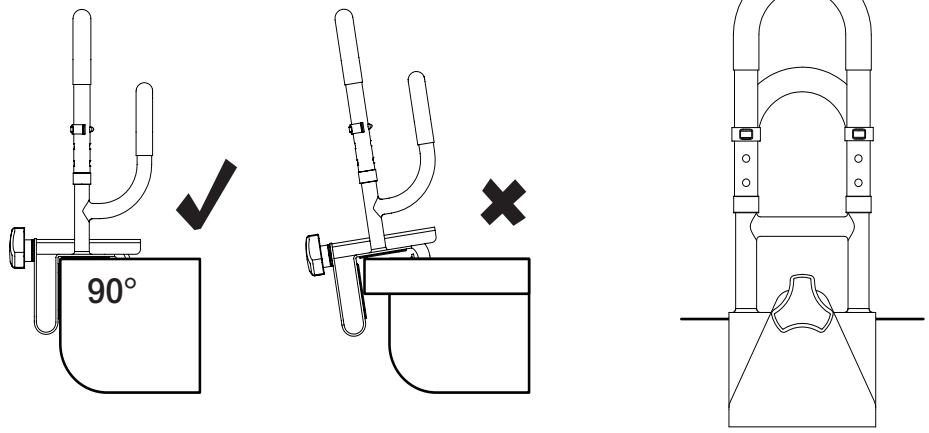


SAFETY INSTRUCTIONS

1. Please follow the instructions in this manual to assemble and use the product.
2. Never apply force over the weight capacity on the product.
3. Keep the parts of this product away from children for safety concerns.
4. Before use, ensure that all parts are connected well and securely.


PARTS LIST

1

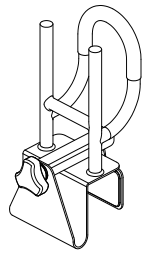


1PC

2

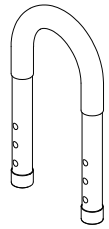

1PC

3

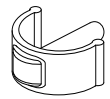


2PCS

ASSEMBLY

STEP 1

1 ×1


2 ×1

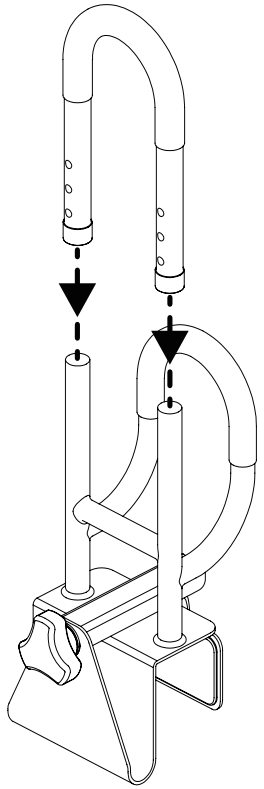

Install the handle on the base as shown.


STEP 2

3 ×1

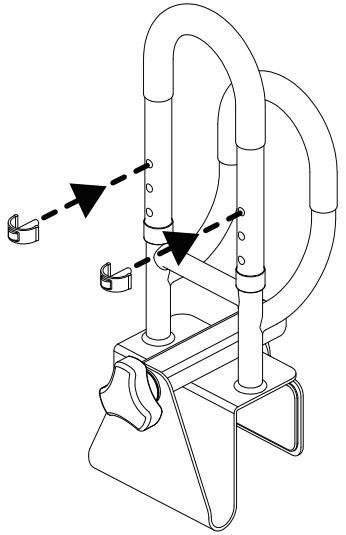

Secure with the two brackets after the high of the handle is decided.

STEP 3

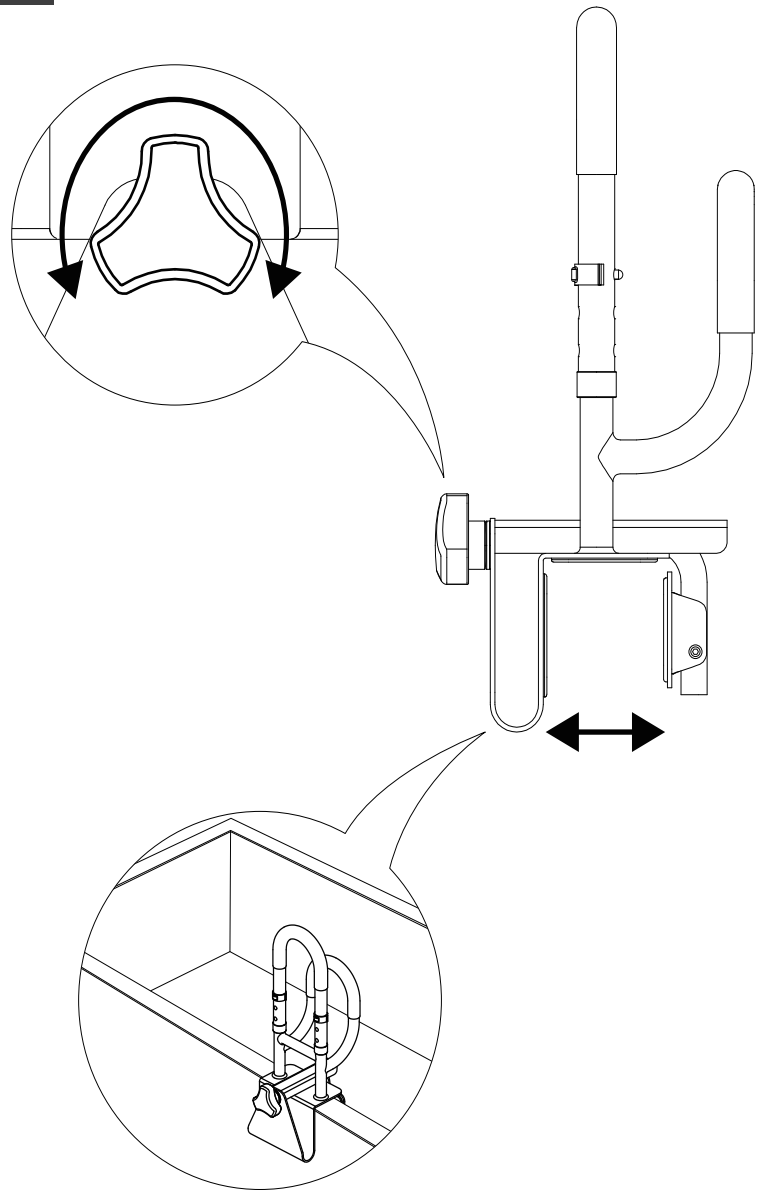

Loosen & fasten the knob to adjust the width of the groove. After the bar is install on the side your bathtub, ensure the knob is fasten enough for securing.