VEVOR

Upgrade · The Home Creator Way


THRESHOLD RAMP

USER MANUAL

VEVOR
Upgrade·The Home Creator Way

THRESHOLD RAMP

Model: LD9YV0 /LYIMV0 /Z2UAV0 /GRBRV0 /OVUSV0 /SNPTV0

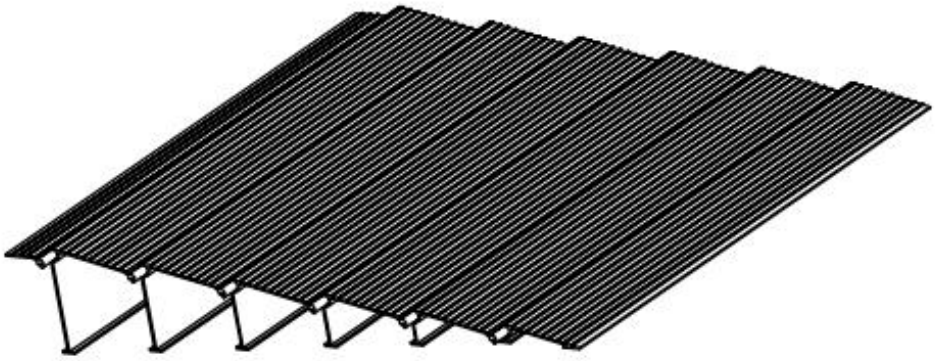This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

1. Safety Instructions

WARNING:

Carefully read and understand instructions and procedures before installing and using the threshold ramp.

Do not modify the threshold ramp design.

- Do not move a threshold ramp while standing on it.
- Do not use it in bad weather, such as wind, rain, snow, and so on. It is recommended to use the ladder in clear weather.

1. During installation, please wear gloves to prevent hand injuries.
2. Do not stack items or allow children to play on the ramp.
3. Regularly check the structure and fasteners of the ramp to ensure its stability and safety.
4. If the surface of the ramp is worn or damaged, please replace or repair it in time.

5 The product complies with relevant standards and regulations. Installation of this device must comply with all applicable regulations, and it should only be operated in well-ventilated areas. Please read the user manual before installing and using this equipment.

6. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

7. The method of fixing stated is not to depend on the use of adhesives since they are not considered to be a reliable fixing means.

2. Maintenance Guide

1. Place the product on a sturdy, secure, accessible, and level starting point. The minimum weight capacity of the starting point must be equal to that of the ramp (800 lbs). Do not exceed the maximum weight capacity of the product.

2. Ensure that the product is securely flush against the starting point. If a safe, accessible, and level starting point cannot be found, or if the product does not fit securely against a higher starting point, do not use the ramp.

3.Do not tamper with, attempt to repair, or modify any part of the product. Exercise caution to avoid pinch points during assembly and handling.
4. If using the product with other mobility equipment, refer to the user guide of the equipment to understand the acceptable usage guidelines, including appropriate slope, incline, and chair direction. Do not exceed the recommendations of the mobility equipment.
5. Before each use, carefully inspect the fasteners on the ramp to ensure they are not loose. Tighten any loose fasteners if necessary.
6. When using the ramp, ensure it is securely flush against the starting point to ensure safety.
7. Before each use, clean the ramp surface of ice, snow, dust, leaves, and other debris. Remove accumulated ice and snow, and sweep the surface clean. For dust and debris, use a soft bristle brush and non-alkaline detergent to clean the ramp surface, then rinse thoroughly.
8. The ramp surface may become slippery in wet conditions. Walk with caution and pay attention to your footing to prevent slipping.
9. If the ramp will not be used for an extended period, store it in a dry and well-ventilated area away from direct sunlight.

SAVE THESE INSTRUCTIONS

3. Product Description and Main Features

The threshold interlocking aluminum alloy ramp is made of high-quality aluminum alloy extrusion. This product is durable and provides convenience for people. It is lightweight, portable, stable, and durable and features an interlocking design for easy installation and storage.

- Main Features:

1. Material: High-strength aluminum alloy, extruded.
2. Surface: Anti-slip treatment, deburring.
3. Interlocking design: Each section’s interface is precisely designed; align each section's interface and push them together for installation.
4. Load limit: Up to 800 lbs.
5. Height optional: suitable for thresholds of different heights.

4. Product Parameter

Item	Model	Picture	Height	Load Limit	N.W	G.W
1	LD9YV0		1inch	800lbs	1.5kg	1.9kg
2	LYIMV0		2inch	800lbs	2.7kg	3.2kg
3	Z2UAV0		3inch	800lbs	4.0kg	4.5kg
4	GRBRV0		4inch	800lbs	5.45kg	6.0kg
5	OVUSV0		5inch	800lbs	7.05kg	7.95g
6	SNPTV0		6inch	800lbs	8.6kg	9.56kg

5. Installation Guide

(1) Before starting the installation, please wear gloves to protect your hands.

(2) Choose an appropriate installation location, ensuring that the ground is level and firm.

(3) Align and interlock each section of the ramp together until all parts are even.

(4) Starting with the 1” extension, install the threshold extensions in the order shown (FIG. 1) until the final height is reached. Bring the threshold extensions together at an approximate 45° angle, then rotate them into place after the bulb is engaged in the slot.

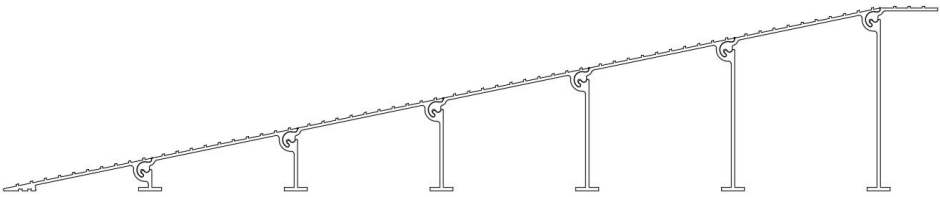
FIG.1



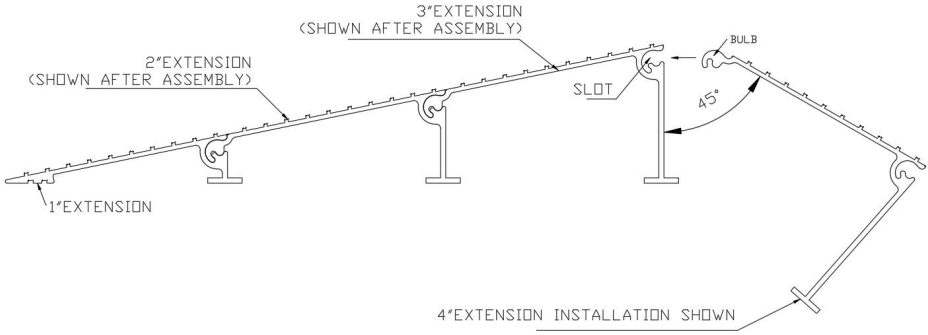
FIG.1

(5)Use the provided screws to tighten each section of the ramp to ensure that the ramp does not slip out or make noise during use.

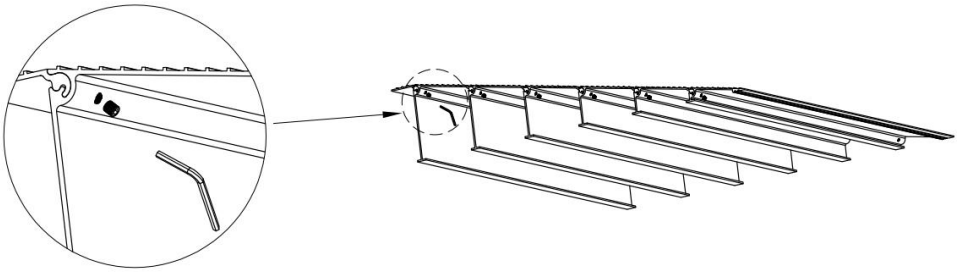(5) After installation, check whether the ramp is stable and tightly attached to the ground.

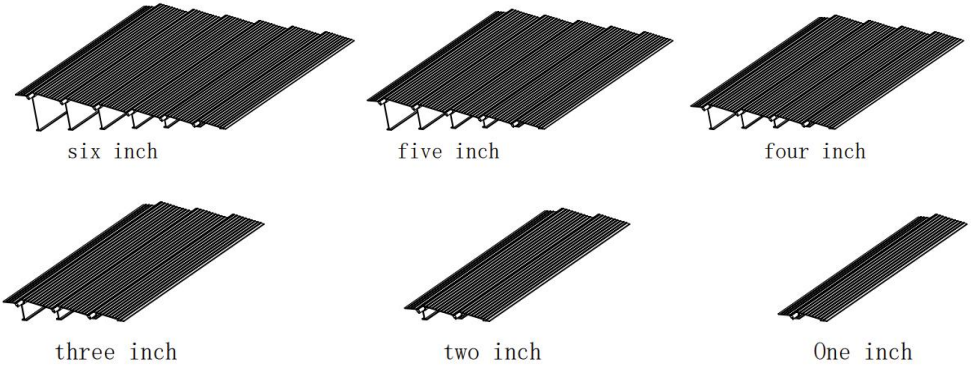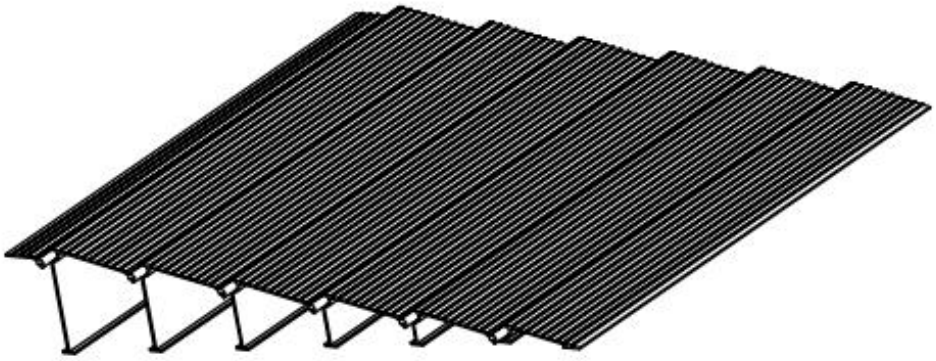
(6) When using the ramp, push the equipment up and down the ramp steadily, avoiding sharp turns or speeding.

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