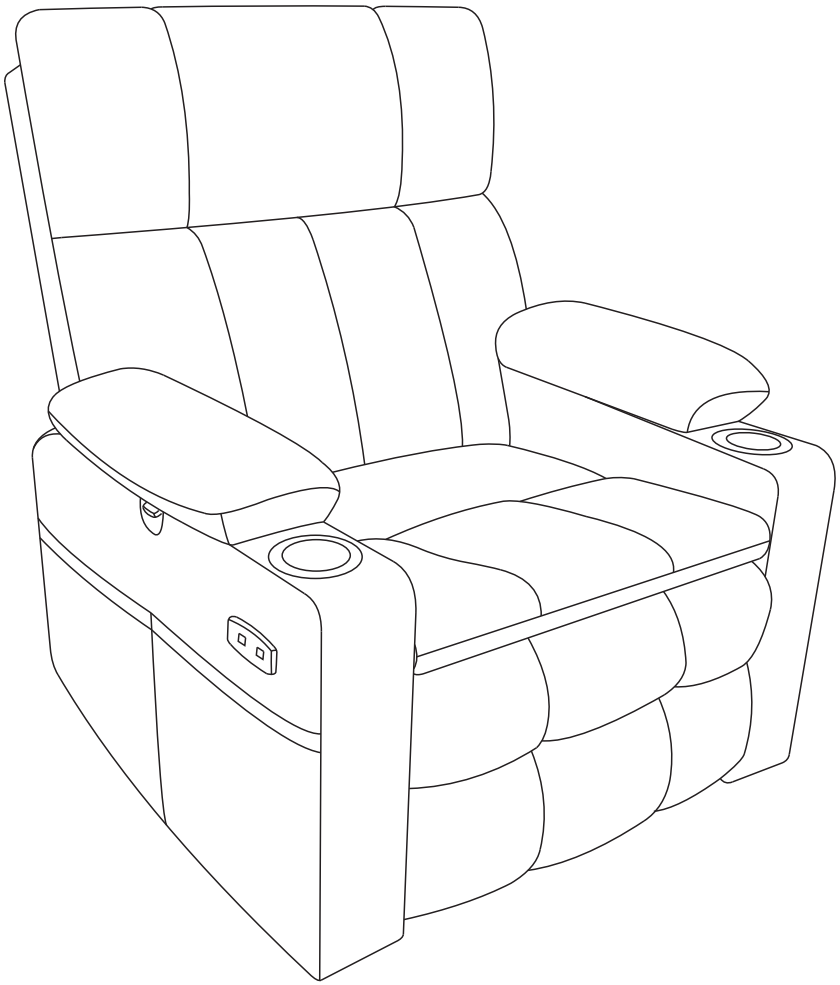
MANUAL

Power Lift Recliner with Adjustable Armrest



1. ASSEMBLE GUIDE

Thank you very much for choosing our product!

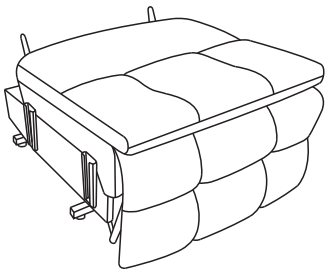
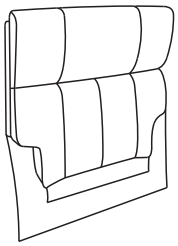
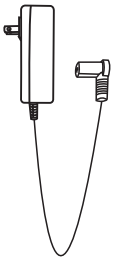
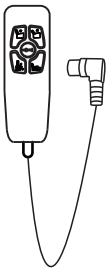
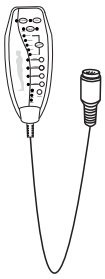
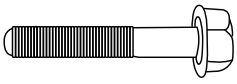
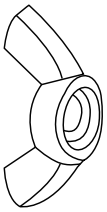
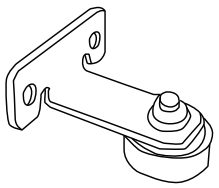
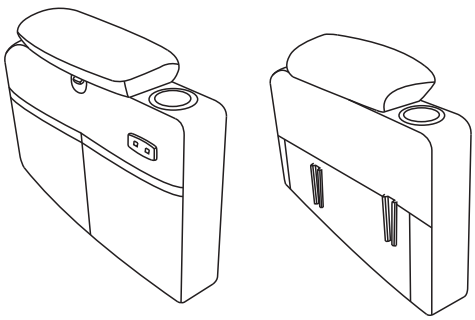
Please take time to read these instructions thoroughly and follow each step carefully for safe and easy operation. Keep this manual for future reference, especially for the safety information.

Also, it is better for you to assemble the chair in cooperation with 2-person, having a second person to help you makes it even easier.

Please read and understand this manual before any assembly. Before beginning assembly of product, make sure all parts are present. Compare parts with packaging contents list.

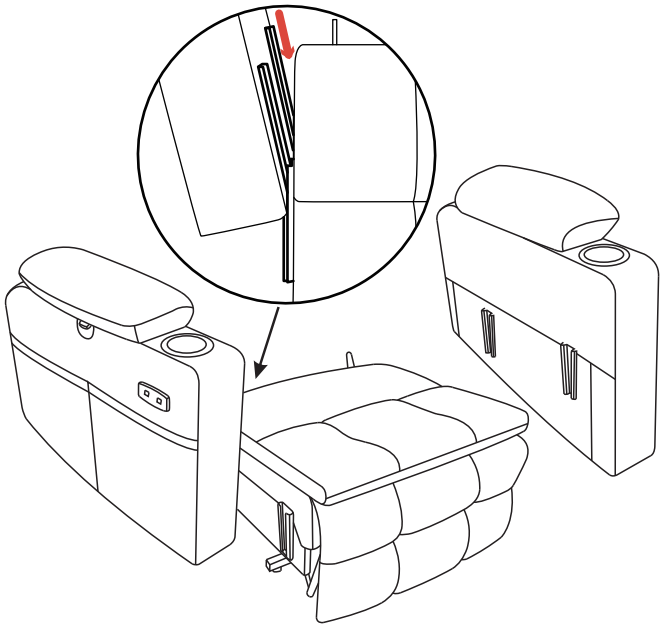
If any part is missing, damaged and any questions, please feel free to contact us.

2. PARTS LIST

BOX 1				
Part	**A**			
Description				
Qty	1PC			
BOX 2				
Part	**B**	**C**	**D**	**E**
Description				
Qty	1PC	1PC	1PC	1PC
Part	**G**	**H**	**I**	
Description				
Qty	4PCS	4PCS	2PCS	
BOX 3				
Part	**F**			
Description				
Qty	2PCS			

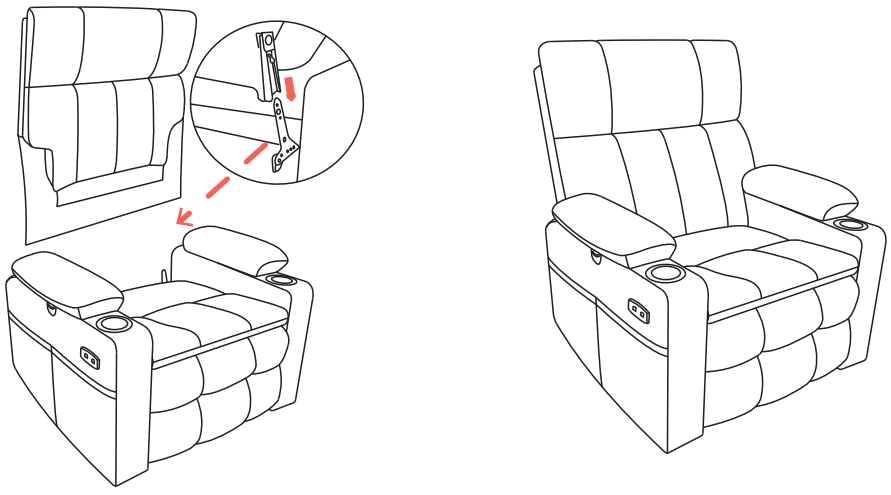
3. ASSEMBLY STEPS

Step1:



Align the two side slots and install the Armrest.

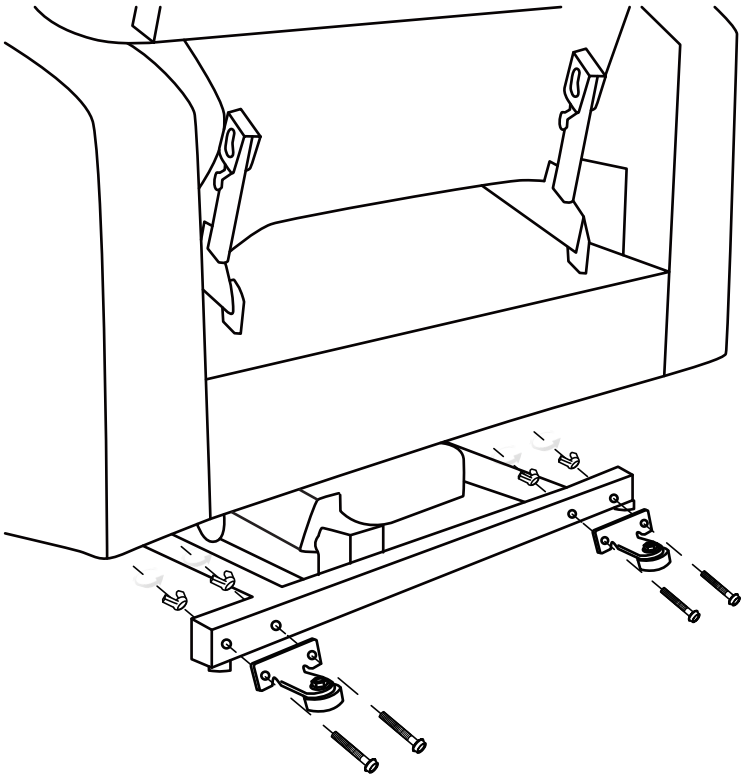
Step2:

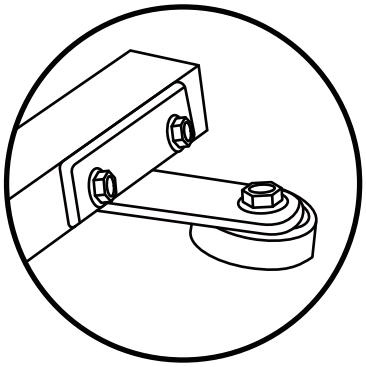


Align the back to the seat base through the slots thento press down backrest, until hear a "click",and then ensure that the backrest is stable and will not be pulled out.

3. ASSEMBLY STEPS

Step3:

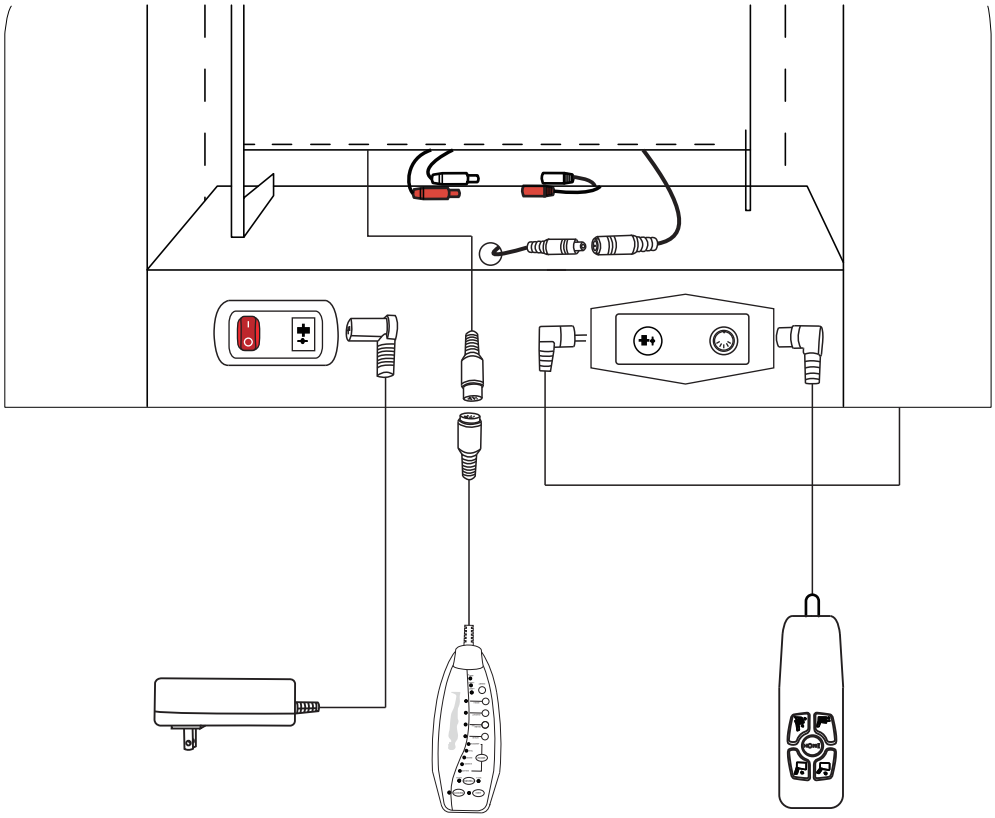




According to the legend, match the screws on the base one by one and install the anti tipping parts

3. ASSEMBLY STEPS

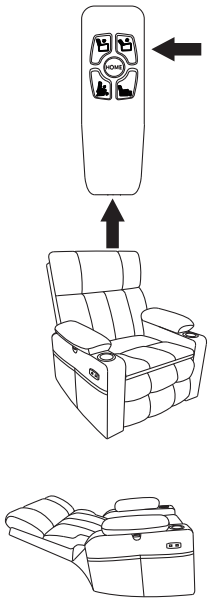
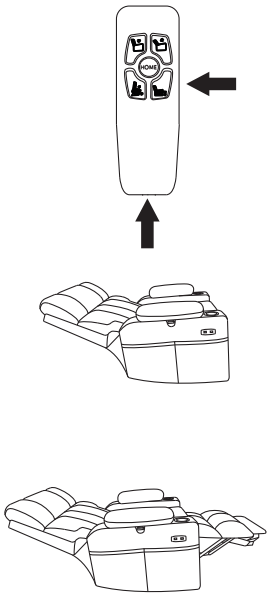
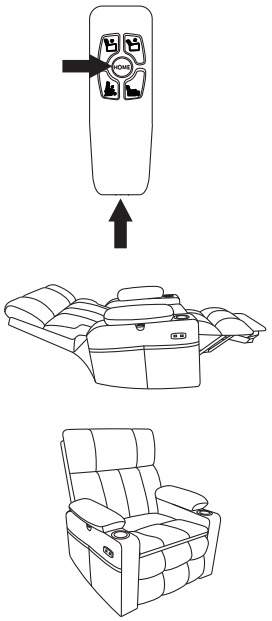
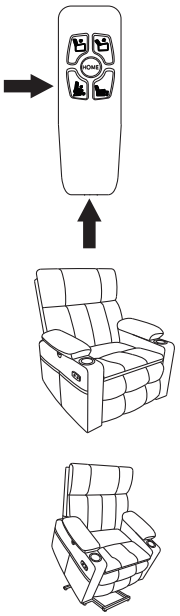
Step4:



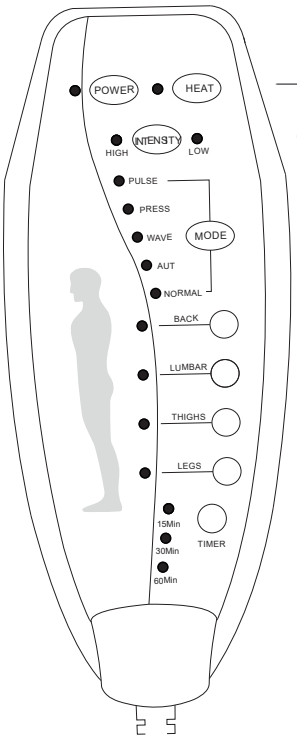
Connect the power cables one by one as shown in the diagram.

Please note:
The USB socket in the bottom of the recliner armrest.

4. OPERATING INSTRUTIONS

STEP 1	STEP 2	STEP 3	STEP 4
			

How to use the remote control for massage?



POWER: Turn on/off remote control.

HEAT: Turn on/off heat function for lumbar part only.

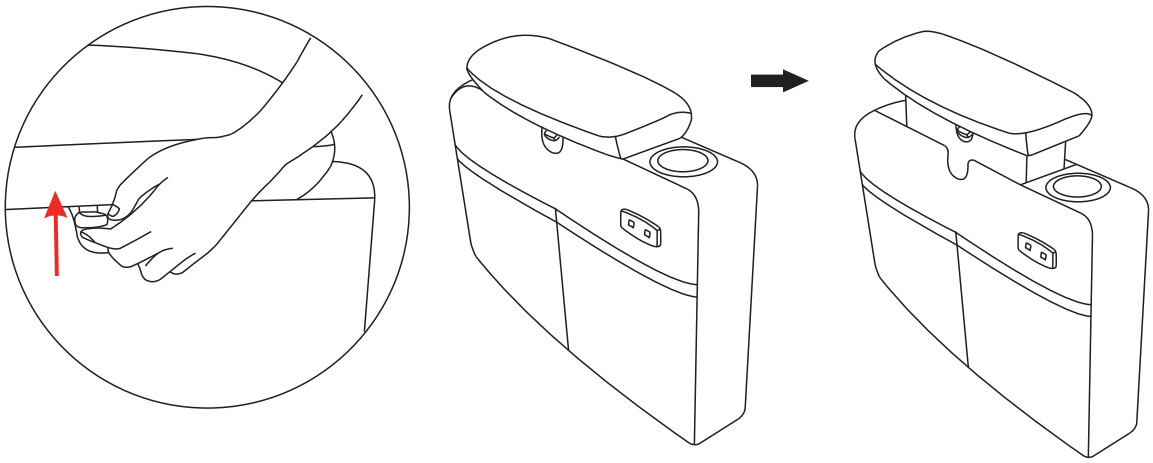
ITENSIITY: Control massage function from low to high.

MODE: Choose one massage mode you like(pulse, press,wave,auto, normal).

You can choose 1-4 massage parts(back, lumbar, thighs, legs)to enjoy.

You can choose 3 time modes as you like: 15min,30min,60min.

4. OPERATING INSTRUTIONS



Press and hold the button on the armrest to control its lifting and lowering.

5. TROUBLESHOOTING

1.What should I do if my shipment is damaged or products received are defective?

We're sorry to hear that! Please contact us and provide adetailed description of the damage, as well as photos (or videos, where appropriate) showing the damage.
We will do our best to resolve the issue as soon as possible.

2.Can replacement parts be ordered?

Just write to our customer service team and let us know your order lD.
and we will be happy to provide further assistance.

3.Who do I contact for repair queries?

Easy! Just email us and provide a detailed description of the fault
as well as photos.

4. How do I claim my warranty?

You can write to us to specify your inquiry, and we will be able to assist you further.

6.What should I do if I want to return the product?

We offer 30-day returnn policy upon product received. To initiate a return request,please message us and with your reason of returning so we canfacilitate the request and assist you further. For more details on the return policy,please refer to Return section.

6. PRODUCT CARE

• Dust regularly with a soft,dry cloth to prevent soil build-up.Keep away from heat to avoid glazing,melting or scorching.Use a professional cleaning service to clean the fabric.

7. WARNING

• This product is not intended for institutional or commercial use.
• FALL HAZARD - To prevent death or serious injury, always keep child within arm's reach.
• Avoid placing this item in direct sunlight to prevent fading.
• Check this product for damaged hardware, loose joints, missing parts or sharp edges before and after assembly and frequently during use.

✱ Customer Service

Any Problem Please Contact Us

• If you have any questions, please provide the order number and related pictures to show the problem through email.
We will reply to you within 24 hours, and solve the problem as soon as possible.