


Massage Seat Cushion

Instruction Manual

AMD-Z01-GY

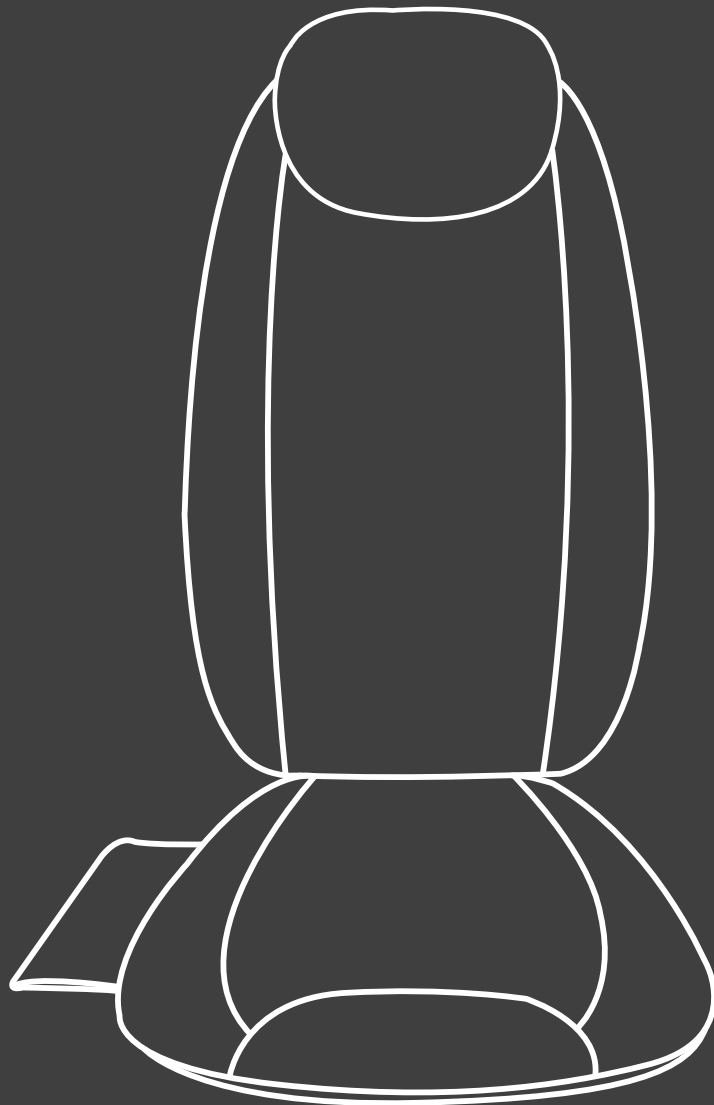

✉ HJ-Customer-Service@outlook.com



☎ (909)637-7665 | (909)637-7593

EST 9:00AM--17:00PM

Congratulations on your purchase of the Massage Seat Cushion. This unit combines soothing heat therapy with four synchronized shiatsu massage nodes to deliver deep, penetrating relief—helping your body feel relaxed and revitalized whether at home or on the go.

CONTENTS

- Massage Seat Cushion
- 12V Home Adapter
- 12V Car Adapter

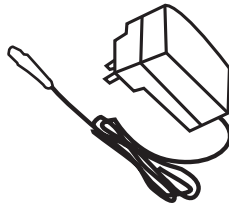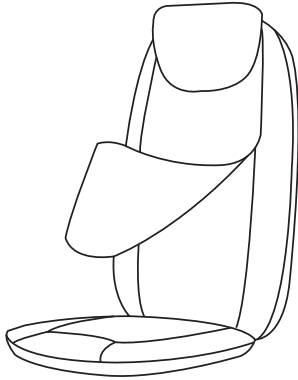

TECHNICAL SPECIFICATIONS

Package Dimension: 26.8*5.9*16.5IN

Model: AMD-Z01-GY

Weight: 6.2LBS

Voltage: 12V 2.5A

Nominal Power: 30W


Preset Operating Time: 5min/10min/15min

THIS DEVICE COMPLIES WITH PART 15 OF THE FCC RULES. OPERATION IS SUBJECT TO THE FOLLOWING TWO CONDITIONS: (1) THIS DEVICE MAY NOT CAUSE HARMFUL INTERFERENCE, AND (2) THIS DEVICE MUST ACCEPT ANY INTERFERENCE RECEIVED, INCLUDING INTERFERENCE THAT MAY CAUSE UNDESIRE

OPERATION. Please read all the safety instructions and warnings in this manual before using this product.

SETUP AND OPERATION

1. Secure the massage seat cushion to a chair using the elastic straps, or place it on another suitable support.
2. Connect the adapter cable to the corresponding input cable on the cushion.
3. Plug the home adapter into a wall outlet, or plug the car adapter into your vehicle's power socket.
4. Turn on the unit using the controller (see instructions on page 4).
5. After use, place the controller into the storage pouch located on the side of the cushion.
6. You may soften the massage pressure by placing a towel between yourself and the unit.

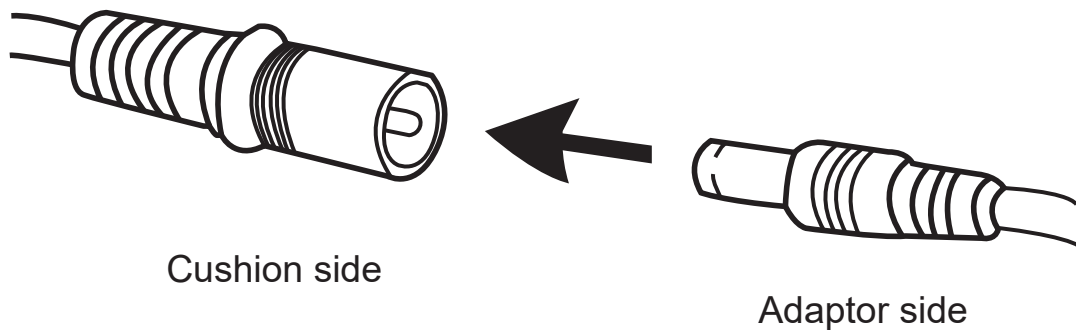

- **Unplug the adapter from the outlet when not in use.**
- **Do not set up or use the massage seat cushion in a bathroom or any other wet / damp areas**
- **Do not use this unit while driving.**

SAFETY INSTRUCTIONS



Please read the following instructions carefully before using your massage seat cushion to ensure trouble-free operation and optimal efficiency. Keep these operating instructions for future reference.

1. The massage seat cushion complies with recognized technical standards and the latest safety regulations.
2. Do not wet the unit, do not use pins, and never remove the cover.
3. This unit is NOT A TOY. Close supervision is required when this appliance is used by, on, or near children or persons with disabilities.
4. This unit should never be left unattended when plugged in.
5. Only use the original power adapter(s) provided with this unit (12V home adapter or 12V car adapter).
6. Do not use any substitute or third-party adapter.
7. Any repairs must only be performed by authorized expert personnel. Improper use and unauthorized repairs are prohibited for safety reasons and will void the warranty.
8. Never touch the power plug with wet hands.
9. Avoid exposing the unit to water, high temperatures, or direct sunlight.
10. Do not use damaged cables, plugs, or loose sockets.
11. Never operate the unit under blankets or any conditions that may block the air vents.
12. If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative, or qualified personnel.
13. In case of malfunction, disconnect the device from the mains immediately.
14. Do not use this unit if you have skin disorders, open wounds, swollen or inflamed areas.
15. Misuse or incorrect use excludes any liability for damage.
16. Do not use this unit while driving.
17. Do not use the unit while sleeping.
18. To prevent overstimulation of muscles and nerves, do not exceed the recommended massage duration of 15 minutes per session.
19. Massage—including hand massage—is contraindicated during pregnancy or if any of the following conditions are present in the massage area: recent injuries, thrombotic diseases, inflammation, swelling, or cancer. Consult a doctor before using the massage seat cushion for treatment of any ailments or illnesses.
20. If you rely on electrical medical devices such as pacemakers, consult your doctor before using this massage seat cushion.
21. The packaging materials supplied are not toys and should be kept away from children.

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

CONTROLLER INSTRUCTIONS


NOTE:

The massage seat cushion is equipped with massage rollers that move to provide relief and comfort.

When the Power button is pressed for the first time, the massage roller will automatically move to its lowest position before starting the massage program.

The massage roller will always return to this lowest position after the device has been turned off.

If the electrical supply is interrupted, the massage roller will revert to the lowest position when power is restored.

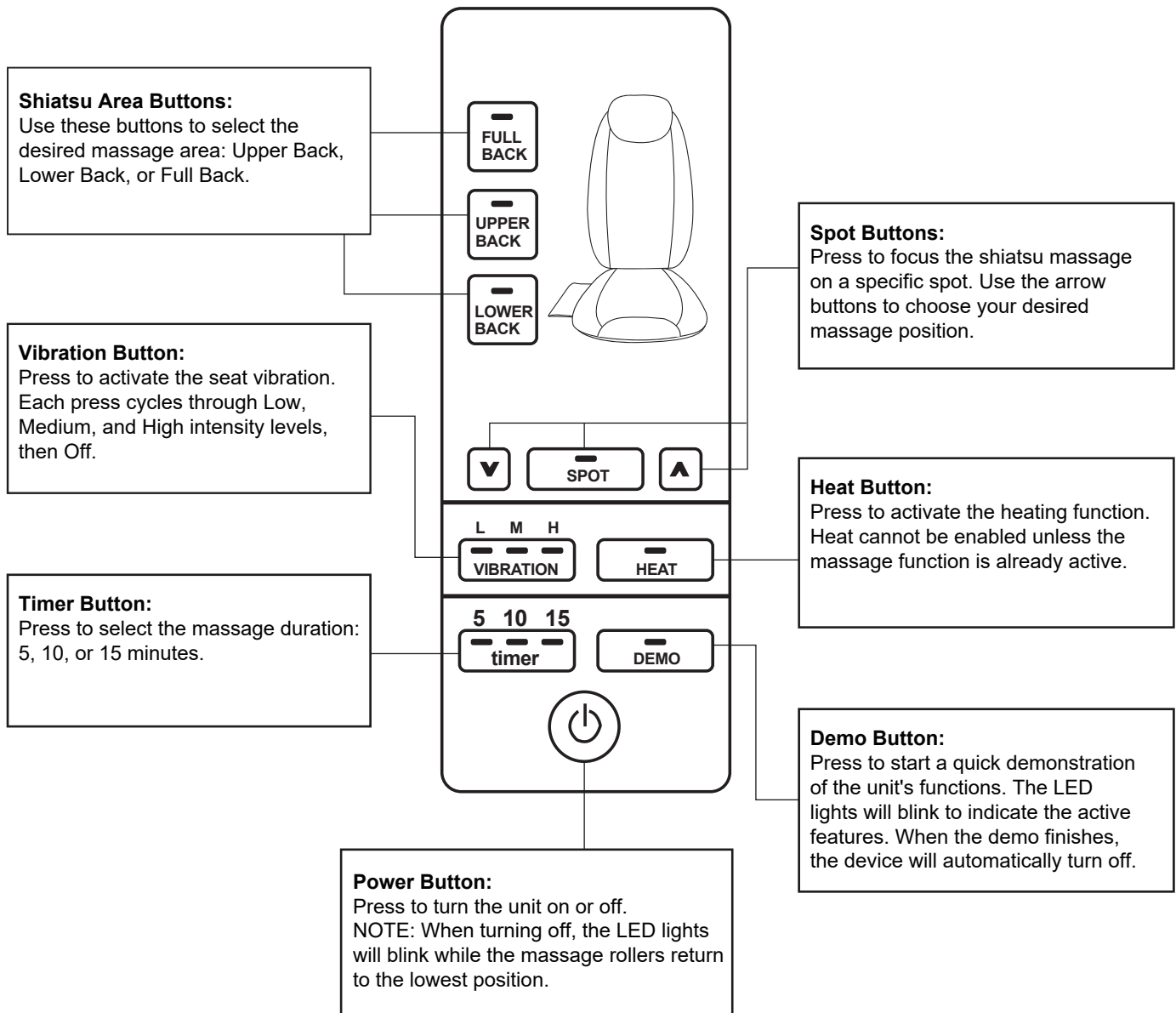


CARE AND MAINTENANCE

1. Do not wash the product or the adapter in a washing machine or immerse them in water.
2. Always switch off the unit and unplug the adapter from the socket with clean, dry hands before cleaning.
3. Wipe the massage seat cushion with a slightly damp cloth and dry it thoroughly. Do not use alcohol, thinner, bleach, or other harsh chemical.
4. Do not overstretch the elastic back straps.
5. Keep the product away from heat sources, open flames, and direct sunlight.
6. Do not apply excessive pressure or impact to the seat cushion or massage rollers. Avoid prolonged or excessive use.
7. Store the massage seat cushion in a cool, dry place, out of reach of children.
8. Avoid contact with sharp edges or pointed objects that may cut or puncture the cushion.

TROUBLESHOOTING

Problem	Solution
Massage heads rotating at reduced speed.	The massage heads may be overloaded. Reduce the pressure and try again.
Massage heads only rotate in the upper or lower region.	The unit may be set to upper or lower back mode. Press the 'Full Back' button to activate both regions.
Unit switches off suddenly.	The unit may have been used continuously for 15 minutes and shut off automatically for safety. Let it rest for 15–30 minutes before restarting.
Unit fails to start.	Ensure the adapter is properly plugged into a working socket. Check the adapter output is connected to the unit's power inlet. If the unit was used for over 15 minutes, allow it to cool down for 15–30 minutes before restarting.
Jerky motion of massage heads / Vibrations from the unit.	Excessive pressure may have been applied onto the unit. Reduce the load and try again.
Infrared lights are on but no motion from massage heads.	The unit may be jammed. Remove the load and press the 'Full Back' button to resume operation.

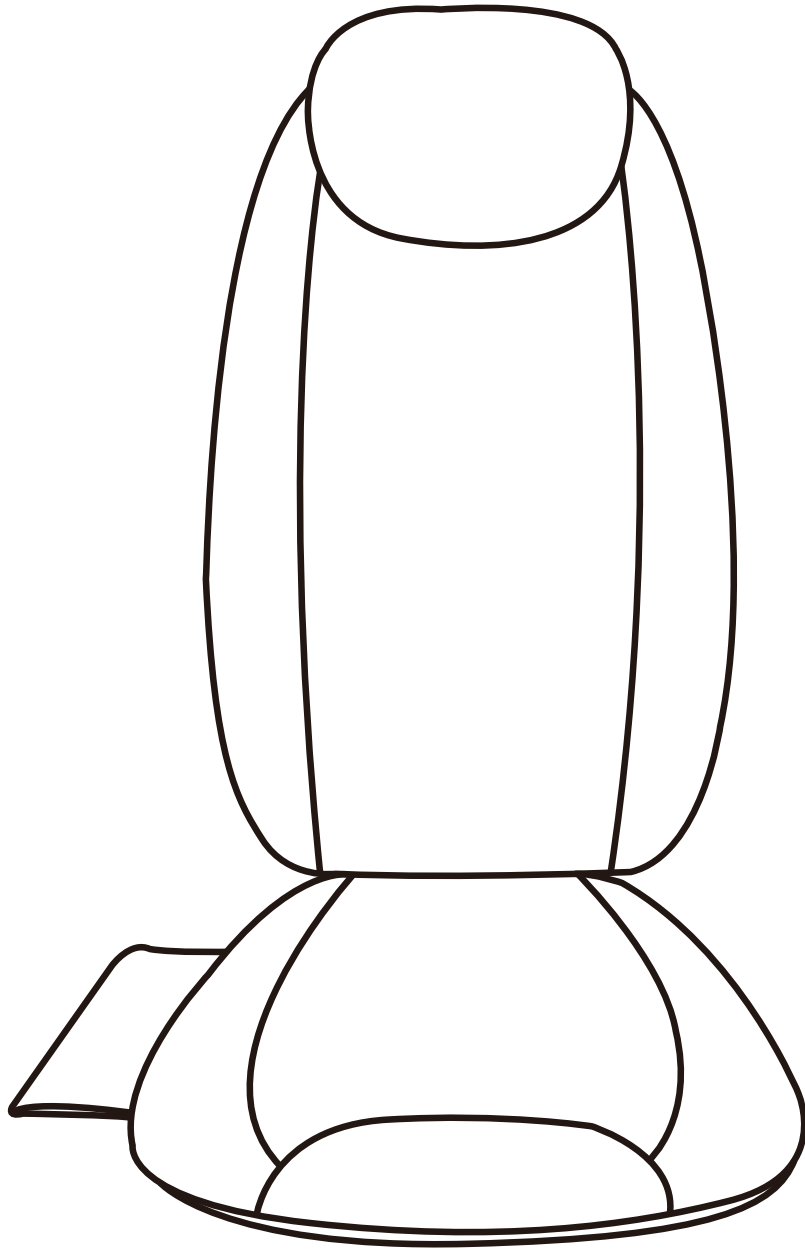

Technical Support

(909)637-7665 | (909)637-7593

Other Info

HJ-Customer-Service@outlook.com

Thank you for choosing this product. We strive to provide you with the best quality products and services in the industry. Should you have any issues, please don't hesitate to contact us.