



3-Section Memory Foam Massage Table TS-3333




CE

USER'S MANUAL


Precautionary Instructions

1. Maintain the table periodically. Make sure all screws are secure and the cable has not deteriorated or frayed. Make sure there are no cracks on the wooden structure.
2. Do not overload the table beyond the recommended weight limit. This could result in serious injury to the user. This product is tested to hold a max working load of 551 lbs (250kgs). Don't load weight over 44 lbs (20kgs) on the headrest cradle and around position of face hole.
3. Be cautious not to tangle the cables with legs or support beams.
4. Do not use alcohol based cleaners when cleaning the vinyl. This could cause damage to the vinyl material. For PU leather, use dry soft cloth to clean the leather or use PU leather friendly cleaning agent to clean and dry it quickly.
5. Make sure when using the table that all four legs are at the same level. Using the legs at different levels can cause the table to be unstable and cause injury.
6. Do not allow children to assemble or disassemble the massage table for you unless under supervision.
7. Be careful when using sharp objects as they can pierce the covering. Keep the table's covering away from being contacted by the sharp items.

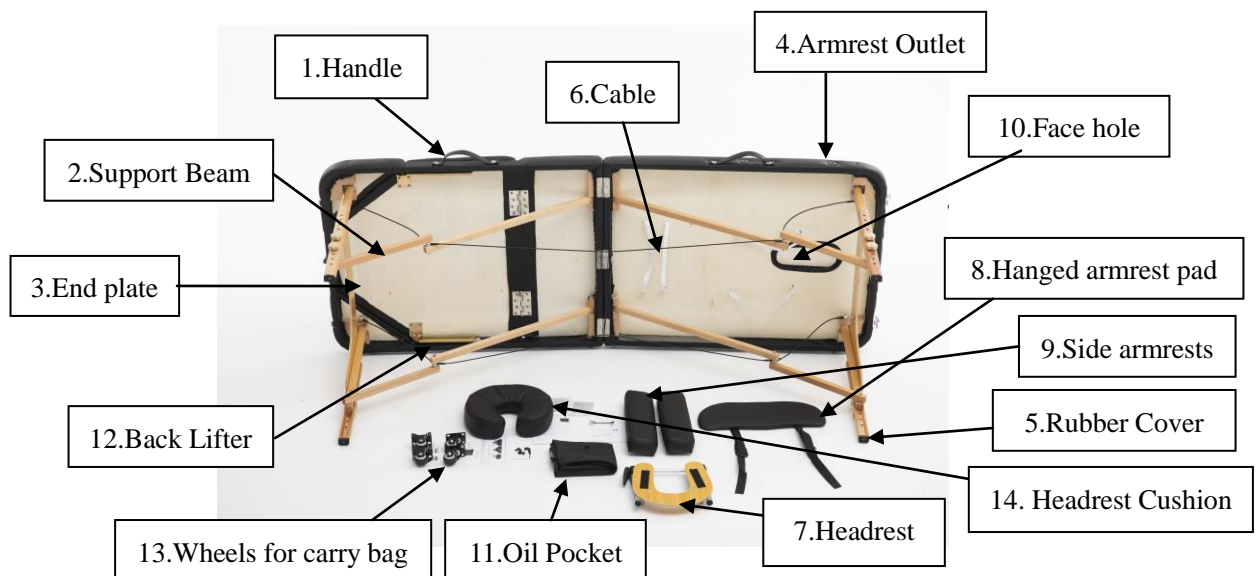

- | | | | | |
|----------------|-----------------|--------------------------|----------------------|-----------------|
| 1. Handle | 2. Support Beam | 3. End plate | 4. Armrest Outlet | 5. Rubber Cover |
| 6. Cable | 7. Headrest | 8. Hanged armrest pad | 9. Side armrests | 10. Face hole |
| 11. Oil Pocket | 12. Back Lifter | 13. Wheels for carry bag | 14. Headrest Cushion | |

Setting Up Your Table

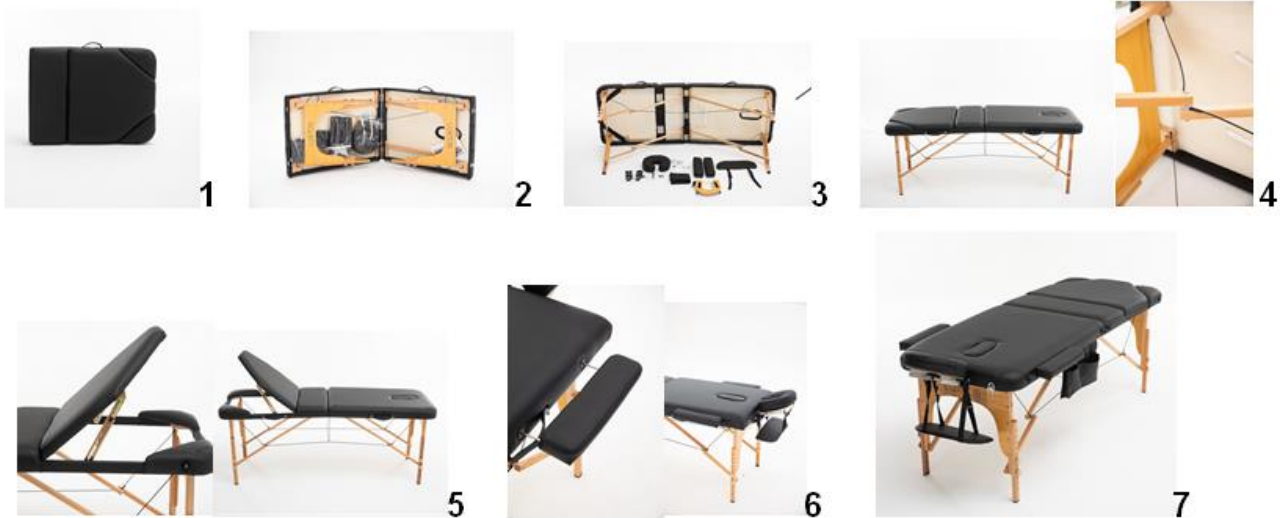

1. Release two buckles, while keep the buckle of back locked.
2. Lay the table on its side with four rubber feet on the floor, open up the table with four legs extend as full as possible.
3. Take out all accessories.
4. Lift the table up until all legs are on the floor, to ensure the cables have all been fully extended, be cautious not to tangle the cables with the legs.
5. If you'd lift up the backrest, release the mid-buckle and lift up the back, you will hear "click" sound in each level of angle, stop at the angle you need, to ensure both sides of backrest are in the same level, no sloping.
6. Put on the side armrests, put on the hanged armrest on the adjustable headrest cradle, and then put on the headrest cushion. Release the side handle of headrest cradle and press it in once the headrest is adjusted to the required position. Remarks: If the headrest is too loose to be locked, press the opposite side of side handle and rotate the handle in clockwise direction to make it firmer and press it in.
7. Now the table is ready for your use. You can also put on table cover if you need.

Re: The height of table is adjustable by turning off the knobs on the leg and take out the movable leg to adjust the height by choosing the holes on the movable leg, put back onto the leg and fasten the knobs firmly. Four legs should be done in the same way one by one and make sure all legs are adjusted into same height level. Don't overload weight of 44 lbs (20kgs) on the position around breast holes.

Remarks: There are 4 wheels for carry bag to be assembled by DIY, follow the guide of assembly in the separate page.

Instructions for Installing Wheels of Carry Bag



Parts Included: **a.** 2 pcs of universal wheels **c.** 9 pcs of screws
 b. 2 pcs of directional wheels **d.** 1 pc of wrench

Assemble all wheels on the bottom of the carry bag. Both Universal wheels should be installed on the same side of the carry bag that has the handle.



To install wheels, place the wheels on the bottom of the carry bag and align the screw with the hole in the bottom of the carry bag. Tighten the screw with a wrench.

Note: Each wheel is to be tightened with two screws.



Note: Avoid using the wheeled carry bag on roads with tough conditions or uneven surfaces

Fold Up and Carry Along



1. Take out all accessories.
2. Lay the table on its side with four rubber feet on the ground.
3. Put all accessories under the elastic strap under the top veneer.
4. Fold the legs on each end at angle towards inside by pressing the corresponding support beams.
5. Fold the table up and lock the buckle.
6. Put table inside the carry bag and zipper up. Carry the table with carry bag with shoulder belt.

Warning: If the frame/legs/table structure is damaged in anyway, then you should not use this table. Using this table when damaged may cause harm or serious injury to you or your patient. You should seek repair immediately.

